

"LEADER IN HOSPITAL CARE PRODUCTS FOR OVER 60 YEARS"



***SSCI***

***SUBURBAN SURGICAL Co., Inc.***

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.  
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

---

***For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061***

***Hours: 8:30am - 5:00pm CST Monday - Friday***

Visit our website at: [www.suburbansurgical.com](http://www.suburbansurgical.com)

# STAINLESS STEEL CARTS

## PRECISION CRAFTED STAINLESS STEEL CARTS DESIGNED FOR DURABILITY AND UTILITY

### Door

Removable double-panel stainless steel doors are designed for easy care: sloped top edges allow water to run off; bottom weep holes ensure effective drainage.

Doors open a full 270 degrees for loading and unloading Cart and latches in open position.

Door lifts off it's hinges for easy removal.

Door knob is self-latching to keep the door(s) latched in the closed position during holding or transit or latched against the cart's side(s) in the open position for loading or cleaning.

### Hinges

.3125-in. (.79cm) diameter stainless steel hinge pins.

Heavy 12-gauge stainless steel die-formed barrel hinges operate smoothly.

### Casters

Four 5" diameter Poly-Loc-Poly Wheels, plate-mounted swivel casters, copper/nickel plated (two with swivel locks). Swivel and Axle bearings have grease fittings and seals. Other casters are available.

### Bumper

Full wrap-around vinyl Bumper is mounted with brackets to provide clearance between the Bumper and Chassis providing adequate drainage. The Bumper extends beyond the external handle(s) and door(s) when they are latched against the cart's side, prevent damage to facility doorways and walls.

### Body

Handle Style A: Recessed Handle

Handle Style B: Recessed Handle installed down from the top of the Case Cart to compensate for the higher height of the cart.

Handle Style C: External Handle

Handle Style D: External Handle installed down from the top of the Case Cart to compensate for the higher height of the cart.

One piece stainless steel top is flanged down to overlap the wrapper's side and back, providing a smooth, crevice-free working surface.

The Top is strengthened with braces spot-welded to the underside allowing it to be used as an instrument table during procedures.

External handle extends the full width of the cart for easier maneuvering.

Integral door jamb on top, bottom and sides of door opening helps seal out environmental contamination.

Cart sides and back are made of a single piece of stainless steel formed with front flanges to create a vertical recessed door jamb.

Cart bottom is spot welded to sides and back for a smooth, easy to clean interior surface. Forward pitch and drain holes in the front corners ensure water can run out when cleaning.

Front flanges of the Top, Bottom and Wrapper are fully heli-arc welded for added strength; welds are ground smooth and polished to a uniform finish.

### Frame

Frame Rail members are heli-arc welded to create a one-piece frame that insures continuous, true alignment for easy, smooth maneuvering.

Frame design prevents any obstruction of water drainage while washing.

One-piece, welded, die-formed frame constructed entirely of stainless steel.

Shelf brackets are designed to eliminate any projecting shelf supports protruding into the cart's interior and allow independent adjustment for all shelf options.

The door Jamb is designed to have an overlapping surface behind the top, bottom and sides of the doors.

Latch striker plate is adjustable to insure tight-closing doors and insure smooth latching of the door.



**SUBURBAN SURGICAL** offers a selection of stainless steel carts for mobilizing supplies and equipment, or an enclosed cart for emergencies or surgical procedures. The high-quality carts are designed to be practical, functional and long-lasting. Our goal is to offer you carts to enhance productivity throughout your facility.



**Case Carts** are ideal for keeping all the supplies needed for individual procedures close at hand. Simply load it with sterile instruments, surgical drapes, sutures and other supplies you need and wheel it to the surgical suite. When the procedure is finished, used items can be placed inside for return to facility's dirty area for disposing of waste and forwarding of reusable items to facility's sanitizing area where instrument's are cleaned and autoclaved. The Case Cart is restocked with sterile instrument packs and items needed for the next surgical procedure or emergency.

The quality of SUBURBAN SURGICAL's Stainless Steel Closed Case Carts quite simply can't be beat. We start with heavy-gauge Type 304 stainless steel, formed, heli-arc welded and polished to a smooth, even finish. We then incorporate a host of construction features shown above designed for industrial-strength durability and exceptionally easy care.



## Single Door Closed Case Cart


Single door is hinged left. Full-width External Handle on right side, Recessed Handle on left side. Cart can be configured with one or two optional shelves.

P/N	DESCRIPTION Handle Style	Cabinet Size			Interior Space (Cu. Ft)	Overall Size			Required Floor Space (Sq. Ft)
		L	W	H		L	W	H	
Recessed Handle									
105487-00	A	24.867"	24.00"	33.00"	9.95	28.625"	27.375"	39.375"	5.44
105482-00	B	24.867"	24.00"	48.00"	14.86	28.625"	27.375"	54.375"	5.44
105483-00	B	24.867"	24.00"	63.00"	19.76	28.625"	27.375"	69.375"	5.44
External Handle									
105587-00	C	24.867"	24.00"	33.00"	9.95	30.625"	27.375"	39.375"	5.82
105582-00	D	24.867"	24.00"	48.00"	14.86	30.625"	27.375"	54.375"	5.82
105583-00	D	24.867"	24.00"	63.00"	19.76	30.625"	27.375"	69.375"	5.82

# Regal Closed Case Carts

## Double Doors with Overlapping Door

When doors are open, cart's full-width optional shelves have full access. Overlapping doors close tightly to help seal out environmental contaminants. When open, large, unobstructed interior space allows easier loading and unloading. Full-width external handles on both sides give these Carts excellent maneuverability. These carts are also available with Recessed Handles.

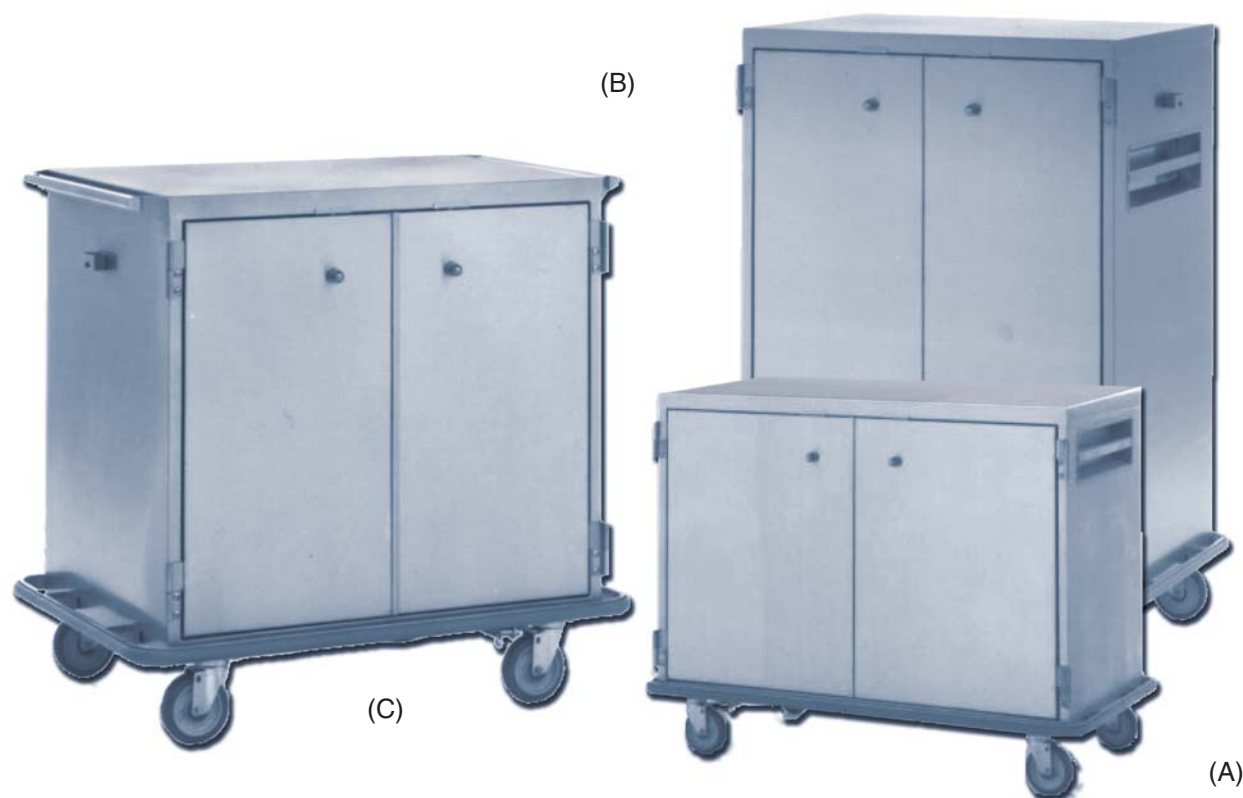

P/N	DESCRIPTION Handle Style	Cabinet Size			Interior Space (Cu. Ft)	Overall Size			Required Floor Space (Sq. Ft)
		L	W	H		L	W	H	
Recessed Handle									
105494-00	A	30"	24"	33"	11.97	33.750"	27.375"	39.375"	6.42
105502-00	A	36"	24"	33"	14.41	39.750"	27.375"	39.375"	7.56
105510-00	A	42"	24"	33"	16.86	45.750"	27.375"	39.375"	8.70
105518-00	A	48"	24"	33"	19.32	51.750"	27.375"	39.375"	9.84
105524-00	B	30"	24"	48"	17.82	35.750"	27.375"	54.375"	6.80
105532-00	B	36"	24"	48"	21.44	41.750"	27.375"	54.375"	7.94
105540-00	B	42"	24"	48"	25.07	47.750"	27.375"	54.375"	9.08
105548-00	B	48"	24"	48"	28.71	53.750"	27.375"	54.375"	10.22
External Handle									
105594-00	C	30"	24"	33"	11.97	35.750"	27.375"	39.375"	6.80
105602-00	C	36"	24"	33"	14.41	41.750"	27.375"	39.375"	7.94
105610-00	C	42"	24"	33"	16.86	47.750"	27.375"	39.375"	9.08
105618-00	C	48"	24"	33"	19.32	53.750"	27.375"	39.375"	10.22



# Regal Closed Case Carts

## Double Doors with Center Post

Allows for cart configuration with shelves both sides of center post. When doors are open, cart's optional shelves have full access. Doors close tightly against center post to help seal out environmental contaminants. Full-width external handles on both sides give these Carts excellent maneuverability.



P/N	DESCRIPTION Handle Style	Cabinet Size			Interior Space (Cu. Ft)	Overall Size			Required Floor Space (Sq. Ft)
		L	W	H		L	W	H	
Recessed Handle									
105495-00	A	30.00"	24.00"	33.00"	11.97	33.750"	27.375"	39.375"	6.42
105503-00	A	36.00"	24.00"	33.00"	14.41	39.750"	27.375"	39.375"	7.56
105511-00	A	42.00"	24.00"	33.00"	16.86	45.750"	27.375"	39.375"	8.70
105519-00	A	48.00"	24.00"	33.00"	19.32	51.750"	27.375"	39.375"	9.84
105525-00	B	30.00"	24.00"	48.00"	17.82	35.750"	27.375"	54.375"	6.80
105533-00	B	36.00"	24.00"	48.00"	21.44	41.750"	27.375"	54.375"	7.94
105541-00	B	42.00"	24.00"	48.00"	25.07	47.750"	27.375"	54.375"	9.08
105549-00	B	48.00"	24.00"	48.00"	28.71	53.750"	27.375"	54.375"	10.22
External Handle									
105595-00	C	30.00"	24.00"	33.00"	11.97	35.750"	27.375"	39.375"	6.80
105603-00	C	36.00"	24.00"	33.00"	14.41	41.750"	27.375"	39.375"	7.94
105611-00	C	42.00"	24.00"	33.00"	16.86	47.750"	27.375"	39.375"	9.08
105619-00	C	48.00"	24.00"	33.00"	19.32	53.750"	27.375"	39.375"	10.22



# Closed Case Carts

## Materials:

### Cabinet:

**Top:** 16 gauge, stainless steel, type 304, #4 finish

**Braces:** 16 gauge, stainless steel, type 304, 2B finish

**Back and Sides:** 18 gauge, stainless steel, type 304, #4 finish

### Handle Styles:

**Styles A & B:** Recessed Push Handle

**Housing:** 20 gauge, stainless steel, type 304, #4 finish

**Push Bar:** 20 gauge, stainless steel, type 304, #4 finish

**Styles C & D:** External Push Handle

**Brackets:** 11 gauge, stainless steel, type 304, #4 finish

**Push Bar:** 18 gauge, 1" dia. tubular stainless steel, type 304, 320 grit finish

**Shelf Brackets:** 18 gauge, stainless steel, type 304, 2B finish

**Bottom:** 18 gauge, stainless steel, type 304, #4 finish

### Door:

**Outer Panel:** 18 gauge, stainless steel, type 304, #4 finish

**Inner Panel:** 20 gauge, stainless steel, type 304, #4 finish

**Latch:** stainless steel, type 304, polished finish

**Latch Striker Plate:** 18 gauge, stainless steel, type 304, 2B finish

**Hinges:** 12 gauge, stainless steel, type 304, polished finish

**Hinge Pin:** .312" dia. rod stainless steel, type 304, polished finish

### Chassis:

**Frame Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Frame Cross Rail/Caster Support:** 11 gauge, stainless steel, type 304, 2B finish

**Casters:** 5" dia. Poly-Loc-Poly wheel plate mounted swivel, copper/nickel plated with grease fittings and seals, 2 with swivel locks mounted diagonally

### Bumper:

**Insert:** non-marring extruded vinyl

**Channel:** 14 gauge, roll-formed stainless steel, "C" shaped, type 304, 2B finish

**Brackets:** 14 gauge, stainless steel, type 304, 2B finish

### Fasteners:

**Machine Screws:** stainless steel, type 18-8

**Lock Nuts:** stainless steel, type 18-8, nylon insert

## Construction:

### Cabinet:

**Top:** one piece with back and both sides flanged down; front is flanged down and formed to create a sloped horizontal recessed door jamb with hemmed and flattened door stop

**Braces:** formed angles; spot-welded to the Top's underside; punched to accommodate drawer runners

**Back and Sides:** one piece wrapper; formed with 90 degree bends; the front flanges of the sides are formed to create a vertical recessed door jamb with hemmed and flattened door stop; wrapper is inserted into the Top with back and sides spot-welded to Top's flanges; Top's front flange with recessed door jamb is fully heli-arc welded to the wrapper's front flanges with recessed door jamb; welds are ground smooth and polished to a uniform finish

### Handles, Styles A & B:

**Housing:** formed cavity with tapered top and bottom and solid ends; spot-welded into Cabinet Sides

**Push Bar:** formed channel with hemmed edges; fastened into Housing with machine screws and lock nuts

### Handles, Styles C & D:

**Bracket:** shaped and angled to inset Push Bar ends from front and back of Top; fully heli-arc welded to Top's corners; welds are polished to a uniform finish

**Push Bar:** ends are fully heli-arc welded to Brackets; welds are ground smooth and polished to a uniform finish

**Shelf Brackets:** punched and formed angles; spot-welded inside Cabinet to wrapper's Back and vertical door stop

**Bottom:** one piece with back and both sides flanged down; front is formed to create a recessed door jamb and stop; Bottom is pitched towards the front door opening with drain holes punched in the front corners; Bottom is inserted into the wrapper with back and sides spot-welded to wrapper to create a smooth, edge-free cart interior; front flange with recessed door jamb is fully heli-arc welded to the wrapper's front flanges with recessed door jamb; welds are ground smooth and polished to a uniform finish; caulked with clear, temperature-resistant silicone sealant to sides and back

### Doors:

**Outer Panel:** one piece formed with top, bottom, and both sides flanged; top flange is sloped

**Inner Panel:** one piece formed with top, bottom, and both sides flanged; Inner Panel is inserted into outer panel; corners are fully heli-arc welded; welds are ground smooth and polished to a uniform finish; top, bottom, and both sides are intermediately heli-arc welded; welds are polished to a uniform finish; drainage holes are incorporated into the door's bottom edge

**Latch:** self-latching latch is incorporated into the door's design

**Latch Striker Plate:** punched and formed; fastened to recessed door jamb with machine screws and lock nuts

**Hinges:** flag-style hinge sets; male portion with pin fastened to the Cabinet with machine screws and lock nuts; female portion heli-arc welded to the door; welds are polished to a uniform finish

### Chassis:

**Frame Rails:** punched and formed angles

**Frame Cross Rail/Caster Support:** punched and formed channel fitted with weldnuts for mounting casters Frame Rails are heli-arc welded creating a one-piece frame; independent chassis is inserted into bottom of Cabinet, recessing the chassis

**Casters:** fastened to frame with machine screws

### Bumper:

**Insert:** slid into Channel

**Channel:** bends are mandrill-bent for wrinkle-free bends

**Brackets:** punched and formed; spot welded to Channel

Bumper is fastened through the Cabinet to Chassis with machine screws and lock nuts



# Accessories For Closed Case Carts



## Wire Shelves

Heavy-duty Type 304 stainless steel wire shelves feature continuous round rod construction welded at all intersections for extra strength. Front and back edges are turned down to create a one-piece shelf designed to handle heavy loads. Plus, there's no front lip which makes it easier to load and unload items on the shelf. Brackets inside cart hold shelves securely, and allow shelf height to be adjusted easily when needed.

The opening is completely recessed to prevent environmental contamination.

Shelf brackets are designed to hold shelves securely in place without projecting into cart's interior.

P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N
<b>Wire Shelves</b>		
220300-00	For Single Door or Half Shelf in Regal Case Cart with center post	105847-00, 105519-00, 105549-00, 105619-00
220301-00	For Double Door Regal Case Cart with overlap	105494-00, 105524-00, 105594-00
220302-00	For Double Door Regal Case Cart with center post	105495-00, 105525-00, 105595-00
220303-00	For Double Door Regal Case Cart with overlap	105502-00, 105532-00, 105602-00
220304-00	For Double Door Regal Case Cart with center post	105503-00, 105533-00, 105603-00
220305-00	For Double Door Regal Case Cart with overlap	105510-00, 105540-00, 105610-00
220306-00	For Double Door Regal Case Cart with center post	105511-00, 105541-00, 105611-00
220307-00	For Double Door Regal Case Cart with overlap	105518-00, 105548-00, 105618-00

## Roll-Out Wire Shelves

Heavy-duty Type 304 stainless steel shelves feature continuous round rod construction welded at all intersections for extra strength. Turned-up back edge keeps items from falling off back when shelf is rolled out, and turned-down front edge eliminates a front lip to make it easier to load and unload items. Shelf suspension tracks fit easily into shelf brackets inside cart. Shelf rolls out to 14.75-in. (37.47cm) past Cart face, and positive stop prevents accidental removal. Height of shelves may be adjusted easily.



Stainless steel shelf runner has smooth-rolling .750-in. (1.91cm) diameter polyethylene bearings.

Retaining lip on sides and back keeps supplies from sliding off shelf when pulling shelf forward.



P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N
<b>Roll-Out Wire Shelves</b>		
220330-00	For Single Door or Half Shelf in Regal Case Cart with center post	105847-00, 105519-00, 105549-00, 105619-00
220331-00	For Double Door Regal Case Cart with overlap	105494-00, 105524-00, 105594-00
220332-00	For Double Door Regal Case Cart with center post	105495-00, 105525-00, 105595-00
220333-00	For Double Door Regal Case Cart with overlap	105502-00, 105532-00, 105602-00
220334-00	For Double Door Regal Case Cart with center post	105503-00, 105533-00, 105603-00
220335-00	For Double Door Regal Case Cart with overlap	105510-00, 105540-00, 105610-00
220336-00	For Double Door Regal Case Cart with center post	105511-00, 105541-00, 105611-00
220337-00	For Double Door Regal Case Cart with overlap	105518-00, 105548-00, 105618-00

## Accessories

### Solid Shelves

Heavy-duty Type 304 stainless steel solid shelves. Front and back edges are turned down to create a one-piece shelf designed to handle heavy loads. Plus, there's no front lip which makes it easier to load and unload items on the shelf. Brackets inside cart hold shelves securely, and allow shelf height to be adjusted easily when needed.



P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N
<b>Solid Shelves</b>		
220310-00	For Single Door or Half Shelf in Regal Case Cart with center post	105847-00, 105519-00, 105549-00, 105619-00
220311-00	For Double Door Regal Case Cart with overlap	105494-00, 105524-00, 105594-00
220312-00	For Double Door Regal Case Cart with center post	105495-00, 105525-00, 105595-00
220313-00	For Double Door Regal Case Cart with overlap	105502-00, 105532-00, 105602-00
220314-00	For Double Door Regal Case Cart with center post	105503-00, 105533-00, 105603-00
220315-00	For Double Door Regal Case Cart with overlap	105510-00, 105540-00, 105610-00
220316-00	For Double Door Regal Case Cart with center post	105511-00, 105541-00, 105611-00
220317-00	For Double Door Regal Case Cart with overlap	105518-00, 105548-00, 105618-00

### Roll-out Solid Shelves

The Solid Roll-out style Shelves are a true labor-saver; the efficient design allows these shelves to be smoothly rolled out of the cart for easy loading and unloading. The continuous flow front and back edges are designed to increase the shelf strength. The back edge is flanged up to prevent stored items from falling off the back as the shelf is rolled out. The front edge is flanged down to provide a smooth, waterfall, edge-free surface for supporting the heaviest loads. The Suspension Tracks are adjustable and suspended from the cart's shelf brackets. The shelves ride on polyethylene bearings and, when rolled out, extend 14.75" past the cart's face to a positive stop preventing accidental removal.



P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N
<b>Roll-Out Solid Shelves</b>		
220340-00	For Single Door or Half Shelf in Regal Case Cart with center post	105847-00, 105519-00, 105549-00, 105619-00
220341-00	For Double Door Regal Case Cart with overlap	105494-00, 105524-00, 105594-00
220342-00	For Double Door Regal Case Cart with center post	105495-00, 105525-00, 105595-00
220343-00	For Double Door Regal Case Cart with overlap	105502-00, 105532-00, 105602-00
220344-00	For Double Door Regal Case Cart with center post	105503-00, 105533-00, 105603-00
220345-00	For Double Door Regal Case Cart with overlap	105510-00, 105540-00, 105610-00
220346-00	For Double Door Regal Case Cart with center post	105511-00, 105541-00, 105611-00
220347-00	For Double Door Regal Case Cart with overlap	105518-00, 105548-00, 105618-00





### Perforated Shelves

The Perforated Adjustable Shelves feature continuous design front and back edges which are flanged down to increase the shelf strength and provide a smooth, waterfall, edge-free surface for supporting the heaviest loads. The shelves are adjustable and secured in the cart's shelf brackets.

P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N
<b>Perforated Shelves</b>		
220320-00	For Single Door or Half Shelf in Regal Case Cart with center post	105847-00, 105519-00, 105549-00, 105619-00
220321-00	For Double Door Regal Case Cart with overlap	105494-00, 105524-00, 105594-00
220322-00	For Double Door Regal Case Cart with center post	105495-00, 105525-00, 105595-00
220323-00	For Double Door Regal Case Cart with overlap	105502-00, 105532-00, 105602-00
220324-00	For Double Door Regal Case Cart with center post	105503-00, 105533-00, 105603-00
220325-00	For Double Door Regal Case Cart with overlap	105510-00, 105540-00, 105610-00
220326-00	For Double Door Regal Case Cart with center post	105511-00, 105541-00, 105611-00
220327-00	For Double Door Regal Case Cart with overlap	105518-00, 105548-00, 105618-00

### Roll-out Perforated Shelves

The Perforated Roll-out style Shelves are a true labor-saver; the efficient design allows these shelves to be smoothly rolled out of the cart for easy loading and unloading. The continuous flow front and back edges are designed to increase the shelf strength. The back edge is flanged up to prevent stored items from falling off the back as the shelf is rolled out. The front edge is flanged down to provide a smooth, waterfall, edge-free surface for supporting the heaviest loads. The Suspension Tracks are adjustable and suspended from the cart's shelf brackets. The shelves ride on polyethylene bearings and, when rolled out, extend 14.75" past the cart's face to a positive stop preventing accidental removal.



P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N
<b>Roll-Out Perforated Shelves</b>		
220350-00	For Single Door or Half Shelf in Regal Case Cart with center post	105847-00, 105519-00, 105549-00, 105619-00
220351-00	For Double Door Regal Case Cart with overlap	105494-00, 105524-00, 105594-00
220352-00	For Double Door Regal Case Cart with center post	105495-00, 105525-00, 105595-00
220353-00	For Double Door Regal Case Cart with overlap	105502-00, 105532-00, 105602-00
220354-00	For Double Door Regal Case Cart with center post	105503-00, 105533-00, 105603-00
220355-00	For Double Door Regal Case Cart with overlap	105510-00, 105540-00, 105610-00
220356-00	For Double Door Regal Case Cart with center post	105511-00, 105541-00, 105611-00
220357-00	For Double Door Regal Case Cart with overlap	105518-00, 105548-00, 105618-00

Accessories



Top Extensions

Heavy-duty Type 304 stainless steel Top Extensions are one piece with front and back and one end hemmed. Flanged corners are fully heli-arc welded, and ground smooth, then polished to a uniform finish. Our Top Extensions are self-storing and easily extend the length of a cart's top.

P/N	DESCRIPTION	REGAL CLOSED CASE CART P/N	Cabinet Top Length	Extends Top Length, per End
Top Extensions				
220400-00	Top Extension for Carts	105518-00, 105519-00, 105548-00, 105549-00, 105847-00	24.867"	14.125"
220401-00	Top Extension for Carts	105494-00, 105495-00 , 105524-00, 105525-00	30.00"	5.125"
220402-00	Top Extension for Carts	105502-00, 105503-00, 105532-00, 105533-00	36.00"	8.125"
220403-00	Top Extension for Carts	105510-00, 105511-00, 105540-00	42.00"	11.125"
220410-00	Top Extension for Carts	105518-00, 105619-00, 105582-00, 105583-00, 105587-00	48.00"	14.125"
220411-00	Top Extension for Carts	105594-00, 105595-00	30.00"	5.125"
220412-00	Top Extension for Carts	105602-00, 105603-00	36.00"	8.125"
220413-00	Top Extension for Carts	105610-00, 105611-00	42.00"	11.125"

Card Holders and Engraved Plates

Card Holders and Engraved Plates must be specified when ordering Case Carts.

\*\* Custom-sized Card Holders available.

P/N	DESCRIPTION	SIZE	
		W	H
Card Holders**			
600013	Open Top Edge, 5"W	5"	3"
603673	Open Top Edge, 6"W	6"	4"
601495	Open Top Edge, 8"W	8"	5"
Engraved Number/Signage Plates			
617375	2 digits, .750"H	2"	1.125"
610794	"THIS END ENTERS LIFT FIRST"	6"	3"
617525	"OTHER END ENTERS LIFT FIRST"	6"	3"





## Stainless Steel Drawer Runner Sets

These Drawer Runner Sets are made of Type 304 Stainless Steel, and can be used in any of our carts.

The drawers are available for use with any of our Closed Case Carts, Open Case Carts, Exchange/Supply Systemizer Carts, or Exchange/Supply Carri Carts. ABS Drawers are available either Ribbed, with vertical channels to accommodate Divider Strips, or Non-ribbed.

P/N	DESCRIPTION
<b>Runners for ABS Drawers</b>	
220430-240	Under the Top of Closed Case Carts
220431-240	Closed Case Carts, Wire Shelves
220432-240	Closed Case Carts, Solid and Perforated Shelves
220433-240	Closed Case Carts, Roll-out Wire Shelves
<b>Runners for Fiberglass Reinforced Plastic Drawer</b>	
220435-240	Under the Top of Closed Case Carts
220436-240	Closed Case Carts, Wire Shelves
220437-240	Closed Case Carts, Solid and Perforated Shelves
220438-240	Closed Case Carts, Roll-out Wire Shelves



Fiberglass Reinforced  
Plastic Drawer



ABS Drawer, Ribbed



ABS Drawer, Non-Ribbed

## ABS and Fiberglass Reinforced Plastic Drawers

The drawers are available for use with any of our Closed Case Carts, Open Case Carts, Exchange/Supply Systemizer Carts, or Exchange/Supply Carri Carts. ABS Drawers are available either Ribbed, with vertical channels to accommodate Divider Strips, or Non-ribbed.

Call your SUBURBAN SURGICAL Representative for more information.

# TOTE-AL Open Case Carts

## Design

Our Open Case Carts are designed with a flat top and with or without a 3-sided Retaining Rail.

These designs are also available with or without External Push Handle(s).

Handle Style A: No Push Handles

Handle Style B: External Push Handles

## Body

Top is one piece with front, back, and both ends hemmed; hemmed edges are flanged down creating a unitized top and double-thick angle frame. This creates a crevice-free, smooth working surface. The top is strengthened with braces spot-welded to the underside.

Open Case Carts are available with or without a 3-sided Retaining Rail which prevents stored items from falling off the top while transporting.

Open Case Carts are available with or without External Push Handles for your convenience.



## Uprights

Uprights are fully heli-arc welded to the unitized Top and double-thick angle frame in each corner for added strength.

Uprights are fully heli-arc welded to the Chassis' notched corners for added strength.

## Bumper

Vinyl inserted, total wraparound design Bumper is mounted with brackets to provide clearance between the Bumper and Chassis providing adequate drainage. The Bumper prevents damage to facility walls.

## Frame

Shelf Supports are stainless steel rods fully heli-arc circumference-welded, spaced 4.330" apart. (Optional spacings available upon request.)

Bottom Deck is one-piece creating a smooth, crevice-free surface.

Frame Rail and Caster Support Rail members are heli-arc welded creating a one-piece frame that insures continuous, true alignment for easy, smooth maneuvering. Frame design prevents any obstruction of water drainage while washing.



**SUBURBAN SURGICAL** offers a selection of stainless steel carts for mobilizing supplies and equipment, or an enclosed cart for emergencies or surgical procedures. The high-quality carts are designed to be practical, functional and long-lasting. Our goal is to offer you carts to enhance productivity throughout your facility.

## TOTE-AL Open Case Carts

Open Case Carts are also available with a 3-sided Retaining Rail around the lower deck. Contact SUBURBAN SURGICAL for Product Numbers.



Style A No Handle



Style B External Handle

Cart Number	Handle Style	L	Top Size W	H	Shelf Positions	L	Overall Size W	H	Required Floor Space (Sq. Ft)
<b>TOTE-AL Open Case Carts</b>									
105431-00	A	27.250"	25.750"	39.375"	6	29.375"	27.875"	39.375"	5.69
105433-00	A	33.250"	25.750"	39.375"	6	35.375"	27.875"	39.375"	6.84
105435-00	A	39.250"	25.750"	39.375"	6	41.375"	27.875"	39.375"	8.01
105437-00	A	45.250"	25.750"	39.375"	6	47.375"	27.875"	39.375"	9.17
105439-00	A	51.250"	25.750"	39.375"	6	53.375"	27.875"	39.375"	10.33
105441-00	B	27.250"	25.750"	39.375"	6	31.500"	27.875"	39.375"	6.10
105443-00	B	33.250"	25.750"	39.375"	6	39.625"	27.875"	39.375"	7.67
105445-00	B	39.250"	25.750"	39.375"	6	45.625"	27.875"	39.375"	8.83
105447-00	B	45.250"	25.750"	39.375"	6	51.625"	27.875"	39.375"	10.00
105449-00	B	51.250"	25.750"	39.375"	6	57.625"	27.875"	39.375"	11.15



# TOTE-AL Open Case Carts

## TOTE-AL Open Case Carts with 3-Sided Retaining Rail

Open Case Carts are also available with a 3-sided Retaining Rail around the lower deck. Contact SUBURBAN SURGICAL for Product Numbers.



Style A No Handle



Style B External Handle

Cart Number	Handle Style	L	Top Size W	H	Shelf Positions	L	Overall Size W	H	Required Floor Space (Sq. Ft)
<b>TOTE-AL Open Case Carts with Retaining Rail</b>									
105431-01	A	27.250"	25.750"	41"	6	29.375"	27.875"	41"	5.69
105433-01	A	33.250"	25.750"	41"	6	35.375"	27.875"	41"	6.84
105435-01	A	39.250"	25.750"	41"	6	41.375"	27.875"	41"	8.01
105437-01	A	45.250"	25.750"	41"	6	47.375"	27.875"	41"	9.17
105439-01	A	51.250"	25.750"	41"	6	53.375"	27.875"	41"	10.33
105441-01	B	27.250"	25.750"	41"	6	31.500"	27.875"	41"	6.10
105443-01	B	33.250"	25.750"	41"	6	39.625"	27.875"	41"	7.67
105445-01	B	39.250"	25.750"	41"	6	45.625"	27.875"	41"	8.83
105447-01	B	45.250"	25.750"	41"	6	51.625"	27.875"	41"	10.00
105449-01	B	51.250"	25.750"	41"	6	57.625"	27.875"	41"	11.15

## Materials:

### Support Structure:

**Top:** 18 gauge, stainless steel, type 304, #4 finish

**Brace:** 18 gauge, stainless steel, type 304, #4 finish

### Handle Styles:

**Style A:** No Push Handle

**Style B:** External Push Handle

**Brackets:** 11 gauge, stainless steel, type 304, #4 finish

**Push Bar:** 18 gauge, 1" dia. tubular stainless steel, type 304, 320 grit finish

### Retaining Rail:

**-00 Models:** No Rail

**-01 Models:** .312" dia. rod stainless steel, type 304, bright buffed finish

**Uprights:** 1" dia. tubular stainless steel, type 304, 320 grit finish

**Shelf Supports:** .312" dia. rod stainless steel, type 304, bright buffed finish

### Chassis:

**Bottom Deck:** 18 gauge, stainless steel, type 304, #4 finish

**Frame Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Frame Cross Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Caster Support Rails:** 12 gauge, stainless steel, type 304, 2B finish

### Casters:

5" dia. Poly-Loc-Poly wheel plate mounted swivel, copper/nickel plated with grease fittings and seals, 2 with swivel locks mounted diagonally

### Bumper:

**Insert:** non-marring extruded vinyl

**Channel:** 14 gauge, roll-formed stainless steel, "C" shaped, type 304, 2B finish

**Brackets:** 14 gauge, stainless steel, type 304, 2B finish

### Fasteners:

**Machine Screws:** stainless steel, type 18-8

**Lock Nuts:** stainless steel, type 18-8, nylon insert

## Construction:

### Support Structure:

**Top:** one piece with front, back, and both ends hemmed; hemmed edges are flanged down creating a unitized top and double-thick angle frame; corners are fully heli-arc welded; welds are ground smooth and polished to a uniform finish

**Brace:** formed channel with hemmed edges; spot-welded to the Top's underside

### External Push Handles:

**Brackets:** shaped and angled to inset Push Bar ends from front and back of Top; fully heli-arc welded to Top's corners; welds are polished to a uniform finish

**Push Bar:** ends are fully heli-arc welded to Brackets; welds are ground smooth and polished to a uniform finish

**Retaining Rail:** one piece, 3-sided, with formed support ends; intermediate supports are fully heli-arc circumference-welded to Rail; Rail is fully heli-arc circumference-welded to Top; welds are buffed to a smooth, bright polished finish

**Uprights:** punched on 4.330" centers to accommodate the Shelf Supports

**Shelf Supports:** inserted into Uprights and fully heli-arc circumference-welded; welds are polished to a uniform finish  
Uprights are fully heli-arc welded to the unitized Top and double-thick angle frame in each corner (welds are concealed)

### Chassis:

**Bottom Deck:** one piece with front, back, and both ends flanged down; corners are notched to fit the Upright's contour

**Frame Rails:** formed angles; ends are notched to fit the Upright's contour

**Frame Cross Rails:** formed angles; ends are notched to fit the Upright's contour

**Caster Support Rails:** formed channel with flanges; flanges are punched to accommodate casters

Frame Rails, Frame Cross Rails, and Caster Support Rails are fully heli-arc welded creating a one piece frame; frame is inserted into Bottom Deck to complete Chassis; Uprights are fully heli-arc welded to Chassis' notched corners; exposed welds are polished to a uniform finish creating a cart with a one piece Support Structure and Chassis

### Casters:

fastened to Caster Support Rails using machine screws and lock nuts

### Bumper:

**Insert:** slid into Channel

**Channel:** bends are mandrill-bent for wrinkle-free bends

**Brackets:** punched and formed; spot welded to Channel

Bumper is fastened to Chassis with machine screws and lock nuts

## Finishing Accessories

### Adjustable Shelves for Open Case Carts

Our adjustable and removable Shelves are available in four styles: Solid, Solid with Retaining Rail, Wire, and Wire with Retaining Rail.



For Cart Number	Solid Shelf	Solid Shelf with 4-Sided Retaining Rail	Wire Shelf	Wire Shelf with 4-Sided Retaining Rail
<b>Adjustable Shelves for Open Case Carts</b>				
105431-00	220600-00	220610-00	220620-00	220630-00
105433-00	220601-00	220611-00	220621-00	220631-00
105435-00	220602-00	220612-00	220622-00	220632-00
105437-00	220603-00	220613-00	220623-00	220633-00
105439-00	220604-00	220614-00	220624-00	220634-00
105441-00	220600-00	220610-00	220620-00	220630-00
105443-00	220601-00	220611-00	220621-00	220631-00
105445-00	220602-00	220612-00	220622-00	220632-00
105447-00	220603-00	220613-00	220623-00	220633-00
105449-00	220604-00	220614-00	220624-00	220634-00

### Adjustable Shelves for Open Case Carts with Retaining Rail



For Cart Number	Solid Shelf	Solid Shelf with 4-Sided Retaining Rail	Wire Shelf	Wire Shelf with 4-Sided Retaining Rail
<b>Adjustable Shelves for Open Carts with Retaining Rails</b>				
105431-01	220600-00	220610-00	220620-00	220630-00
105433-01	220601-00	220611-00	220621-00	220631-00
105435-01	220602-00	220612-00	220622-00	220632-00
105437-01	220603-00	220613-00	220623-00	220633-00
105439-01	220604-00	220614-00	220624-00	220634-00
105441-01	220600-00	220610-00	220620-00	220630-00
105443-01	220601-00	220611-00	220621-00	220631-00
105445-01	220602-00	220612-00	220622-00	220632-00
105447-01	220603-00	220613-00	220623-00	220633-00
105449-01	220604-00	220614-00	220624-00	220634-00

**All of our** Open Case Cart Options and Accessories are designed to further enhance the efficiency and convenience of our carts. Each product is crafted with the same high-quality manufacturing standards as our main product line.

## Drawer Runner Sets, Stainless Steel for ABS Drawers\*

P/N	DESCRIPTION
<b>For ABS Drawers</b>	
220660-257	For under the Top of "Tote-al" Open Case Carts
220661-257	For "Tote-al" Open Case Carts, Solid Shelves
220662-257	For "Tote-al" Open Case Carts, Wire Shelves
<b>For Fiberglass Reinforced Plastic Drawers*</b>	
220665-257	For under the Top of "Tote-al" Open Case Carts
220666-257	For "Tote-al" Open Case Carts, Solid Shelves
220667-257	For "Tote-al" Open Case Carts, Wire Shelves



## Card Holders, Stainless Steel

P/N	DESCRIPTION
<b>Card Holders, Stainless Steel</b>	
105420-00	Hang-on 5"W 3"H, Open on 5" side
105421-00	Hang-on 6"W 4"H, Open on 6" side
105422-00	Hang-on 8"W 5"H, Open on 8" side



## Linen Hamper Ring, Stainless Steel

P/N	DESCRIPTION
<b>Linen Hamper Ring, Stainless Steel</b>	
221075-00	Fold-down Linen Hamper Ring for "Tote-al" Open Case Carts, fits 18" bag



## SUBURBAN SURGICAL

Exchange/Supply Carts are designed in two different styles: the Systemizer Cart is designed to support Solid or Wire Shelves, and the Carri Cart is designed to support Solid or Wire Shelves and Drawers.

### Body

Support Structure Frame Members are mandrill-bent for wrinkle-free bends.

Shelf Supports are stainless steel rods fully heli-arc circumference welded, spaced 4" apart. (Optional spacings available upon request.)

Support Structure Frame Members and Chassis are fully heli-arc welded for added strength.

### Casters

Four 5" or 8" diameter\* Poly-Loc-Poly wheels, plate-mounted swivel casters, copper/nickel plated (two with swivel locks). Swivel and Axle bearings have grease fittings and seals. (Other casters are available.)



### Frame

Frame Rail and Caster Support Rail members are heli-arc welded creating a one-piece frame that insures continuous, true alignment for easy, smooth maneuvering. Frame design prevents any obstruction of water drainage while washing.

Vinyl inserted, total wraparound design bumper is mounted with brackets to provide clearance between the Bumper and Chassis providing adequate drainage. The Bumper prevents damage to facility walls.

Bottom Deck is one-piece creating a smooth, crevice-free surface.

## Materials:

### Support Structure, Systemizer Cart:

#### Style A:

**Frame Members:** 1" dia. tubular stainless steel, type 304, 320 grit finish

**Shelf Supports:** .312" dia. rod stainless steel, type 304, bright buffed finish

#### Style B:

**Frame Members:** 1.250" dia. tubular stainless steel, type 304, 320 grit finish

**Shelf Supports:** .312" dia. rod stainless steel, type 304, bright buffed finish

### Support Structure, Carri Cart:

#### Style B:

**Frame Members:** 1.250" dia. tubular stainless steel, type 304, 320 grit finish

**Vertical Frame Member:** 1.250" dia. tubular stainless steel, type 304, 320 grit finish

**Shelf Supports:** .312" dia. rod stainless steel, type 304, bright buffed finish

**Drawer Runners:** 16 gauge, stainless steel, type 304, #4 finish

### Chassis:

**Bottom Deck:** 18 gauge, stainless steel, type 304, #4 finish

**Frame Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Frame Cross Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Caster Support Rails:** 12 gauge, stainless steel, type 304, 2B finish

### Casters:

**Style A:** 5" dia. Poly-Loc-Poly wheel plate mounted swivel, copper/nickel plated with grease fittings and seals, 2 with swivel locks mounted diagonally

**Style B:** 8" dia. Poly-Loc-Poly wheel plate mounted swivel, copper/nickel plated with grease fittings and seals, 2 with swivel locks mounted diagonally

### Bumper:

**Insert:** non-marring extruded vinyl

**Channel:** 14 gauge, roll formed stainless steel, "C" shaped, type 304, 2B finish

**Brackets:** 14 gauge, stainless steel, type 304, 2B finish


### Fasteners:

**Machine Screws:** stainless steel, type 18-8

**Lock Nuts:** stainless steel, type 18-8, nylon insert



# TOTE-AL Systemizer Carri Exchange/Supply Carts

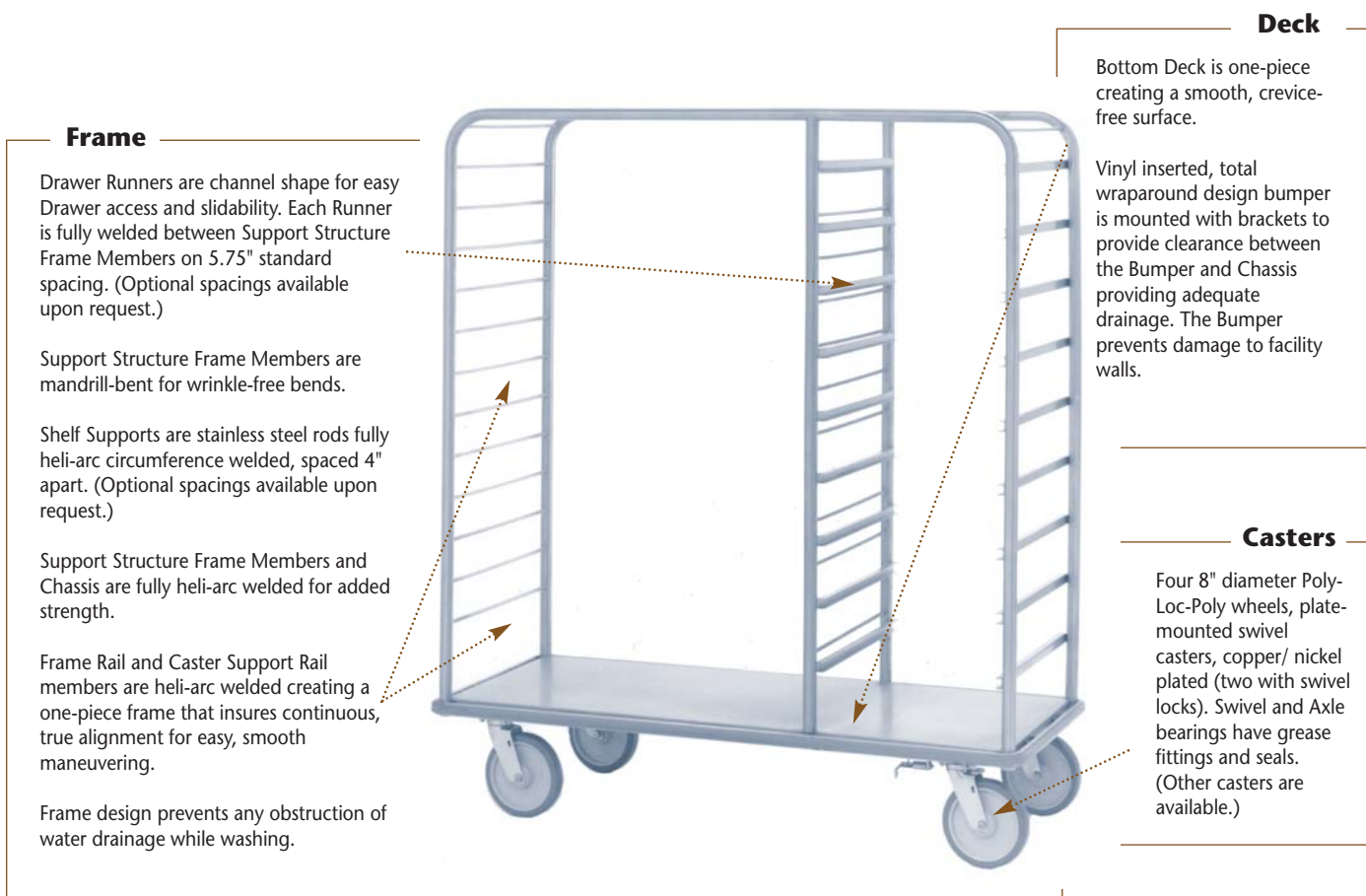

\* Cart Style A has 5" casters, Cart Style B has 8" casters

## Construction:

### Support Structure, Systemizer Cart, A & B Styles:

**Frame Members:** one piece "U"-shaped; mandrill-bent for wrinkle-free bends with holes punched on 4" centers to accommodate Shelf Supports

**Shelf Supports:** inserted into Frame Members and fully heli-arc circumference-welded; welds are polished to a uniform finish

### Support Structure, Carri Cart, B Style:

**Frame Members:** one piece "U"-shaped; mandrill-bent for wrinkle-free bends with holes punched on 4" centers to accommodate Shelf Supports

**Vertical Frame Members:** punched on 4" centers to accommodate Shelf Supports; top end coped and butted to Frame Member; bottom end fitted with a bracket; both ends are fully heli-arc welded; welds are polished to a uniform finish

**Shelf Supports:** inserted into Frame Members and Vertical Frame Members and fully heli-arc circumference-welded; welds are polished to a uniform finish

**Drawer Runners:** formed channels; fully heli-arc welded to Frame Members and Vertical Frame Members on 5.75" centers; welds are polished to a uniform finish

### Chassis:

**Bottom Deck:** one piece with front, back, and both ends flanged down; corners are notched to fit the Frame Member's contour

**Frame Rails:** formed angles; ends are notched to fit the Member's contour

**Frame Cross Rails:** formed angles; ends are notched to fit the Frame Member's contour

**Caster Support Rails:** formed channel with flanges; flanges are punched to accommodate casters

Frame Rails, Frame Cross Rails, and Caster Support Rails are fully heli-arc welded creating a one piece frame; frame is inserted into Bottom Deck to complete Chassis; Frame Members are fully heli-arc welded to Chassis' notched corners; exposed welds are polished to a uniform finish creating a cart with a one piece Support Structure and Chassis

**Casters:** are fastened to Caster Support Rails using machine screws and lock nuts

### Bumper:

**Insert:** slid into Channel

**Channel:** bends are mandrill-bent for wrinkle-free bends

**Brackets:** punched and formed; spot welded to Channel Bumper is fastened to Chassis with machine screws and lock nuts

# TOTE-AL Systemizer Carts



Cart Style A



Cart Style B

Cart Number	Cart Style	Shelf Positions	Drawer Positions	Caster Size	L	Overall Size W	H	Required Floor Space (Sq. Ft)
<b>TOTE-AL Systemizer Carts</b>								
105080-00	A	7	N/A	5"	24"	25.625"	42"	4.27
105081-00	A	10	N/A	5"	24"	25.625"	55"	4.27
105082-00	A	12	N/A	5"	24"	25.625"	64.75"	4.27
105090-00	A	7	N/A	5"	30"	25.625"	42"	5.34
105091-00	A	10	N/A	5"	30"	25.625"	55"	5.34
105092-00	A	12	N/A	5"	30"	25.625"	64.75"	5.34
105100-00	A	7	N/A	5"	36"	25.625"	42"	6.41
105101-00	A	10	N/A	5"	36"	25.625"	55"	6.41
105102-00	A	12	N/A	5"	36"	25.625"	64.75"	6.41
105104-00	B	10	N/A	8"	42"	25.625"	58.875"	7.74
105105-00	B	12	N/A	8"	42"	25.625"	68"	7.74
105110-00	B	10	N/A	8"	48"	25.625"	58.875"	8.54
105111-00	B	12	N/A	8"	48"	25.625"	68"	8.54
105120-00	B	10	N/A	8"	60"	25.625"	58.875"	10.68
105121-00	B	12	N/A	8"	60"	25.625"	68"	10.68

# Accessories For Systemizer Cart



## ABS DRAWERS

ABS Drawers are available either Ribbed, with vertical channels to accommodate Divider Strips, or Non-ribbed.



ABS Drawer, Ribbed



ABS Drawer, Non-Ribbed.

Call your SUBURBAN SURGICAL Rep for available sizes & colors.

## FIBERGLASS REINFORCED PLASTIC DRAWERS

Certain sizes of Fiberglass Reinforced Plastic Drawers are available to accommodate Divider Strips. Fiberglass Reinforced Plastic Drawers can be autoclaved at a temperature of 250° constant.



Call your SUBURBAN SURGICAL Rep for available sizes & colors.

For Cart Number	Solid Shelf	Wire Shelf	Single Retaining Rod	Double Retaining Rod	Back Panel	Covers
<b>TOTE-AL Systemizer Cart Accessories</b>						
105080-00	105168-00	105380-00	105300-00	105320-00	207230-00	105030-
105081-00	105168-00	105380-00	105300-00	105320-00	207231-00	105031-
105082-00	105168-00	105380-00	105300-00	105320-00	207232-00	105032-
105090-00	105169-00	105381-00	105301-00	105321-00	207233-00	105036-
105091-00	105169-00	105381-00	105301-00	105321-00	207234-00	105037-
105092-00	105169-00	105381-00	105301-00	105321-00	207235-00	105038-
105100-00	105170-00	105382-00	105302-00	105322-00	207236-00	105042-
105101-00	105170-00	105382-00	105302-00	105322-00	207237-00	105043-
105102-00	105170-00	105382-00	105302-00	105322-00	207238-00	105044-
105104-00	105171-00	105383-00	105303-00	105323-00	207239-00	105048-
105105-00	105171-00	105383-00	105303-00	105323-00	207240-00	105049-
105110-00	105172-00	105384-00	105304-00	105324-00	207241-00	105053-
105111-00	105172-00	105384-00	105304-00	105324-00	207242-00	105054-
105120-00	105174-00	105386-00	105306-00	105326-00	207245-00	105063-
105121-00	105174-00	105386-00	105306-00	105326-00	207246-00	105064*

See page 39 for the types of Cover Materials available.

**For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061**  
**Hours: 8:30am - 5:00pm CST Monday - Friday**  
 Visit our website at: [www.suburbansurgical.com](http://www.suburbansurgical.com)

# TOTE-AL Carri Carts




Cart Style B

Cart Number	Cart Style	Shelf Positions	Drawer Positions	Caster Positions	L	Overall Size W	H	Required Floor Space (Sq. Ft)
<b>TOTE-AL Carri Carts</b>								
105107-00	B	10	7	8"	42"	25.625"	58.875"	7.74
105108-00	B	12	9	8"	42"	25.625"	68.000"	7.74
105113-00	B	10	7	8"	48"	25.625"	58.875"	8.54
105114-00	B	12	9	8"	48"	25.625"	68.000"	8.54
105125-00	B	10	7	8"	60"	25.625"	58.875"	10.68
105126-00	B	12	9	8"	60"	25.625"	68.000"	10.68





## Accessories for Carri Cart

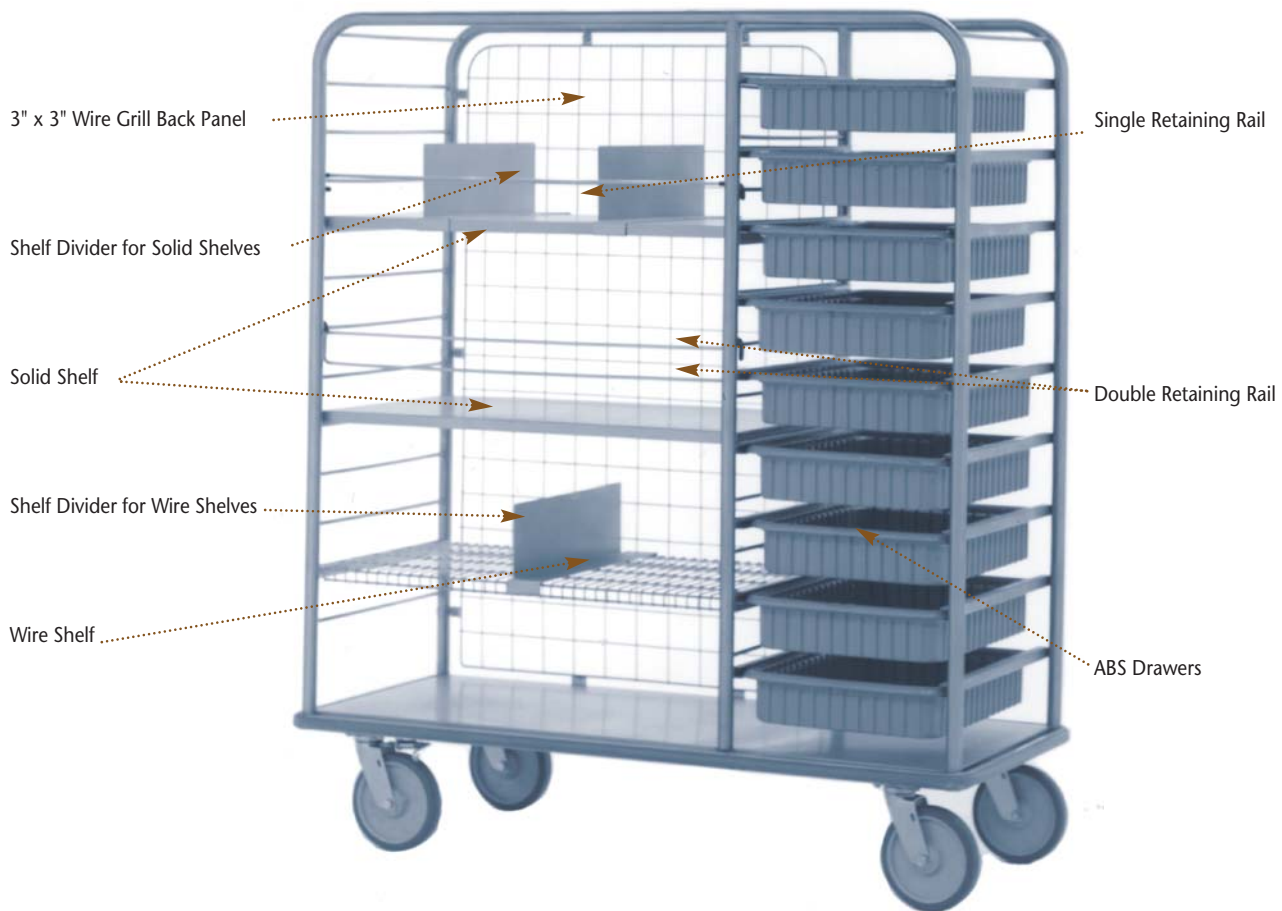

### ABS DRAWERS

ABS Drawers are available either Ribbed, with vertical channels to accommodate Divider Strips, or Non-ribbed.



ABS Drawer, Ribbed



ABS Drawer, Non-Ribbed.

Call your SUBURBAN SURGICAL Rep for available sizes & colors.

### FIBERGLASS REINFORCED PLASTIC DRAWERS

Certain sizes of Fiberglass Reinforced Plastic Drawers are available to accommodate Divider Strips. Fiberglass Reinforced Plastic Drawers can be autoclaved at a temperature of 250° constant.



Call your SUBURBAN SURGICAL Rep for available sizes & colors.

For Cart Number	Solid Shelf	Wire Shelf	Single Retaining Rod	Double Retaining Rod	Back Panel	Covers*
<b>Accessories for TOTE-AL Carri Carts</b>						
105107-00	105176-00	105387-00	105309-00	105329-00	207251-00	105048-
105108-00	105176-00	105387-00	105309-00	105329-00	207252-00	105049-
105113-00	105177-00	105388-00	105310-00	105330-00	207253-00	105053-
105114-00	105177-00	105388-00	105310-00	105330-00	207254-00	105054-
105125-00	105179-00	105390-00	105312-00	105332-00	207257-00	105063-
105126-00	105179-00	105390-00	105312-00	105332-00	207258-00	105064-

See page 25 for the types of Cover Materials available.



# TOTE-AL Linen Carts

## TOTE-AL Linen Carts

Our "Tote-al" Linen Carts feature Slide-out/Fold-down Shelves that support clean linen and then fold down to create a bin for gathering soiled linen.

### Design

Support Structure Frame Members are mandrill-bent for wrinkle-free bends.

3" x 3" Wire Grill Back and Side Panels feature stainless steel Spokes resistance-welded at each junction for strength within a stainless steel Frame. The entire Panel is electro-polished for long lasting durability.

Full-width External Push Handle is fully heli-arc welded to Support Structure Frame Members for strength and easy maneuverability.

Slide-out/Fold-down Shelf Supports are channel shape for easy Shelf slidability and fully heli-arc welded to Support Structure Frame Members.

Slide-out/Fold-down Shelves feature stainless steel Spokes resistance-welded at each junction for strength within a stainless steel Frame. The entire Shelf is electro-polished for long lasting durability.



Cart Style A

### Body

Support Structure Frame Members and Chassis are fully heli-arc welded for added strength.

Vinyl inserted, total wraparound design bumper is mounted with brackets to provide clearance between the Bumper and Chassis providing adequate drainage. The Bumper prevents damage to facility walls.

Four 5" or 8" diameter\* Poly-Loc-Poly wheels, plate-mounted swivel casters, copper/nickel plated (two with swivel locks). Swivel and Axle bearings have grease fittings and seals. (Other casters are available.)

Frame Rail and Caster Support Rail members are heli-arc welded creating a one-piece frame that insures continuous, true alignment for easy, smooth maneuvering. Frame design prevents any obstruction of water drainage while washing.

Bottom Deck is one-piece creating a smooth, crevice-free surface.

Cart Number	Cart Style	Caster Diameter	L	Overall Size W	H	Required Floor Space (Sq. Ft)	Covers
<b>TOTE-AL Linen Carts</b>							
109484-00	A	5"	36"	25.625"	64.750"	6.41	105044-
109485-00	B	8"	48"	25.625"	68"	8.54	105054-
109486-00	B	8"	60"	25.625"	68"	10.67	105064*

\* Cart Style A has 5" casters, Cart Style B has 8" casters.

**All of our** Systemizer and Carri Cart Options and Accessories are designed to further enhance the efficiency and convenience of our carts. Each product is crafted with the same high-quality manufacturing standards as our main product line.



### Drawer Runner Sets, Stainless Steel

P/N	DESCRIPTION
<b>for ABS Drawers*</b>	
221630-256	Fits Under Wire Shelves for Style A Carts
221631-256	Fits Under Wire Shelves for Style B Carts
221632-256	Fits Under Solid Shelves for Style A Carts
221633-256	Fits Under Solid Shelves for Style B Carts
<b>for Fiberglass Reinforced Plastic Drawers*</b>	
221640-256	Fits Under Wire Shelves for Style A Carts
221641-256	Fits Under Wire Shelves for Style B Carts
221642-256	Fits Under Solid Shelves for Style A Carts
221643-256	Fits Under Solid Shelves for Style B Carts

### Shelf Dividers, Stainless Steel

P/N	DESCRIPTIONS
<b>Fits Solid Shelves for Style A Carts</b>	
105358-04	4"H Divider
105358-06	6"H Divider
105358-08	8"H Divider
<b>Fits Wire Shelves for Style A Carts</b>	
105362-04	4"H Divider
105362-06	6"H Divider
105362-08	8"H Divider
<b>Fits Solid Shelves for Style B Carts</b>	
105359-04	4"H Divider
105359-06	6"H Divider
105359-08	8"H Divider
<b>Fits Wire Shelves for Style B Carts</b>	
105363-04	4"H Divider
105363-06	6"H Divider
105363-08	8"H Divider



### Card Holders, Stainless Steel

P/N	DESCRIPTION
105420-00	Hang-on 5"W 3"H, Open on 5" side
105421-00	Hang-on 6"W 4"H, Open on 6" side
105422-00	Hang-on 8"W 5"H, Open on 8" side

### Covers

Covers are custom-designed per application. They incorporate either a Velcro Closure or a Zipper Closure, information pockets, etc. Covers are made of the following materials:

Clear Vinyl	Vinyl-coated Nylon
Staph-Chek	Poly-cotton
200 Denier Nylon	400 Denier Nylon

Call your SUBURBAN SURGICAL Factory Direct Representative for more information.



## Materials:

### Support Structure:

#### Style A:

**Frame Members:** 1" dia. tubular stainless steel, type 304, 320 grit finish

**Slide-out/Fold-down Shelf Supports:** 16 gauge, stainless steel, type 304, 2B finish

#### Style B:

**Frame Members:** 1.250" dia. tubular stainless steel, type 304, 320 grit finish

**Slide-out/Fold-down Shelf Supports:** 16 gauge, stainless steel, type 304, 2B finish

#### External Push Handle:

**Brackets:** 11 gauge, stainless steel, type 304, #4 finish

**Push Bar:** 18 gauge, 1" dia. tubular stainless steel, type 304, 320 grit finish

#### Wire Grill Back and Side Panels:

**Clamps:** 16 gauge, stainless steel, type 304, 2B finish

**Frame:** .2375" dia. rod stainless steel, type 304, bright buffed finish

**Spokes:** .162" dia. rod stainless steel, type 304, bright buffed finish

**Mounting Tabs:** 12 gauge, stainless steel, type 304, 2B finish

#### Slide-out/Fold-down Shelves:

**Frame:** .312" dia. rod stainless steel, type 304, bright buffed finish

**Spokes:** .162" dia. rod stainless steel, type 304, bright buffed finish

#### Chassis:

**Bottom Deck:** 18 gauge, stainless steel, type 304, #4 finish

**Frame Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Frame Cross Rails:** 12 gauge, stainless steel, type 304, 2B finish

**Caster Support Rails:** 12 gauge, stainless steel, type 304, 2B finish

#### Casters:

**Style A:** 5" dia. Poly-Loc-Poly wheel plate mounted swivel, copper/nickel plated with grease fittings and seals, 2 with swivel locks mounted diagonally

**Style B:** 8" dia. Poly-Loc-Poly wheel plate mounted swivel, copper/nickel plated with grease fittings and seals, 2 with swivel locks mounted diagonally

#### Bumper:

**Insert:** non-marring extruded vinyl

**Channel:** 14 gauge, roll formed stainless steel, "C" shaped, type 304, 2B finish

**Brackets:** 14 gauge, stainless steel, type 304, 2B finish

#### Fasteners:

**Machine Screws:** stainless steel, type 18-8

**Lock Nuts:** stainless steel, type 18-8, nylon insert

## Construction:

### Support Structure, A & B Styles:

**Frame Members:** one piece "U"-shaped; mandrill-bent for wrinkle-free bends

**Slide-out/Fold-down Shelf Supports:** formed channels; fully heli-arc welded to Frame Members; welds are polished to a uniform finish

#### External Push Handle:

**Brackets:** shaped; fully heli-arc welded to Frame Member; welds are polished to a uniform finish

**Push Bar:** ends are fully heli-arc welded to Brackets; welds are ground smooth and polished to a uniform finish

#### Wire Grill Back and Side Panels:

**Clamps:** formed to wraparound Support Structure's Frame Members

**Frame and Spokes:** resistance welded at each junction with all ends flame-deburred

**Mounting Tabs:** fully heli-arc welded to Frame

Electro-polished finish

Fastened to the Cart's Support Structure through Mounting Tabs using Clamps, machine screws, and lock nuts

#### Slide-out/Fold-down Shelves:

**Frame and Spokes:** resistance welded at each junction with all ends flame-deburred

Electro-polished finish

### Chassis:

**Bottom Deck:** one piece with front, back, and both ends flanged down; corners are notched to fit the Frame Member's contour

**Frame Rails:** formed angles; ends are notched to fit the Frame Member's contour

**Frame Cross Rails:** formed angles; ends are notched to fit the Frame Member's contour

**Caster Support Rails:** formed channel with flanges; flanges are punched to accommodate casters

Frame Rails, Frame Cross Rails, and Caster Support Rails are fully heli-arc welded creating a one piece frame; frame is inserted into Bottom Deck to complete Chassis; Frame Members are fully heli-arc welded to Chassis' notched corners; exposed welds are polished to a uniform finish creating a cart with a one piece Support Structure and Chassis

**Casters:** fastened to Caster Support Rails using machine screws and lock nuts

### Bumper:

**Insert:** slid into Channel

**Channel:** bends are mandrill-bent for wrinkle-free bends

**Brackets:** punched and formed; spot welded to Channel

Bumper is fastened to Chassis with machine screws and lock nuts

## Stainless Steel Mobile Container Cart

features a 2.5" high raised side to prevent stored containers from falling off the Cart while transporting. The Container Cart is designed to hold containers up to 23.125"L x 14.875"W.

P/N	DESCRIPTION
108075-00	Stainless Steel Mobile Container Cart, Raised Side Supports 24"L 16"W 36"H Overall



This Stainless Steel Container Cart can hold either a Stainless Steel or ABS Container

## Materials:

### Support Structure:

#### Container Cart:

**Frame Members:** 16 gauge, 1" dia. tubular stainless steel, type 304, bright buffed finish

**Cross Members:** .375" dia. rod stainless steel, type 304, bright buffed finish

**Retaining Rail:** 11 gauge, stainless steel, type 304, #4 finish

#### Container Supports:

**Upper:** .375" dia. rod stainless steel, type 304, bright buffed finish

**Lower:** .375" dia. rod stainless steel, type 304, bright buffed finish

**Casters:** 3" dia. non-marking, swivel ball bearing, zinc plated

## Construction:

### Container Cart:

**Frame Members:** one piece "U"-shaped; bends are mandrill-bent for wrinkle-free bends

**Cross Members:** fully heli-arc circumference-welded to Frame Members; welds are buffed to a smooth, bright polished finish

**Retaining Rails:** formed ends with 1" radius bends; heli-arc welded to Frame Members; welds are buffed to a smooth, bright polished finish

#### Container Supports:

**Upper:** butted to Retaining Rails; fully heli-arc circumference-welded; welds are buffed to a smooth, bright polished finish

**Lower:** heli-arc welded to Cross Members; ends are flame-deburred; welds are buffed to a smooth, bright polished finish

**Casters:** stem-mounted into caster socket with caster socket inserted into ends of Frame Member



# Dressing Carts

## Stainless Steel Dressing Cart

Each Dressing Cart features a 36"L 20"W Top that stands 34"H from the floor and one Shelf. There is 18.75" of space between the Top and the Shelf. Both the Top and the Shelf have a 4-sided Retaining Rail.



FOR CART NUMBER	DESCRIPTION	L	OVERALL W	H
108080-00	Stainless Steel Dressing Cart	36.75"	20.75"	36"
108081-00	Stainless Steel Dressing Cart with Pail and Basin, 13 Quart Pail and 7 Quart Basin	51.50"	20.75"	36"
108081-01	Stainless Steel Dressing Cart with 13 Quart Pail, 7 Quart Basin, and Mayo Tray	51.50"	20.75"	
	Mayo Tray:	21.25"	16.25"W	
	Tray adjusts from			42.5" - 63"

## Materials:

### Chassis:

**Top and Shelf:** 18 gauge, stainless steel, type 304, #4 finish

**Brace:** 18 gauge, stainless steel, type 304, #4 finish

**Retaining Rail:** .375" dia. rod stainless steel, type 304, bright buffed finish

**Legs:** 16 gauge, 1" dia. tubular stainless steel, type 304, bright buffed finish

**Collet:** .750" I.D., stainless steel, type 304, bright buffed finish

**Knob:** triangular-shaped, black plastic, .375"-16 threads

**Pail and Basin Supports:** 11 gauge, stainless steel, type 304, #4 finish

### Tray Holder:

**Frame:** 11 gauge, stainless steel, type 304, #4 finish

**Brace:** 18 gauge, 1" square tubular stainless steel, type 304, 320 grit finish

**Telescoping Tube:** 18 gauge, .750" dia. tubular stainless steel, bright buffed finish

**Casters:** 5" Poly-Loc-Poly Wheel, Expandable Stem-Mounted Swivel Casters, Copper/Nickel plated w/Grease Fittings & Seals

thick angle frame; corners are notched to fit the leg contour; unitized top and double-thick angle frame's vertical flanges are fully heli-arc welded to the Legs (welds are concealed)

**Brace:** formed channel with hemmed edges; spot-welded to the Top's and Shelf's under sides

**Retaining Rail:** one piece; total wraparound design; fully heli-arc welded to legs; welds are buffed to a smooth, bright polished finish

### Legs:

**Collet:** machined and tapped one Leg is fitted with Collet;

Collet is fully heli-arc circumference-welded to Leg; welds are buffed to a smooth, bright polished finish; three Legs have top end capped and fully heli-arc welded; welds are ground smooth and buffed to a smooth, bright polished finish

**Knob:** threaded into Collet

**Pail and Basin Supports:** one piece formed ring with brackets; fully heli-arc welded to the Legs; welds are buffed to a smooth, bright polished finish

### Tray Holder:

**Frame:** one piece formed into oblong shape; ends are butted and fully heli-arc welded; welds are ground smooth and polished to a uniform finish

**Brace:** supports the Frame from end to end; fully heli-arc welded to the Frame; welds are buffed to a smooth, bright polished finish  
**Telescoping Tube:** fully heli-arc welded to Tray Holder Frame; welds are buffed to a smooth, bright polished finish

**Casters:** inserted into bottom ends of Legs

## Construction:


### Chassis:

**Top and Shelf:** one piece with front, back, and both ends hemmed; hemmed edges are flanged down creating a unitized top and double-



## Stainless Steel Utility/Equipment Cart

Each Utility/Equipment Cart features a 36"L x 20"W Top that stands 34"H from the floor and one or two shelves. Both the Top and Shelves have a 3-sided Retaining Lip.

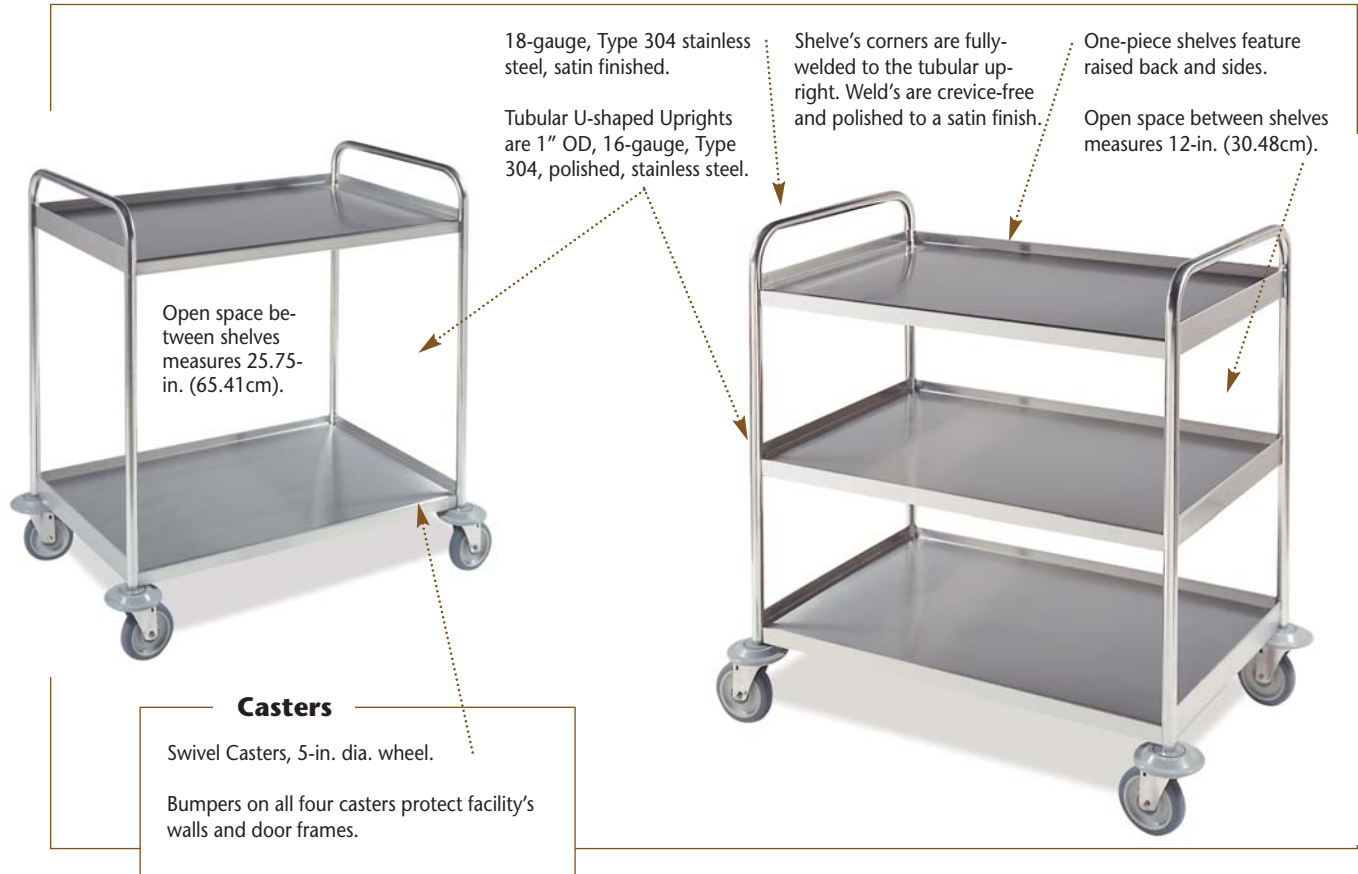

FOR CART NUMBER	DESCRIPTION	SPACE BETWEEN SHELVES	L	OVERALL W	H
108110-00	Stainless Steel Utility/Equipment Cart with Shelf	25.75" of open space between shelves	39.809"	27.779"	40"
108110-01	Stainless Steel Utility/Equipment Cart with 2 Shelves	12" of open space between shelves	39.809"	27.779"	40"

## Materials:

### Utility/Equipment Cart:

**Frame Members:** 16 gauge, 1" dia. tubular stainless steel, type 304, bright buffed finish

**Top and Shelves:** 18 gauge, stainless steel, type 304, #4 finish

**Brace:** 18 gauge, stainless steel, type 304, #4 finish

**Bumpers:** 5" dia., vinyl, rolling type

**Casters:** 5" Poly-Loc-Poly Wheel, Expandable Stem-Mounted Swivel Casters, Copper/Nickel plated w/Grease Fittings & Seals

## Construction:

### Utility/Equipment Cart:

**Frame Members:** one piece "U"-shaped; bends are mandrill-bent for wrinkle-free bends

**Top and Shelves:** one piece with retaining lip around back and both ends; front, back, and both ends hemmed; hemmed edges are flanged down creating a unitized top and double-thick angle frame; corners are notched to fit the Frame Members contour; unitized top and double-thick angle frame's vertical flanges are fully heli-arc welded to the Frame Members (welds are concealed)

**Brace:** formed channel with hemmed edges; spot-welded to the Top's and Shelf's under sides

**Bumpers:** inserted onto bottom ends of Frame Members

**Casters:** inserted into bottom ends of Frame Members