

"LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS"



SSCI

SUBURBAN SURGICAL Co., Inc.

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

STAINLESS STEEL RACKS

CLASSIC ADJUSTABLE MOBILE SHELF RACKS

Solid Shelves

SOLID Shelves are fabricated 18 gauge stainless steel with a 2B finish in various lengths and widths.

The front, back and ends of each shelf are flanged down with the lower edge flanged back, forming a channel for strength.

The front, back and ends of each shelf are flanged down with the lower edge flanged back, forming a channel for strength.

The shelves are adjustable to any position along the full height of the tubular legs.

Wire Shelves

WIRE shelves are fabricated grills made of .162" diameter stainless steel wire spokes.

.237" diameter cross members resistance welded at each junction.

Each grill has 14 gauge stainless steel die-formed channels heli-arc welded to the front, back and ends for strength.

The entire shelf is electropolished.

Uprights

Made from 16 gauge stainless steel 1.110" diameter tubular uprights'.

The upper ends are capped.

The lower ends are fitted with caster couplers.

Both ends are welded and ground smooth.

The legs are polished to a satin finish, removing any discoloration.

Casters

Chrome Plated Swivel Casters, 5-in. dia. polyolefin wheels

Four wheel brakes.

Swivel and axle bearings have seals and zerk grease fittings.



Stainless Steel Classic Adjustable Mobile Racks

These Racks are available with either solid or wire shelves. The shelves are adjustable to any position along the full height of the tubular legs.

Solid shelves are fabricated of 18 gauge stainless steel with a2B finish in various lengths and widths. The front, back and ends of each shelf are flanged down with the lower edge flanged back, forming a channel for strength.

Wire shelves are fabricated grills made of .162" diameter stainless steel wire spokes on 1" centers, with .237" diameter cross members resistance welded at each junction. Each grill has 14 gauge stainless steel die-formed channels heli-arc welded to the front, back and ends for strength. The entire shelf is electropolished.

11 gauge stainless steel corner clamps are fastened to each corner of the shelves with four fasteners, two on each side of the tubular legs. This assures a "slip-free" fit which secures the shelves in position.

Each set of 16 gauge stainless steel 1.110" diameter tubular up-rights' upper ends are capped, the lower ends are fitted with caster couplers, and both ends are welded and ground smooth. The legs are polished to a satin finish, removing any discolorations. The legs are then fitted with chrome plated swivel casters with threaded stems and 5" polyolefin wheels, four with wheel brakes. Swivel and axle bearings have seals and zerK grease fittings.

Options:

Stainless Steel Casters

Shelf Dividers



P/N	OVERALL HEIGHT
Upright Sets with 5" diameter Polyolefin Wheels, Chrome Plated Swivel Casters, 4 with Wheel Brakes	
105440-6000	66.125"
105440-6300	69.125"
105440-6600	72.125"
105440-6900	75.125"
105440-7200	78.125"
Upright Sets with 5" diameter Polyolefin Wheels, Stainless Steel Swivel Casters, 4 with Wheel Brakes	
105442-6000	66.125"
105442-6300	69.125"
105442-6600	72.125"
105442-6900	75.125"
105442-7200	78.125"
Upright Sets with 5" diameter Polyurethane Wheels, Stainless Steel Swivel Casters, 4 with Wheel Brakes	
105443-6000	66.125"
105443-6300	69.125"
105443-6600	72.125"
105443-6900	75.125"
105443-7200	78.125"

P/N SOLID SHELVES	P/N WIRE SHELVES	OVERALL SIZE	
		Length	Depth
105430-3615	105434-3615	36"	15"
105430-3618	105434-3618	36"	18"
105430-3624	105434-3624	36"	24"
105430-4815	105434-4815	48"	15"
105430-4818	105434-4818	48"	18"
105430-4824	105434-4824	48"	24"
105430-6015	105434-6015	60"	15"
105430-6018	105434-6018	60"	18"
105430-6024	105434-6024	60"	24"
105430-6030	105434-6030	60"	30"
105430-7215	105434-7215	72"	15"
105430-7218	105434-7218	72"	18"
105430-7224	105434-7224	72"	24"
105430-7230	105434-7230	72"	30"

Stainless Steel Racks

Stainless Steel Wash Racks with Fold Down Basket Shelves

The Fold Down Basket Shelves are fabricated from plates, die-formed angles and rod. The plates, angles and rods are butted at all joints and welded for strength. The rack's frame is fabricated from die-formed channels and angles. The channel uprights, angle cross members and channel cross members are butted at all joints and welded for strength.


The entire unit has a glass beaded satin finish, removing any discoloration.

The rack is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

Options:

Stainless Steel Casters

P/N	DESCRIPTION	Inches
109263	Rack	68.125"W x 32"D x 79.375"H
	w/Caster Overhang	75.625"W x 32"D x 79.375"H

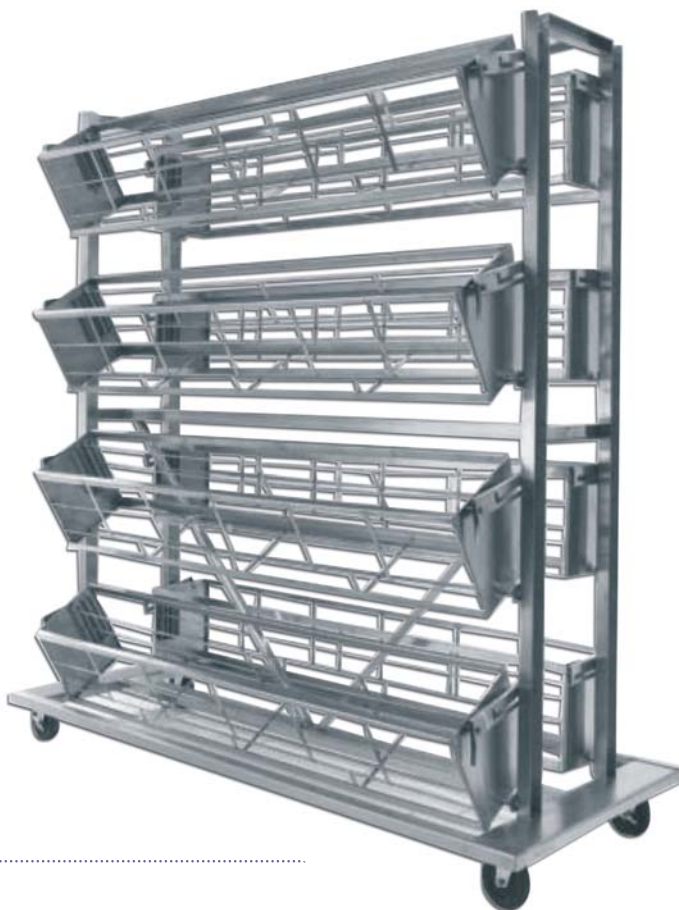

Stainless Steel Angular Racks for Washing Collection Pans

This Stainless Steel Angular Rack for Washing Collection Pans is fabricated from die-formed angles and channels. Rod dividers hold pans securely in place and allow thorough washing and drying.

The angular frame is fabricated from 11 gauge type 304 stainless steel with a 2B finish.

Each pan support is .312" diameter wire, type 304 stainless steel.

The rack has a glass beaded satin finish, removing any discoloration.

The rack is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

Options:

Stainless Steel Casters

P/N	DESCRIPTION	Inches
107363	Angular Rack for Collection Pans	68.00"L 38.00"W 53.75"H





**Stainless Steel Angular
Collection Pan Storage Racks**

The Stainless Steel Angular Collection Pan Storage Rack is fabricated from die-formed angles. The angle uprights' lower ends have caster couplers. The angle uprights and angle cross members, forming the rack's frame, are butted at all joints and welded for strength.

Each pan is held securely in place by a pair of full-length Pan Runners that are welded to the rack's uprights.

The rack's frame has a glass beaded stain finish, removing any discoloration.

The rack is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

Racks are fabricated in custom sizes to fit collection pans.

Options:

Stainless Steel Casters

P/N	DESCRIPTION
109649-00	Stainless Steel Angular Collection Pan Storage Rack

**Stainless Steel Tubular Collection Pan
Storage Racks**

The Stainless Steel Tubular Collection Pan Storage Rack support structure frame members are man drill bent round tubing; angle cross members and channel cross members are butted at all joints and welded for strength, forming the rack's frame.

Each pan is held securely in place by a pair of full-length Pan Runners that are welded to the rack's uprights.

The rack's frame has a glass beaded stain finish, removing any discoloration.

The rack is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

Racks are fabricated in custom sizes to fit collection pans.

Options:

Stainless Steel Casters

P/N	DESCRIPTION
107364-00	Stainless Steel Tubular Collection Pan Storage Rack

