

"LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS"



***SSCI***

***SUBURBAN SURGICAL Co., Inc.***

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.  
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

# STAINLESS STEEL CARTS

## PRECISION CRAFTED STAINLESS STEEL CARTS DESIGNED FOR DURABILITY AND UTILITY

**SSCI** offers a selection of stainless steel carts for mobilizing supplies and equipment, or an enclosed cart for emergencies or surgical procedures. The high-quality carts

are designed to be practical, functional and long-lasting. Our goal is to offer you carts to enhance productivity throughout your facility.

### Shelves

18-gauge, Type 304 stainless steel, satin finished.

One-piece shelves feature raised back and sides.

Shelves' corners are fully-welded to the tubular upright. Weld's are crevice-free and polished to a satin finish.

Open space between shelves measures 12-in. (30.48cm).

### Uprights

Tubular U-shaped Uprights are 1" OD, 16-gauge, Type 304, polished, stainless steel.

### Casters

Swivel Casters, 5-in. dia. wheel.

Bumpers on all four casters protect facility's walls and door frames.



## Classic Utility/Equipment Carts


Carts are designed to easily move supplies or equipment throughout your facility. Sturdy Type 304 stainless steel construction allows for sanitary food preparation and distribution. Tubular uprights also function as a push-handle for both ends of the cart. Carts roll easily on 5-in. (12.7cm) diameter wheel swivel casters to let you move supplies or equipment where you need them. Shelves measure 34-in. wide x 22.779-in. deep (86.36cmW x 57.86cmD) and feature raised sides.

Open space between shelves measures 25.75-in. (65.41cm).



P/N	DESCRIPTION	OVERALL SIZE	
		Inches	Centimeters
108110-00	Utility/Equipment Cart with Two Shelves	39.809Wx40Hx27.779D	101.11Wx101.60Hx70.56D
108110-01	Utility/Equipment Cart with Three Shelves	39.809Wx40Hx27.779D	101.11Wx101.60Hx70.56D

## ALL OF OUR REGAL CASE CARTS ARE DESIGNED WITH THESE QUALITY FEATURES:

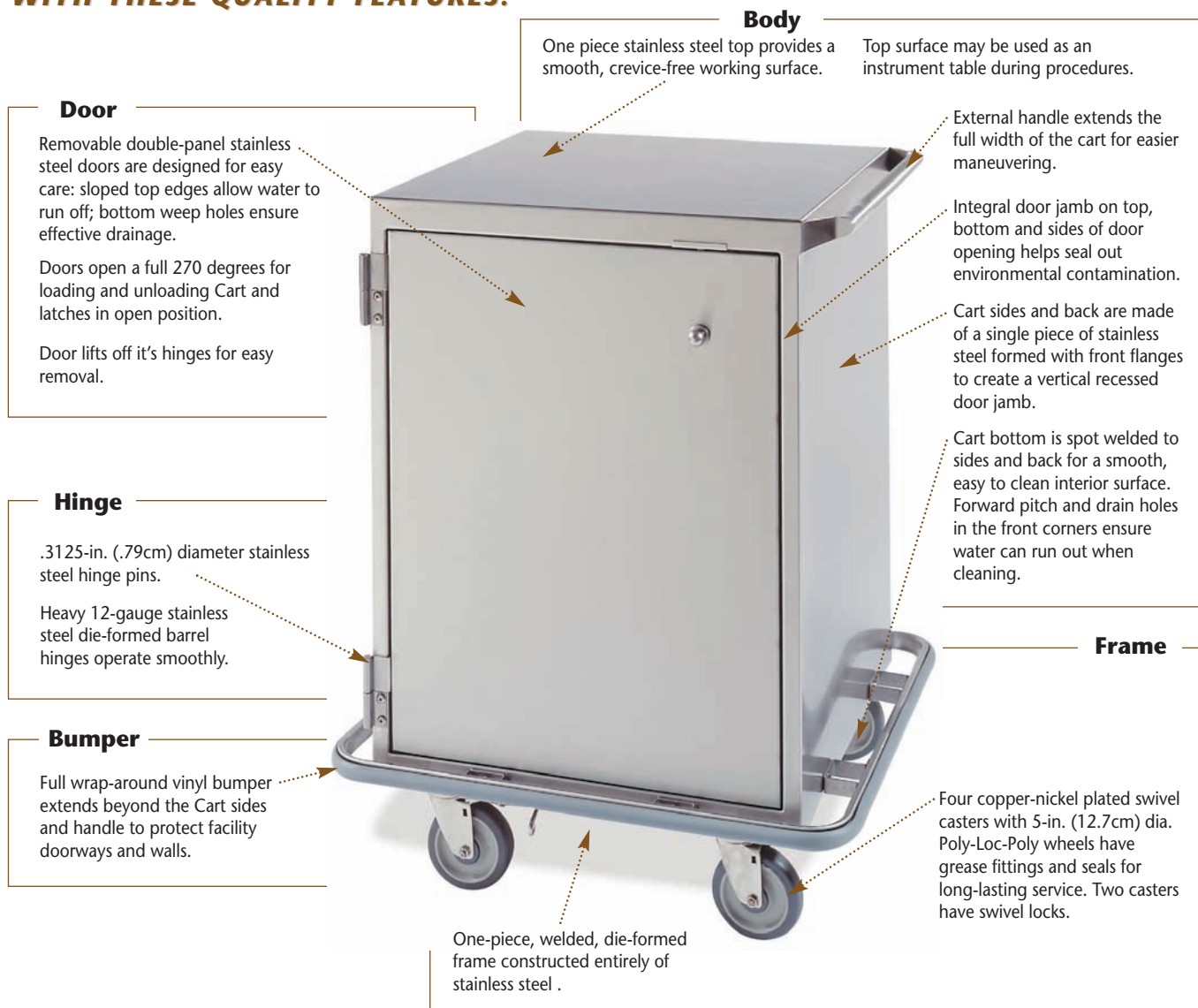

**Case Carts** are ideal for keeping all the supplies needed for individual procedures close at hand. Simply load it with sterile instruments, surgical drapes, sutures and other supplies you need and wheel it to the surgical suite. When the procedure is finished, used items can be placed inside for return to facility's dirty area for disposing of waste and forwarding of reusable items to facility's sanitizing area where instrument's are cleaned and autoclaved. The Case Cart is restocked with sterile instrument packs and items needed for the next surgical procedure or emergency.

The quality of SSCI Stainless Steel Closed Case Carts quite simply can't be beat. We start with heavy-gauge Type 304 stainless steel, formed, heli-arc welded and polished to a smooth, even finish. We then incorporate a host of construction features shown above designed for industrial-strength durability and exceptionally easy care.



# Regal Case Carts

## REGAL CASE CARTS: CHOOSE FROM THREE STYLES

### Single Door Case Cart

Single door is hinged left. Full-width external handle on right side. Cart can be configured with one or two optional shelves.

P/N	DESCRIPTION	
	Inches	Centimeters
105587-00 Cart	30.625Wx39.375Hx27.375D	77.79Wx100.01Hx69.53D
105582-00 Cart	30.625Wx54.375Hx27.375D	77.79Wx138.11Hx69.53D



### Double Doors with Overlapping Door

When doors are open, cart's full-width optional shelves have full access. Overlapping doors close tightly to help seal out environmental contaminants. When open, large, unobstructed interior space allows easier loading and unloading. Full-width external handles on both sides give these Carts excellent maneuverability.

P/N	DESCRIPTION	
	Inches	Centimeters
105602-00 Cart	41.750Wx39.375Hx27.375D	106.05Wx100.01Hx69.53D

### Double Doors with Center Post

Allows for cart configuration with shelves both sides of center post. When doors are open, cart's optional shelves have full access. Doors close tightly against center post to help seal out environmental contaminants. Full-width external handles on both sides give these Carts excellent maneuverability.

P/N	DESCRIPTION	
	Inches	Centimeters
105618-00 Cart	53.750Wx39.375Hx27.375D	136.53Wx100.01Hx69.53D





# Finishing Accessories For Regal Case Cart

## CHOOSE FROM TWO STYLES OF SHELVES TO COMPLETE YOUR CLOSED CASE CART:




Shelf brackets are designed to hold shelves securely in place without projecting into cart's interior.

### Wire Shelves

Heavy-duty Type 304 stainless steel wire shelves feature continuous round rod construction welded at all intersections for extra strength. Front and back edges are turned down to create a one-piece shelf designed to handle heavy loads. Plus, there's no front lip which makes it easier to load and unload items on the shelf. Brackets inside cart hold shelves securely, and allow shelf height to be adjusted easily when needed.

P/N	DESCRIPTION
220300-00	For Single Door Regal Case Cart
220303-00	For Double Door Regal Case Cart with overlap
220307-00	For Double Door Regal Case Cart with center post

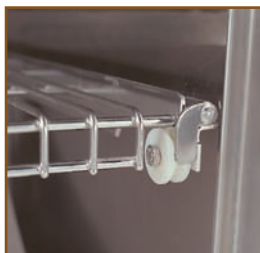

Stainless steel shelf runner has smooth-rolling .750-in. (1.91cm) diameter polyethylene bearings.



Retaining lip on sides and back keeps supplies from sliding off shelf when pulling shelf forward.

### Pull-out Wire Shelves

Heavy-duty Type 304 stainless steel shelves feature continuous round rod construction welded at all intersections for extra strength. Turned-up back edge keeps items from falling off back when shelf is rolled out, and turned-down front edge eliminates a front lip to make it easier to load and unload items. Shelf suspension tracks fit easily into shelf brackets inside cart. Shelf rolls out to 14.75-in. (37.47cm) past Cart face, and positive stop prevents accidental removal. Height of shelves may be adjusted easily.

P/N	DESCRIPTION
220330-00	For Single Door Regal Case Cart
220333-00	For Double Door Regal Case Cart with overlap
220337-00	For Double Door Regal Case Cart with center post

# Platform Carts

This cart is designed to quietly and efficiently carry loads up to 3,000 pounds.



## “TOTE-AL” Aluminum Platform Cart

The “Tote-al” Aluminum Platform Cart, with a removable stainless steel tubular handle, is designed to quietly and efficiently carry loads up to 3000 lbs.

The one-piece diamond plate deck with welded corners has sides and ends formed down into a channel and is crevice-free.

Handle sockets, cross members and reinforced caster mounts are welded for strength.

Each cart is equipped with a vinyl bumper mounted on each corner of the deck. Also included with the cart are two (2) rigid and two (2) ball-bearing swivel heavy-duty casters with polyurethane tread molded onto aluminum hubs with roller bearings. Swivel and axle bearings have seals and zerk grease fittings.

## Options:

Stainless Steel Casters


P/N	DESCRIPTION	DECK SIZE
109110	“TOTE-AL” Aluminum Platform Cart	48"L x 24"W x 10"H
109111	“TOTE-AL” Aluminum Platform Cart	60"L x 30"W x 10"H

# Classic Patient Transport Carts

**Classic** Gurneys are designed to transport patients to and from various departments within your facility. Classic Gurneys are designed to withstand the jars and jolts of everyday use to perform like new for years to come. The strong tubular U-shaped frame virtually eliminates bends and breaks so common in straight-frame rolling carts. What's more, with our one-piece welded tubular frame, all the joints and connection points are welded together for strength, there's no mechanical fasteners which could loosen or fail.

Stainless steel platform top and lower shelf are flanged down and hemmed to help eliminate any sharp edges.

Carts roll easily on Swivel Casters with 5-in. (12.70cm) diameter wheels; two have wheel brakes for safety.

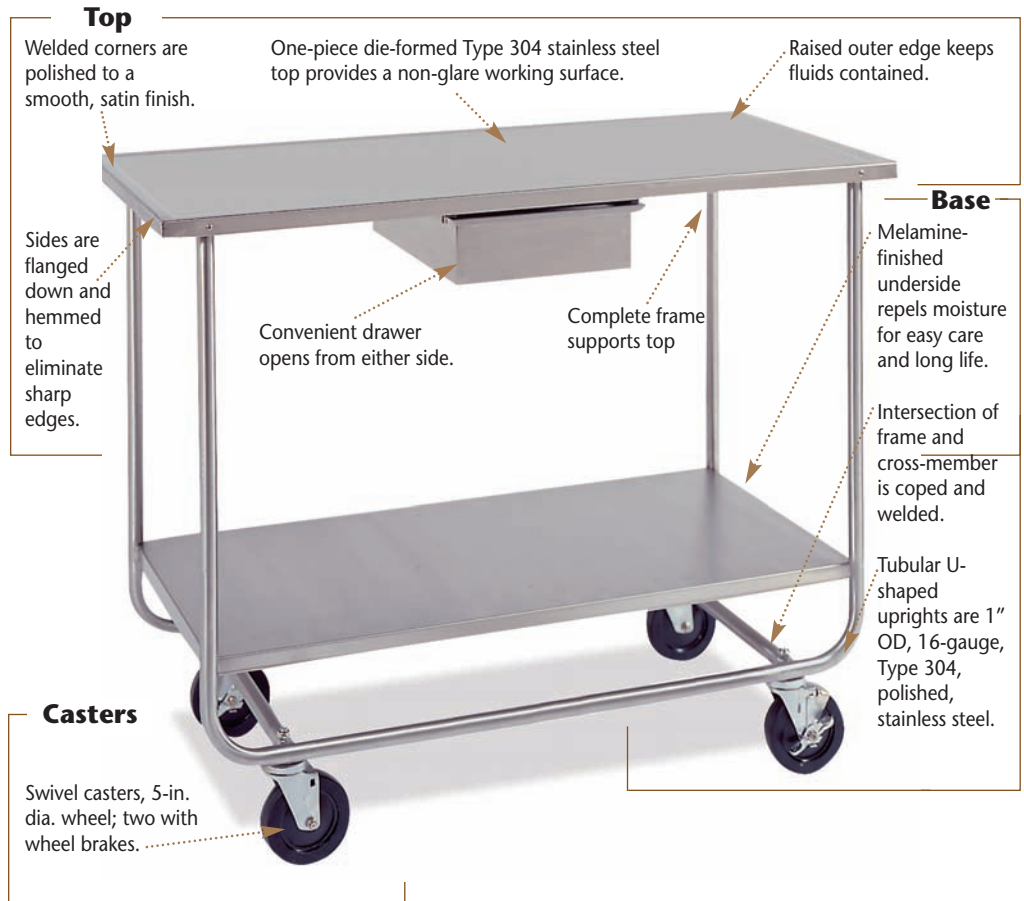

## Classic Patient Transport Cart

Stainless steel drawer slides open easily from either side to keep supplies close at hand.

P/N	DESCRIPTION	Inches	Centimeters
12511-00-GYFQDG	Classic Patient Transport Cart	45.221W x 36.513H x 21.284D	114.86W x 92.74H x 54.06D



## Classic Patient Transport Cart with Regal 300 Electronic Scale

Features an Regal 300 Electronic Platform Scale with black vinyl mat. Allows you to weigh patients anywhere in your facility. Ideal for obtaining an accurate weight of an animal just prior to surgery. See Regal 300 Electronic Scale for electronic scale information.

P/N	DESCRIPTION	Inches	Centimeters
12527-00-GZFS DH	Classic Patient Transport Cart with Regal 300 Electronic Scale	45.438W x 37.750H x 21.438D	115.41W x 95.89Hx54.37D



# Classic Vet-Mate Gurney Lift Table

**A MUST FOR EVERY VETERINARY HOSPITAL. PROTECTS PATIENTS AND YOUR STAFF BY SAFELY LIFTING AND TRANSPORTING PATIENTS THROUGHOUT THE FACILITY.**

## Lift Table • Transport Gurney • Exam/Treatment Table "ALL STAINLESS STEEL CONSTRUCTION"

### Table Top

Adjustable Tie-down system lets you safely position animals.

Stainless Steel Examination Treatment Top with Fluid Retention Lip -- 45.221" x 21.284" wide.

### THREE TABLES IN ONE:

- Lift Table
- Transport Gurney
- Exam/Treatment Table

### Frame

Fully-welded chassis and frame channel are heavy gauge tubular stainless steel --

**A SUBURBAN SURGICAL EXCLUSIVE.**

Stainless steel ball bearings roll thru the channel frame.

Bronze bushings at all pivot points ensures long wear.

5-in. Corrosion resistant locking casters for easier movement throughout your facility.

### Lift Mechanism

Powerful actuator lifts up to 200 lbs smoothly from 13" to 38" high.

Rechargeable 12 v power pack with plug-in cord and charger, gives you lifting capacity when and where you need it.

Foot operated chassis mounted UP/DOWN lift mechanism is located at the base of the table and can be mounted on either the End or the Side.





## Classic Vet-Mate Gurney Lift Table

It's a lift table, a transport gurney, and a treatment/exam table. The Classic Vet-Mate Gurney Lift Table is all this and more!

Use it to move anesthetized patients to and from holding cages, or to lift reluctant patients of just about any size up to the right height for treatment. You can even take your Vet-Mate Table outside to meet critical patients as they arrive and quickly transport them inside for treatment. Eliminates the back-breaking lifting often necessary when treating large, heavy or reluctant animals.

This versatile table lowers to just 13-in. (33.02cm) above the floor to allow you to maneuver even the heaviest patients onto it with minimal lifting. Or the Vet-Mate Table can be set at cage height, car seat height, treatment height or any height you need up to 38-in. (96.52cm) high to allow you to move patients on and off with the least amount of effort. The scissor-action electric lifting mechanism is operated by a rechargeable battery to smoothly and quietly lift the table no matter where you are inside or outside of your facility. The battery recharges quickly with the built-in charger.

Heavy-gauge tubular stainless steel channel construction stands up to rugged daily use. All joints are heli-arc welded for superior strength.




**The adjustable Classic VET-MATE Lift Gurney Table can be set anywhere from**

**13-in. (33.02cm) to 38-in. (96.52cm)**

**high to give you the optimal working height for your needs. Lowers to just 13-in. (33.02cm) for easier maneuvering and storage.**

**Foot operated chassis mounted UP/DOWN lift mechanism is located at the base of the table and can be mounted on either the End or the Side.**

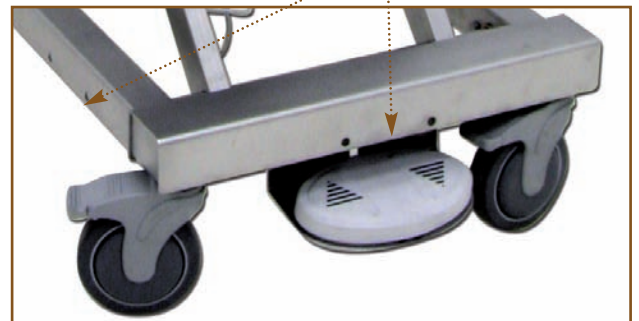

P/N	DESCRIPTION	VOLTS	OVERALL SIZE	LOWEST HEIGHT:	MAXIMUM HEIGHT:
12761-01-IGFVDL	Classic VET-MATE Lift Gurney	120V AC▲	53.63-in. wide x 22.5-in. deep (136.22 cmH x 57.15 cmW)	13-in. (33.02cm)	38-in. high (96.52cm)

▲ 220V requires step-down transformer. Contact SSCI.

# Regal Vet-Mate Gurney Lift Table

**A MUST FOR EVERY VETERINARY HOSPITAL. PROTECTS PATIENTS AND YOUR STAFF BY SAFELY LIFTING AND TRANSPORTING PATIENTS THROUGHOUT THE FACILITY.**

## Table

Adjustable Tie-down system lets you safely position animals.

Satin-finish Table Top surface is easy on the eyes even under bright surgery lights.

One piece, die-formed Type 304 18-gauge stainless steel Flat Operating Table Top is designed for long life.

Drain trough with no-spill raised outer lip channels fluids to the collection hole.

## FOUR TABLES IN ONE:

- Lift Table
- Transport Gurney
- Exam/Treatment Table
- Surgery Table

## Lift Mechanism

Scissor-action lifting mechanism has smooth, quiet bronze bushings for long-lasting service.

12 v. 26 amp rechargeable battery gives you lifting capacity when and where you need it.

12 v. DC motor-driven screw lift mechanism will raise and lower animals up to 300 lbs. (136 kgs.). Scissors Housing is a one-piece, die-formed 10-ga. stainless steel.

Scissor Frame is a one-piece, fully-welded 1-1/2" x 2", 7-ga. stainless steel tubing.

Top supported by one-piece, fully-welded frame of 2" square, 7-ga. stainless steel tubing.

Scissor-action is supported by bearings riding in a milled groove in the base.

## Base

Angle joints are cut, mitered and welded.

Casters recess under chassis, out of the way for safety.

Base is a one-piece, fully-welded frame of 2" square, 7-ga. stainless steel tubing.

1-pc., fully-welded, tubular stainless steel frame with a 7-gauge wall thickness has a uniform, non-glare finish.

Swivel Casters feature 5-in. dia. wheels and wheel brakes to keep table in place when needed.

Smooth-rolling, non-marring 5-in. (12.70cm) casters rotate 360° for incredible maneuverability.



# Regal Vet-Mate Gurney Lift Table

## It's a transport table,

a prep table, a treatment table, a lift table, a gurney, even an operating table. The Regal Vet-Mate Gurney Lift Table is all this and more!

Use it to move anesthetized patients to and from holding cages, or to lift reluctant patients of just about any size up to the right height for treatment. You can even take your Vet-Mate Table outside to meet critical patients as they arrive and quickly transport them inside for treatment. Eliminates the back-breaking lifting often necessary when treating large, heavy or reluctant animals.

This versatile table lowers to just 10.5-in. (26.67cm) above the floor to allow you to maneuver even the heaviest patients onto it with minimal lifting. Or the Vet-Mate Table can be set at cage height, car seat height, treatment height or any height you need up to 42-in. (106.68cm) high to allow you to move patients on and off with the least amount of effort. The scissor-action electric lifting mechanism is operated by a rechargeable battery to smoothly and quietly lift the table no matter where you are inside or outside of your facility. The battery recharges quickly with the built-in charger.

Heavy-gauge tubular stainless steel construction stands up to rugged daily use. All joints are heli-arc welded for superior strength.



**Infinitely adjustable Vet-Mate Gurney Table can be set anywhere from 10.5-in. (26.67cm) to 42-in. (106.68cm) high to give you the optimal working height for your needs. Lowers to just 10.5-in. (26.67cm) for easier maneuvering and storage.**



**Lift mechanism operates by convenient switches located at both ends of the table top.**

P/N	DESCRIPTION	OVERALL SIZE	LOWEST HEIGHT:	MAXIMUM HEIGHT:
12760-00-ITGLDK	Regal Vet-Mate Gurney Lift Table	57-in. wide x 22.25-in. deep (144.78cmH x 56.51cmW)	10.5-in. (26.67cm).	42-in. high (106.68cm).
<b>Accessories</b>				
854480	12 volt Re-chargeable Battery			



# Bulk Cart

## STAINLESS STEEL BULK CARTS: CHOOSE FROM A VARIETY OF ACCESSORIES

### FRAME

Type 304 Stainless Steel, 1.25" square tube frame with 3/8" rods.

Rated to hold up to 1500 pounds.

Heavy duty construction for durability.

Enclosed on 3 sides by heavy duty stainless steel rods.

Solid Metal Floor are Sloped to the back for drainage.

Front columns have hooks for optional full or half gates.

6" Caster with brakes.

Autoclavable.

### GRILL

Ideal for keeping handling cages and bulky packages on the cart as you travel from one area to another.

#### GRILL OPTIONS

Full Front Grill  
Half Front Grill

### SHELVES

Ideal for handling cages and bulky packages.

Large cage stacking surface.

Smooth surfaces for secure sterilization.

#### SHELF OPTIONS

- Solid Shelf  
- Wire Shelf





## A MUST FOR EVERY FACILITY.

**PROTECTS PATIENTS AND YOUR STAFF BY SAFELY  
LIFTING AND TRANSPORTING ITEMS THROUGHOUT THE FACILITY.**

**SSCI** offers a selection of stainless steel carts for mobilizing supplies and equipment. These high-quality carts are designed to be practical, functional and long-lasting. Our goal is to offer you carts to enhance productivity throughout your facility.

Our Bulk Carts are made this same way. Type 304 stainless steel, all welded, modular construction allows for lightweight functionality and maneuverability. Access is limited to one side, adjustable shelving is standard, and the polished finish provides aesthetic appeal. They are easy to clean and decontaminate and protect patients and staff.

### Standard Features:

- Rated to hold up to 1500 pounds
- Heavy duty construction for durability
- Type 304 Stainless Steel, 1.25" square tube frame with 3/8" rods
- Solid Metal Floor, Optional Solid Metal Shelves, Sloped to the back for drainage
- Smooth surfaces for secure sterilization
- Large cage stacking surface
- 6" Caster with brakes
- Autoclavable

P/N	DESCRIPTION	OVERALL SIZE
58000-00-HJKHDR	Stainless Steel Bulk Cart HD - Solid Bottom	48W X 67H X 24D
58000-00-HJKHEX	Stainless Steel Bulk Cart HD - Solid Bottom	48W X 67H X 32D
58000-00-JFKHDR	Stainless Steel Bulk Cart HD - Solid Bottom	60W X 67H X 24D
58000-00-JFKHEX	Stainless Steel Bulk Cart HD - Solid Bottom	60W X 67H X 32D
58001-00-HJKHDR	Stainless Steel Bulk Cart HD - Wire Bottom	48W X 67H X 24D
58001-00-HJKHEX	Stainless Steel Bulk Cart HD - Wire Bottom	48W X 67H X 32D
58001-00-JFKHDR	Stainless Steel Bulk Cart HD - Wire Bottom	60W X 67H X 24D
58001-00-JFKHEX	Stainless Steel Bulk Cart HD - Wire Bottom	60W X 67H X 32D

P/N	DESCRIPTION	OVERALL SIZE
<b>Accessories</b>		
58025-00-HJAADR	Stainless Steel Shelf - Wire HD	48W X 24D
58025-00-HJAAEX	Stainless Steel Shelf - Wire HD	48W X 32D
58025-00-JFAADR	Stainless Steel Shelf - Wire HD	60W X 24D
58025-00-JFAAEX	Stainless Steel Shelf - Wire HD	60W X 32D
58026-00-HJAADR	Stainless Steel Shelf - Wire STD	48W X 24D
58026-00-HJAAEX	Stainless Steel Shelf - Wire STD	48W X 32D
58026-00-JFAADR	Stainless Steel Shelf - Wire STD	60W X 24D
58026-00-JFAAEX	Stainless Steel Shelf - Wire STD	60W X 32D
58040-00-HJAADR	Stainless Steel Shelf - Solid	48W X 24D
58040-00-HJAAEX	Stainless Steel Shelf - Solid	48W X 32D
58040-00-JFAADR	Stainless Steel Shelf - Solid	60W X 24D
58040-00-JFAAEX	Stainless Steel Shelf - Solid	60W X 32D
58055-00-HJKHAA	Stainless Steel Front Grill - Full	48W X 67D
58055-00-JFKHAA	Stainless Steel Front Grill - Full	60W X 67D
58056-00-HJKHAA	Stainless Steel Front Grill - Half	48W X 67D
58056-00-JFKHAA	Stainless Steel Front Grill - Half	60W X 67D