

"LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS"



***SSCI***

***SUBURBAN SURGICAL Co., Inc.***

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.  
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

***For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061***

***Hours: 8:30am - 5:00pm CST Monday - Friday***

Visit our website at: [www.suburbansurgical.com](http://www.suburbansurgical.com)

Catalog No. 702488 01/04/06


# SURGEON SCRUB SINKS

## ***SURGEON SCRUB SINKS HAVE THESE QUALITY FEATURES:***

**Our** 16-gauge, Type 304 Stainless Steel Surgeon Scrub Sinks are fully heli-arc welded into one piece to create the strongest Surgeon Scrub Sink available. Welds are continuous, crevice-free and ground smooth. Coved corners eliminate areas where dirt and debris accumulate. Radius bends on bottom and all sides and a uniform, non-glare finish results in a tub interior that is

easy to clean. Scrub Sinks feature a 10.50-in. (26.67cm) deep basin, a 12-in. (30.48cm) high back splash. The Surgeon Scrub Sink is designed to accept a back splash-mounted goose neck faucet and 3.5-in. (8.89cm) drain. Stainless Steel Wall-Mounted Brackets are included.

### **One-Piece Welded Surgeon scrub Sink**

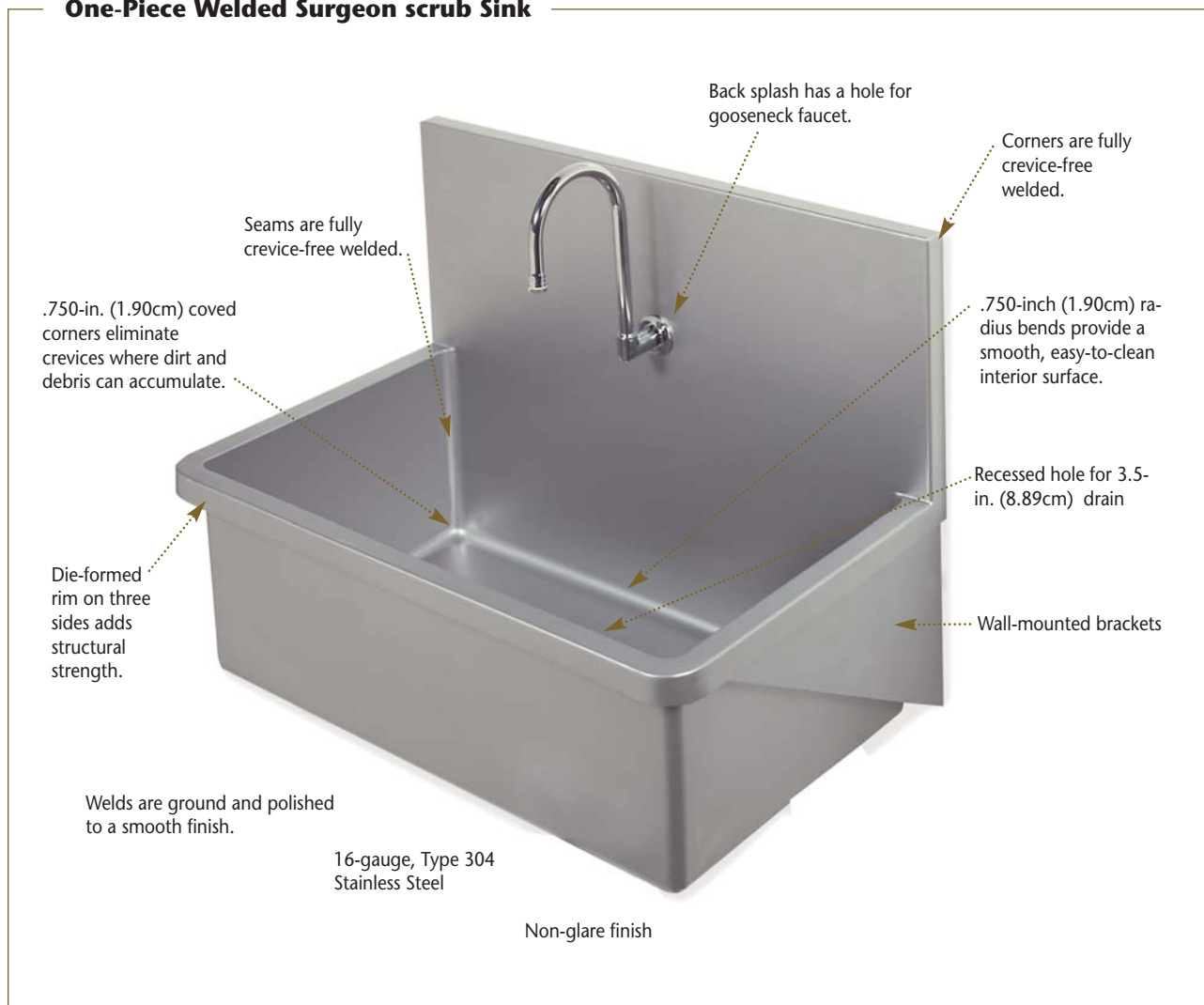

**A clean** surgical scrub area is an essential part of sanitary, bacteria-free procedures. We make it easy for you to scrub up with these top-quality all stainless steel Surgeon Scrub Sinks.

Both sizes of Surgeon Scrub Sinks are designed to accept a gooseneck faucet. Faucets can be controlled by three types of hands-free valves; foot-operated, knee-operated and electronic-motion sensor (see Section 18).



Single-Station



Dual-Station

Surgeon Scrub Sinks

P/N	Description	OVERALL SIZE	
		Inches	Centimeters
Surgeon Scrub Sinks			
13320-00-EPDMCX	Single-Station, Gooseneck Compatible	30W x 22.56H x 19D	76.20W x 57.30H x 48.26D
13320-01-EPDMCX	Single-Station, Gooseneck & Eyewash Compatible	30W x 22.56H x 19D	76.20W x 57.30H x 48.26D
13321-00-HJDMCX	Dual-Station, Gooseneck Compatible	48W x 22.56H x 19D	121.92W x 57.30H x 48.26D
13321-01-HJDMCX	Dual-Station, Gooseneck & Eyewash Compatible	48W x 22.56H x 19D	121.92W x 57.30H x 48.26D