

“LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS”



**SSCI**

**SUBURBAN SURGICAL Co., INC.**

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI. Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

**Catalog No. 702488-0308**

**For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061**  
**Hours: 8:30am - 5:00pm CST Monday - Friday**  
Visit our website at: [www.suburbansurgical.com](http://www.suburbansurgical.com)

# ELECTRONIC SCALES

## REGAL 300: YOU WON'T FIND A MORE ACCURATE SCALE!

Console Bracket mounts easily on vertical and horizontal surfaces.

Platform's 12-ft. (3.66m) long data cable connects console to platform.

Console Mounting Column is available as an accessory.

One-piece die-formed Type 304 stainless steel top provides a non-glare working surface.

The data from each of the 4 load cells is transferred to a summing panel and uploaded to the console.

Four load cells (one to each corner) for a more accurate weight reading for all size animals.

### Console

Easy-to-read LCD displays weights to the nearest 0.1 lb. or kg.

Display indicates whether pounds or kilograms has been selected.

Display console tilts for the best viewing angle.

Select unit of measure between kilograms or pounds.

Zeros out weight and turns off scale.

HOLD button locks in the weight for restless animals.

Console Mounting Bracket

TARE button compensates for the weight of a carrier. Simply weigh the carrier first and push the TARE button to reset the console to zero. Then place the animal and carrier on the scale to find the true weight of the animal.

Turns on scale and prints weight.

Welded corners are polished to a smooth, satin finish.

Melamine-finished underside repels moisture for easy care and long life.

Stable, secure platform remains steady for accurate readings even with anxious, fidgety animals.

Heavy, reinforced, 1-pc., fully-welded frame creates the strongest, most rigid platform scale available.

Sides are flanged down and hemmed to eliminate sharp edges.

Raised outer edge keeps fluids contained.

Non-skid rubber mat removes easily for cleaning.

### THE REGAL 300 ELECTRONIC SCALE IS AVAILABLE ON A WIDE VARIETY OF TABLES:



Regal Wall-Mount Lateral and Longitudinal Lift Tables  
Regal Floor-standing Lift Tables are also available.



Elite Floor-standing Lateral and Longitudinal Lift Tables  
Elite Wall-Mount Lift Tables are also available.



Premier Exam Table

**Accuracy.** It's the first thing to look for in a scale. So much depends on knowing the accurate weight of your patients: anesthesia, medications, inoculations, even feeding amounts. That's why we've made sure you can depend on the accuracy of our Regal 300 Electronic Scale.

A reinforced, 1-piece, fully-welded stainless steel-covered frame creates the strongest, most rigid platform scale available.

Stainless Steel motionless platform is supported by a load cell in each corner, with self-aligning floor pads to insure accuracy on most floors. State-of-the-art electronics take load readings from all four corner-mounted load cells to get the most accurate weight even if the animal is placed off center. Platform's 12-ft. (3.66m) long data cable plugs into the console. Many less expensive scales take readings from only one or two load cells and the platform is not motionless.

Hold button on the console allows you to lock in an electronically-averaged weight for restless animals to keep the reading from "jumping around." The Regal 300 also features a Tare function to compensate for the extra weight when an animal is weighed inside a carrier or held by its owner.

Precision Regal 300 display console can measure weights in pounds or kilograms with the push of a button, has a 300 lb. (136 kg) capacity. The easy-to-read LCD screen displays weights to the nearest 0.1 lbs. or 0.1 kg. with  $\pm 1\%$  accuracy. The scale is touch-activated with auto shut-off. The mounting bracket allows the console to be mounted to a horizontal or vertical surface and tilted to make LCD Screen easier to see.

The Regal 300 Scale is truly portable while operating with its' built-in 6-v. rechargeable battery or it can continuously powered by leaving the battery's charger plugged in.



**Pedestal Exam/  
Treatment Table**




**Transport/Prep Cart**

### Regal 300 Electronic Stainless Steel Platform Scale

Non-slip rubber mat gives anxious animals better traction, and removes easily for cleaning. Display console may be mounted on the optional mounting column, on a wall or cabinet above the scale, or in a remote location. A 12-foot (3.66m) long cable is included.

**Console Bracket mounts easily on vertical or horizontal surfaces. Platform's 12-ft. (3.66m) long data cable allows positioning Console in remote location.**

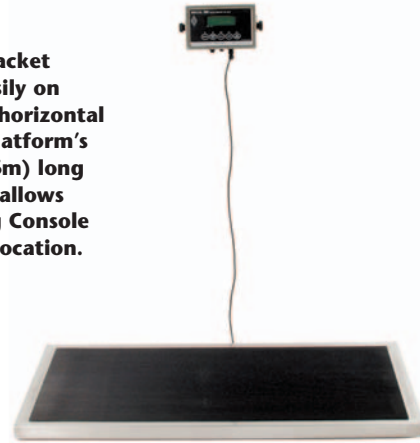

P/N	Description	Volts	Inches	DIMENSIONS	
				Inches	Centimeters
12452-00-GZAHDH	Floor Model	120V AC ▲	45.438W x 2H x 21.438D	115.41W x 5.08H x 54.45D	
<b>Accessories</b>					
12470-00-AJEPAA	Mounting Column		2.5W x 30H	6.35W x 76.20H	

▲ 220V requires step-down transformer. Contact SSCI.

### Auto 44 Electronic Stainless Steel Table Scale

With a capacity of 44 lbs (20 kg) and readability of .5oz (10g), the Auto 44 scale will accurately weigh small animals. The scale can operate on AC power (adapter included) or a single 9 volt DC battery. Hold button on the console allows you to lock in an electronically-averaged weight for restless animals to keep the reading from "jumping around."



P/N	Description	Volts	Inches	DIMENSIONS	
				Inches	Centimeters
12455-00-DDBTAT	Table Scale	120V AC ▲	20.5W x 4.875H x 11.5D	52.07W x 12.388H x 29.21D	

▲ 220V requires step-down transformer. Contact SSCI.