

“LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS”



SSCI

SUBURBAN SURGICAL Co., Inc.

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

STAINLESS STEEL GURNEYS

SSCI OFFERS SEVERAL STYLES OF GURNEYS DESIGNED TO MEET ALL OF YOUR FACILITY'S NEEDS:



Regal Vet-Mate Lift-Gurney



Classic Vet-Mate Lift-Gurney



Classic Gurney




**Regal Transport Cage,
see Section 1**

Classic Gurneys are designed to transport patients to and from various departments within your facility. Classic Gurneys are designed to withstand the jars and jolts of everyday use to perform like new for years to come. The strong tubular U-shaped frame virtually eliminates bends and breaks so common in straight-frame rolling carts. What's more, with our one-piece welded tubular frame, all the joints and connection points are welded together for strength, there's no mechanical fasteners which could loosen or fail.

Stainless steel platform top and lower shelf are flanged down and hemmed to help eliminate any sharp edges.

Carts roll easily on Swivel Casters with 5-in. (12.70cm) diameter wheels; two have wheel brakes for safety.

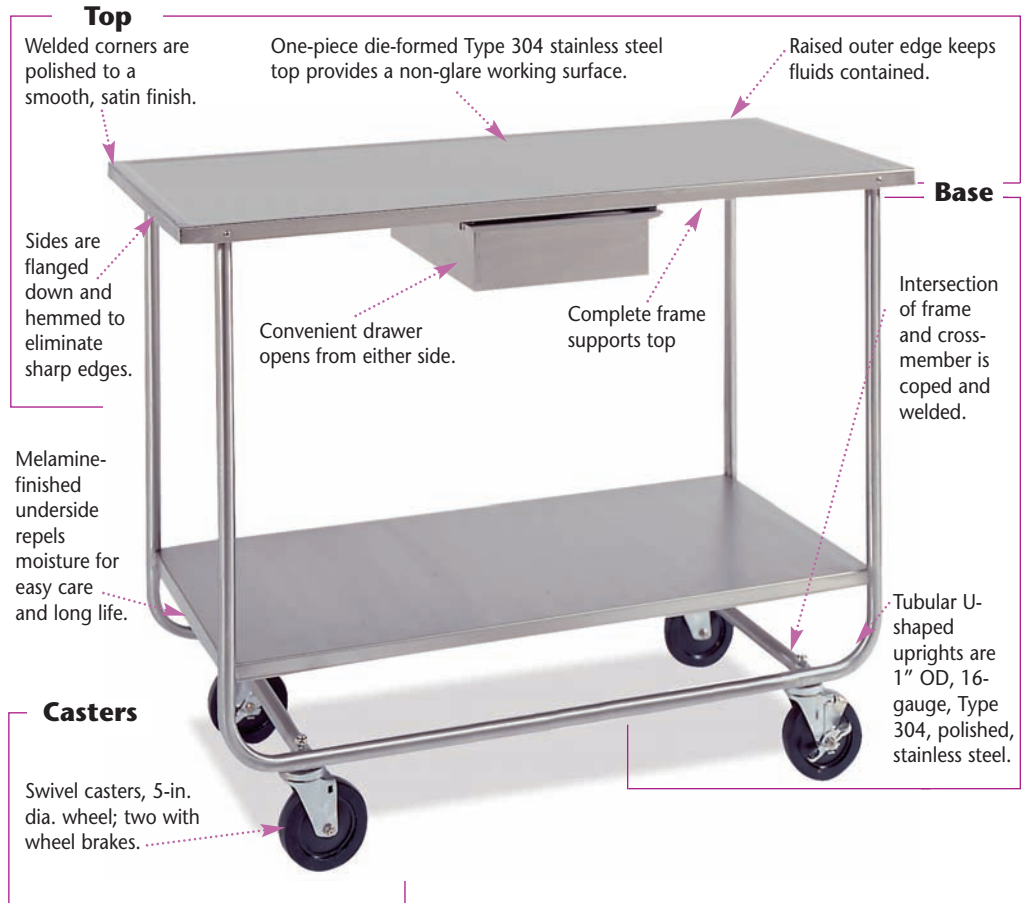

Classic Patient Transport Cart

Stainless steel drawer slides open easily from either side to keep supplies close at hand.

P/N	DESCRIPTION	
	Inches	Centimeters
12511-00-GYFQDG	Classic Patient Transport Cart 45.221W x 36.513H x 21.284D	114.86W x 92.74H x 54.06D



Classic Patient Transport Cart with Regal 300 Electronic Scale

Features an Regal 300 Electronic Platform Scale with black vinyl mat. Allows you to weigh patients anywhere in your facility. Ideal for obtaining an accurate weight of an animal just prior to surgery. See Regal 300 Electronic Scale for electronic scale information.

P/N	DESCRIPTION	
	Inches	Centimeters
12527-00-GZFS DH	Classic Patient Transport Cart with Regal 300 Electronic Scale* 45.438W x 37.750H x 21.438D	115.41W x 95.89Hx54.37D

*See Electronic Scales Section for information on the Regal 300

CLASSIC VET MATE LIFT GURNEY

A MUST FOR EVERY VETERINARY HOSPITAL. PROTECTS PATIENTS AND YOUR STAFF BY SAFELY LIFTING AND TRANSPORTING PATIENTS THROUGHOUT THE FACILITY.

Lift Table • Transport Gurney • Exam/Treatment Table
"ALL STAINLESS STEEL CONSTRUCTION"

Table Top

Adjustable Tie-down system lets you safely position animals.

Stainless Steel Examination Treatment Top with Fluid Retention Lip --
45.221" x 21.284" wide.

THREE TABLES IN ONE:

- Lift Table
- Transport Gurney
- Exam/Treatment Table

Frame

Fully-welded chassis and frame channel are heavy gauge tubular stainless steel --

A SUBURBAN SURGICAL EXCLUSIVE.

Stainless steel ball bearings roll thru the channel frame.

Bronze bushings at all pivot points ensures long wear.

5-in. Corrosion resistant locking casters for easier movement throughout your facility.

Lift Mechanism

Powerful actuator lifts up to 200 lbs smoothly from 13" to 38" high.

Rechargeable 12 v power pack with plug-in cord and charger, gives you lifting capacity when and where you need it.

Foot operated chassis mounted UP/DOWN lift mechanism is located at the base of the table and can be mounted on either the End or the Side.



Classic Vet-Mate Gurney Lift Table

It's a lift table, a transport gurney, and a treatment/exam table. The Classic Vet-Mate Gurney Lift Table is all this and more! Use it to move anesthetized patients to and from holding cages, or to lift reluctant patients of just about any size up to the right height for treatment. You can even take your Vet-Mate Table outside to meet critical patients as they arrive and quickly transport them inside for treatment. Eliminates the back-breaking lifting often necessary when treating large, heavy or reluctant animals.

This versatile table lowers to just 13-in. (33.02cm) above the floor to allow you to maneuver even the heaviest patients onto it with minimal lifting. Or the Vet-Mate Table can be set at cage height, car seat height, treatment height or any height you need up to 38-in. (96.52cm) high to allow you to move patients on and off with the least amount of effort. The scissor-action electric lifting mechanism is operated by a rechargeable battery to smoothly and quietly lift the table no matter where you are inside or outside of your facility. The battery recharges quickly with the built-in charger.

Heavy-gauge tubular stainless steel channel construction stands up to rugged daily use. All joints are heli-arc welded for superior strength.



The adjustable Classic VET-MATE Lift Gurney Table can be set anywhere from

13-in. (33.02cm) to 38-in. (96.52cm)

high to give you the optimal working height for your needs. Lowers to just 13-in. (33.02cm) for easier maneuvering and storage.

Foot operated chassis mounted UP/DOWN lift mechanism is located at the base of the table and can be mounted on either the End or the Side.



P/N	DESCRIPTION	VOLTS	OVERALL SIZE	LOWEST HEIGHT:	MAXIMUM HEIGHT:
12761-01-IGFVDL	Classic VET-MATE Lift Gurney	120V AC▲	53.63-in. wide x 22.5-in. deep (136.22 cmH x 57.15 cmW)	13-in. (33.02cm)	38-in. high (96.52cm)

▲ 220V requires step-down transformer. Contact SSCI.

REGAL VET MATE LIFT GURNEY

A MUST FOR EVERY VETERINARY HOSPITAL. PROTECTS PATIENTS AND YOUR STAFF BY SAFELY LIFTING AND TRANSPORTING PATIENTS THROUGHOUT THE FACILITY.

Table

Adjustable Tie-down system lets you safely position animals.

Satin-finish Table Top surface is easy on the eyes even under bright surgery lights.

One piece, die-formed Type 304 18-gauge stainless steel Flat Operating Table Top is designed for long life.

Drain trough with no-spill raised outer lip channels fluids to the collection hole.

FOUR TABLES IN ONE:

- **Lift Table**
- **Transport Gurney**
- **Exam/Treatment Table**
- **Surgery Table**

Lift Mechanism

Scissor-action lifting mechanism has smooth, quiet bronze bushings for long-lasting service.

12 v. 26 amp rechargeable battery gives you lifting capacity when and where you need it.

12 v. DC motor-driven screw lift mechanism will raise and lower animals up to 300 lbs. (136 kgs.). Scissors Housing is a one-piece, die-formed 10-ga. stainless steel.

Scissor Frame is a one-piece, fully-welded 1-1/2" x 2", 7-ga. stainless steel tubing.

Top supported by one-piece, fully-welded frame of 2" square, 7-ga. stainless steel tubing.

Scissor-action is supported by bearings riding in a milled groove in the base.

Base

Angle joints are cut, mitered and welded.

Casters recess under chassis, out of the way for safety.

Base is a one-piece, fully-welded frame of 2" square, 7-ga. stainless steel tubing.

1-pc., fully-welded, tubular stainless steel frame with a 7-gauge wall thickness has a uniform, non-glare finish.

Swivel Casters feature 5-in. dia. wheels and wheel brakes to keep table in place when needed.

Smooth-rolling, non-marring 5-in. (12.70cm) casters rotate 360° for incredible maneuverability.

It's a transport table,

a prep table, a treatment table, a lift table, a gurney, even an operating table. The Regal Vet-Mate Gurney Lift Table is all this and more!

Use it to move anesthetized patients to and from holding cages, or to lift reluctant patients of just about any size up to the right height for treatment. You can even take your Vet-Mate Table outside to meet critical patients as they arrive and quickly transport them inside for treatment. Eliminates the back-breaking lifting often necessary when treating large, heavy or reluctant animals.

This versatile table lowers to just 10.5-in. (26.67cm) above the floor to allow you to maneuver even the heaviest patients onto it with minimal lifting. Or the Vet-Mate Table can be set at cage height, car seat height, treatment height or any height you need up to 42-in. (106.68cm) high to allow you to move patients on and off with the least amount of effort. The scissor-action electric lifting mechanism is operated by a rechargeable battery to smoothly and quietly lift the table no matter where you are inside or outside of your facility. The battery recharges quickly with the built-in charger.

Heavy-gauge tubular stainless steel construction stands up to rugged daily use. All joints are heli-arc welded for superior strength.



Infinitely adjustable Vet-Mate Gurney Table can be set anywhere from 10.5-in. (26.67cm) to 42-in. (106.68cm) high to give you the optimal working height for your needs. Lowers to just 10.5-in. (26.67cm) for easier maneuvering and storage.



Lift mechanism operates by convenient switches located at both ends of the table top.

P/N	DESCRIPTION	VOLTS	OVERALL SIZE	LOWEST HEIGHT:	MAXIMUM HEIGHT:
12760-00-ITGLDK	Regal Vet-Mate Gurney Lift Table	120V AC▲	57-in. wide x 22.25-in. deep (144.78cmH x 56.51cmW)	10.5-in. (26.67cm)	42-in. high (106.68cm)
Replacement Parts					
854480	12 volt Re-chargeable Battery				

▲ 220V requires step-down transformer. Contact SSCI.