

“LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS”



SSCI

SUBURBAN SURGICAL Co., INC.

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI. Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

Catalog No. 702488-0308

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
Visit our website at: www.suburbansurgical.com


STAINLESS STEEL LIFT TABLES

SSCI Lift Tables raise and lower animals up to 300 lbs. (136.08kgs) and have the quickest lifting speeds available.

Elite Lift Tables raise working surface from 7.250 in. to 40 in. (18.41cm to 101.60cm) via an electric motor-driven screw lifting mechanism.

Elite models are offered with lateral or longitudinal exam tables fitted with or without an Regal 300 Electronic Scale.

ALL STYLES LIFT TABLES ARE DESIGNED WITH THESE QUALITY FEATURES:

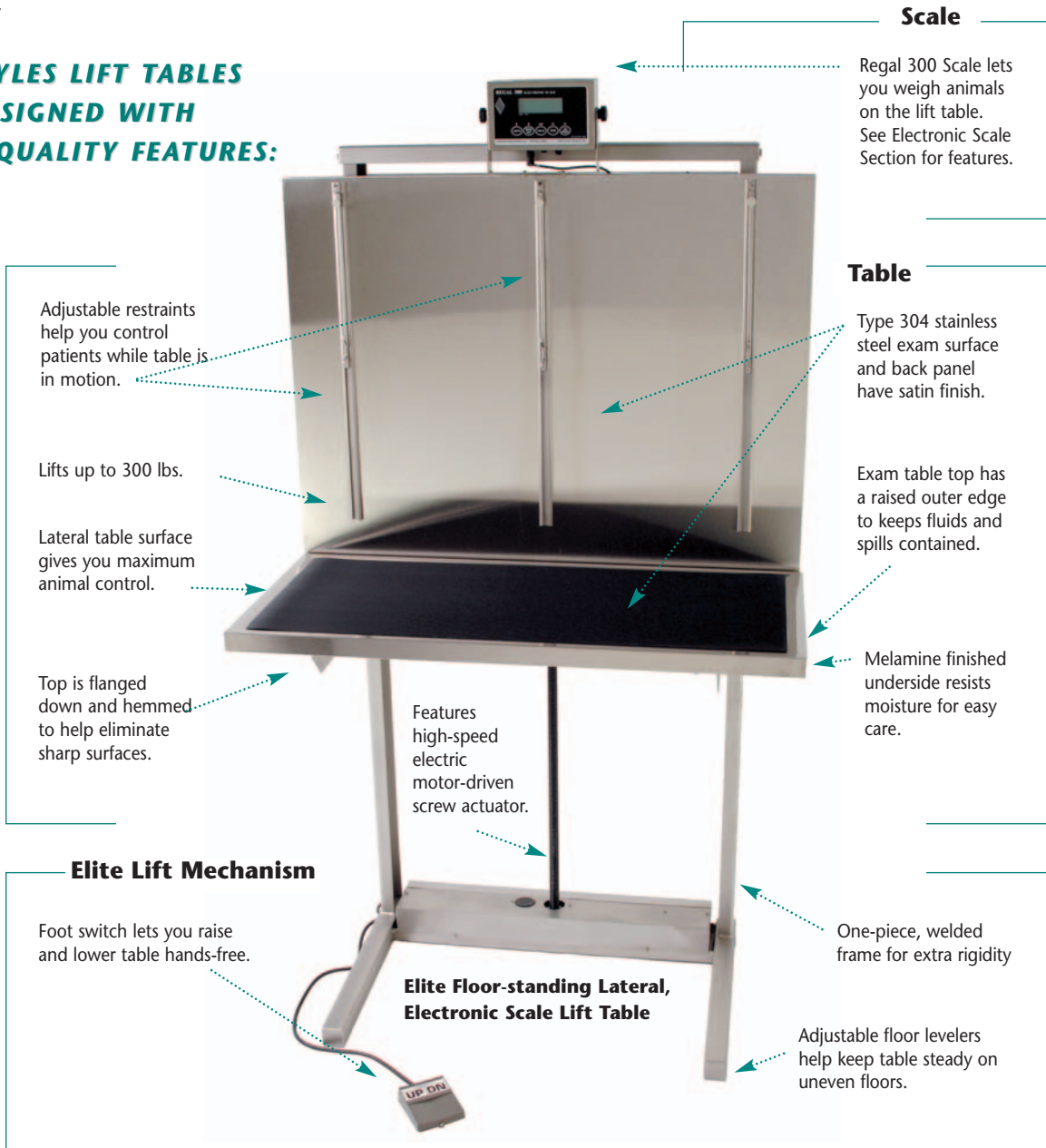

With our stainless steel Lift Tables, you'll be able to examine large patients at the height that's most comfortable for you. No more lifting large animals, and best of all, there's less chance of being harmed by a frightened animal you are trying to lift to the exam table.

All our Lift Tables feature a stainless steel exam table surface and attached back panel that slide smoothly up and down on the stainless steel frame to bring animals to an ideal working height. Foot-operated pedal controls the direction of the table, leaving your hands free to soothe and steady the animal. Back panel has adjustable restraining loops with tie-offs for greater animal control and safety while the table is in motion.

Elite Stainless Steel Lift Tables



Elite Floor-standing Lateral Lift Table



Elite Floor-standing Longitudinal Lift Table

Floor-Standing Lift Tables

Four floor levelers adjust to hold table steady on uneven floors for smooth operation and safety.

P/N	VOLTS	TOP WORKING SURFACE		OVERALL DIMENSIONS	
		INCHES	CENTIMETERS	INCHES	CENTIMETERS
Elite Floor-Standing Lateral Lift Table					
12660-02-GYLWEF	120V AC ▲, 4.6 A, 50/60 Hz	45.221Wx21.284D	114.86Wx54.06D	45.250Wx77.125Hx27D	114.94Wx195.90Hx68.58D
Elite Floor-Standing Longitudinal Lift Table					
12664-02-DPLWIB	120V AC ▲, 4.6 A, 50/60 Hz	21.284W x 45.221D	54.06Wx114.86D	23.500Wx77.125H x 51D	59.69Wx195.90Hx129.54D


▲ 220V requires step-down transformer. Contact SSCI.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com

Elite Lift Tables



Elite Floor-standing Longitudinal, Regal 300 Electronic Scale Lift Table

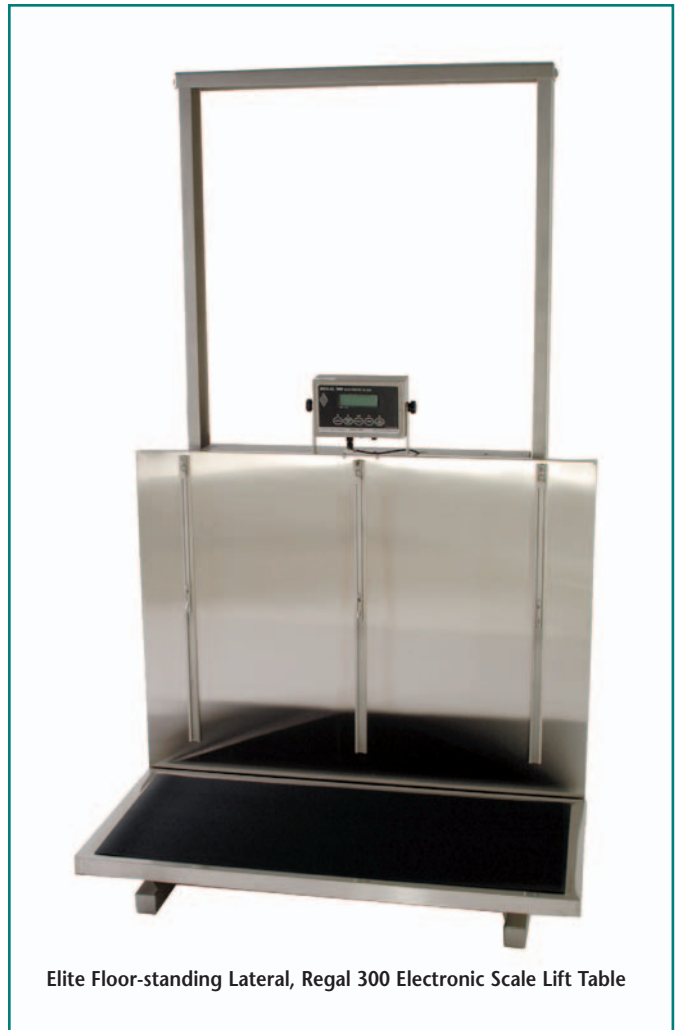

Elite Floor-standing Lateral, Regal 300 Electronic Scale Lift Table

Floor-Standing Lift Tables with Regal 300 Electronic Scale

Features an Regal 300 Electronic Platform Scale for weighing patients in pounds or kilograms up to 300 lbs. (136 kgs.). For Scale specifications; see Electronic Scale Section.

P/N	DESCRIPTION VOLTS	TOP WORKING SURFACE		OVERALL DIMENSIONS	
		INCHES	CENTIMETERS	INCHES	CENTIMETERS
Elite Floor-Standing Lateral, Regal 300 Electronic Scale Lift Table					
12661-02-GZLWEF	120V AC ▲, 4.6 A, 50/60 Hz	45.438Wx21.438D	115.33Wx54.37D	45.438Wx77.125Hx27D	115.41Wx195.90Hx68.58D
Elite Floor-Standing Longitudinal, Regal 300 Electronic Scale Lift Table					
12665-02-DPLWHX	120V AC ▲, 4.6 A, 50/60 Hz	21.438Wx45.438D	54.37Wx115.33D	23.5Wx77.125Hx51.25D	59.69Wx195.90Hx130.18D

▲ 220V requires step-down transformer. Contact SSCI.



Elite Wall-mounted Longitudinal Lift Table



Elite Wall-mounted Lateral Lift Table


Wall-Mount Lift Tables

Mounts to wall to free up floor space under the table. The wall must be structurally sound. Masonry wall or a studded wall with proper backing plates is required. Contact SSCI for construction details.

P/N	DESCRIPTION VOLTS	TOP WORKING SURFACE		OVERALL DIMENSIONS	
		INCHES	CENTIMETERS	INCHES	CENTIMETERS
Elite Lateral Lift Table					
12662-02-GZLWEF	120V AC ▲, 4.6 A, 50/60 Hz	45.221Wx21.284D	114.86Wx54.06D	45.250Wx77.125Hx27D	114.94Wx195.90Hx68.58D
Elite Longitudinal Lift Table					
12666-02-DPLWHX	120V AC ▲, 4.6 A, 50/60 Hz	21.284Wx45.221D	54.06Wx114.86D	23.500Wx77.125H x 51D	59.69Wx195.90Hx129.54D

▲ 220V requires step-down transformer. Contact SSCI.

Elite Lift Tables

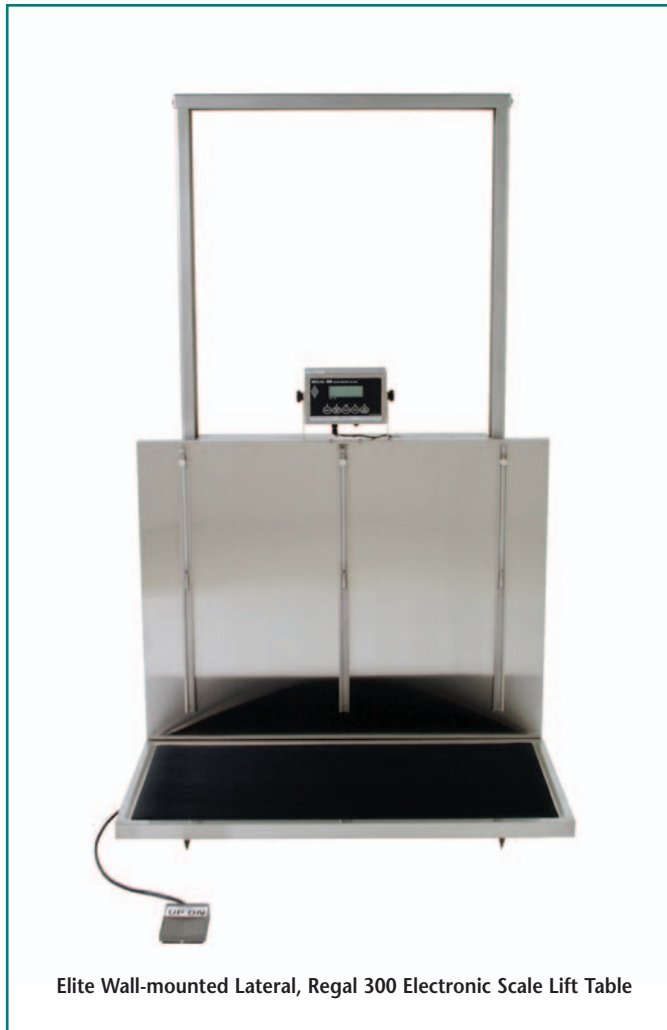

Elite Wall-mounted Lateral, Regal 300 Electronic Scale Lift Table



Elite Wall-mounted Longitudinal, Regal 300 Electronic Scale Lift Table

Wall-Mount Lift Tables with Regal 300 Electronic Scale

Features an Electronic Platform Scale for weighing patients in pounds or kilograms up to 300 lbs. (136 kgs.). For Scale specifications; see Electronic Scale Section.

P/N	DESCRIPTION VOLTS	TOP WORKING SURFACE		OVERALL DIMENSIONS	
		INCHES	CENTIMETERS	INCHES	CENTIMETERS
Elite Electronic Scale Lateral Lift Table					
12663-02-GZLWEF	120V AC ▲, 4.6 A, 50/60 Hz	45.438Wx21.438D	115.33Wx54.37D	45.438Wx77.125Hx27D	115.41Wx195.90Hx68.58D
Elite Electronic Scale Longitudinal Lift Table					
12667-02-DPLWHX	120V AC ▲, 4.6 A, 50/60 Hz	21.438Wx45.438D	54.37Wx115.33D	23.5Wx77.125Hx51.25D	59.69Wx195.90Hx130.18D

▲ 220V requires step-down transformer. Contact SSCI.



Swivel Casters for Lift Tables


Add these Swivel Casters to any of our floor-standing Lift Tables to create a mobile lift table. Swivel Casters with 3-in. (7.62cm) diameter wheel with wheel brake, allow for easy movement of lift table.

P/N	DESCRIPTION
C010-23061001-CH000	Nickel-plated 5-in. (12.70cm) diam. wheel w/Brake

Black Vinyl Ribbed Mat

Rubber vinyl gives animals surer footing and greater security when table is raised and lowered.

P/N	FITS TABLE STYLE
12800-00-GNAACV	Fits Elite Lift Tables

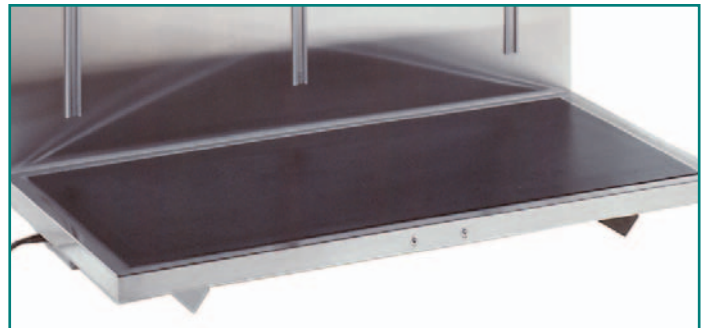

Add-on Drawer for Lift Tables

Stainless steel drawer mounts easily underneath Lift Table to keep supplies close at hand. May restrict lowest position of table's top.

P/N	DESCRIPTION	INCHES	CENTIMETERS
12720-00-CLAPBR	Add-on Stainless Steel Drawer Kit	10.875Wx4Hx15.936D	(27.62Wx10.16Hx40.48D)

▲ 220V requires step-down transformer. Contact SSCI.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com

