

“LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS”



SSCI

SUBURBAN SURGICAL Co., INC.

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI. Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

Catalog No. 702488-0308

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
Visit our website at: www.suburbansurgical.com


SURGEON SCRUB SINKS

SURGEON SCRUB SINKS HAVE THESE QUALITY FEATURES:

Our 16-gauge, Type 304 Stainless Steel Surgeon Scrub Sinks are fully heli-arc welded into one piece to create the strongest Surgeon Scrub Sink available. Welds are continuous, crevice-free and ground smooth. Coved corners eliminate areas where dirt and debris accumulate. Radius bends on bottom and all sides and a uniform, non-glare finish results in a tub interior

that is easy to clean. Scrub Sinks feature a 10.50-in. (26.67cm) deep basin, a 12-in. (30.48cm) high back splash. The Surgeon Scrub Sink is designed to accept a back splash-mounted gooseneck faucet and 3.5-in. (8.89cm) drain. Stainless Steel Wall-Mounted Brackets are included.

One-Piece Welded Surgeon scrub Sink

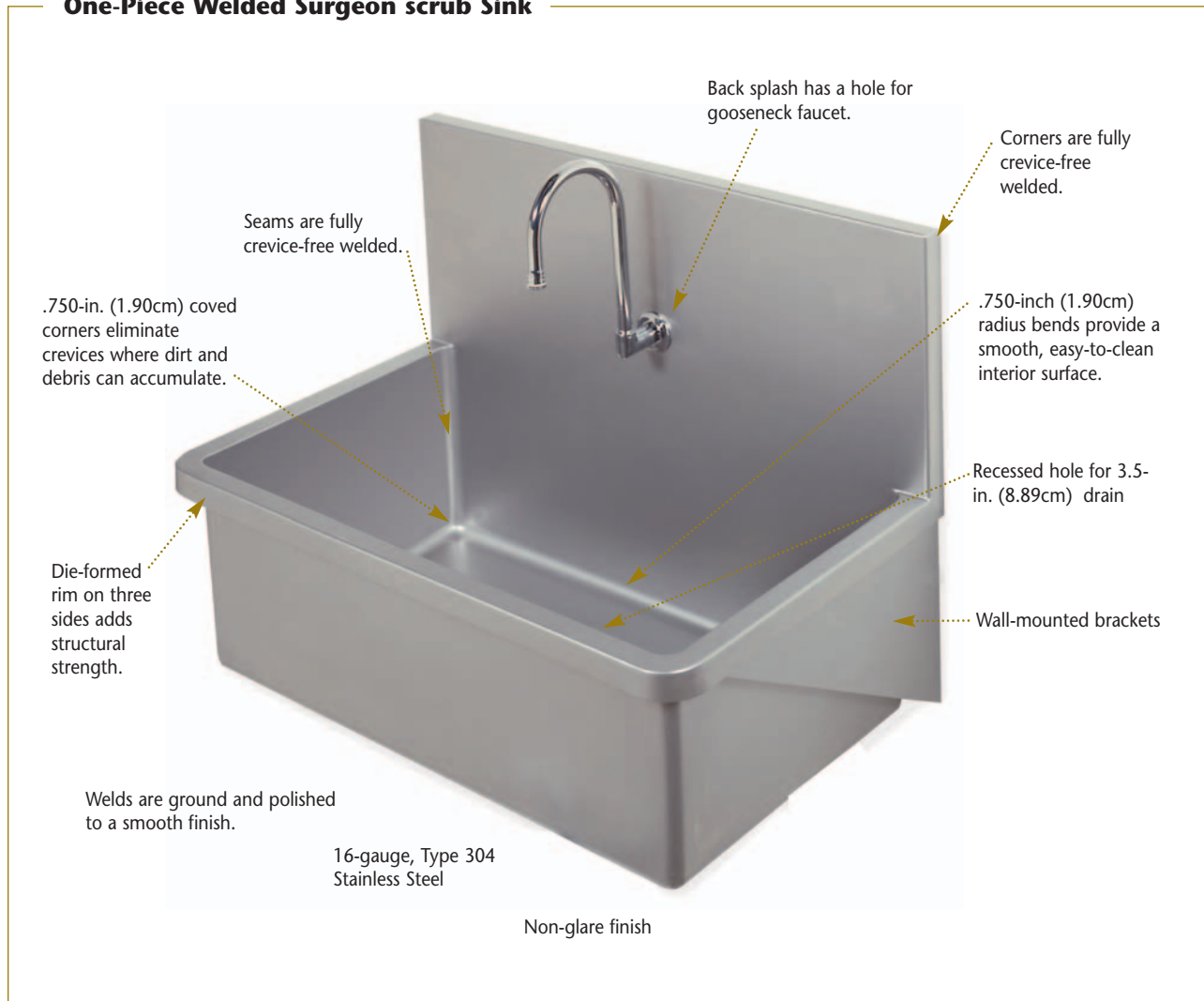

A clean surgical scrub area is an essential part of sanitary, bacteria-free procedures. We make it easy for you to scrub up with these top-quality all stainless steel Surgeon Scrub Sinks.

Both sizes of Surgeon Scrub Sinks are designed to accept a gooseneck faucet. Faucets can be controlled by three types of hands-free valves; foot-operated, knee-operated and electronic-motion sensor (see Section 18).



Single-Station



Dual-Station

Surgeon Scrub Sinks

P/N	Description	OVERALL SIZE	
		Inches	Centimeters
13320-00-EPDMCX	Single-Station	30W x 22.56H x 19D	76.20W x 57.30H x 48.26D
13321-00-HJDMCX	Dual-Station	48W x 22.56H x 19D	121.92W x 57.30H x 48.26D

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com