

“LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS”



SSCI

SUBURBAN SURGICAL Co., Inc.

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
Visit our website at: www.suburbansurgical.com

VETERINARY LIGHTING

A LIGHT FOR EVERY VETERINARY FUNCTION

Selecting the right lighting is a very important element to consider when outfitting your veterinary practice. You need proper illumination without generating a lot of heat or creating color distortion.

At Suburban Surgical, we've taken great care to find the very best medical lighting available. We have surgical lights, treatment/procedure lights and exam diagnostic lights to fit your needs and budget. And most of our lights are available in a wall-

mount, ceiling-mount or portable floor model. Many styles may also be track-mounted for the utmost positioning flexibility.

In addition, we offer the Medi-Vue Magnifying Light as well as a "Woods" Diagnostic Magnifying Light in wall-mount, table-mount and portable floor models to enhance your view of the treatment area.



ExceLED
The New Standard
in Surgery Lights



Table Mount Medi-Vue
Magnifying Light




Single or Dual Trolley ChuttleTrak®
Each trolley accommodates up to 2 lights for a maximum of 4 lights per track.



Ceiling Mount CENTURION
Excel Series Surgery Light



Century "Woods"
Diagnostic Light

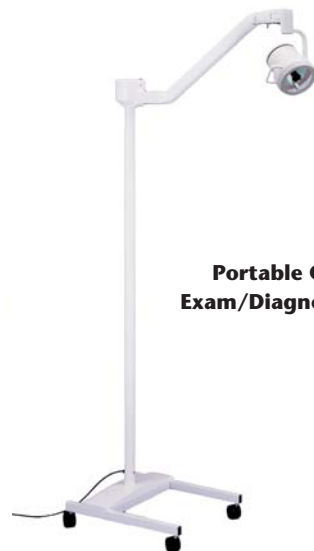

Portable CENTRY
Exam/Diagnostic Light



Wall-Mount CENTURA
Full Variable Spot Light



SYSTEM TWO Metal Halide
TRIO Ceiling Mounted Light




EXCELED A NEW ERA UNFOLDS WITH GREEN TECHNOLOGY

ExceLED represents the new standard in surgery lighting with its green technology LED lights. It features a larger reflector with high efficiency at low heat, a larger pattern and long life LED bulbs.


- Very Low Heat
- High Efficiency — Half the Electric Consumption of Standard Halogen Lights
- Long Life L.E.D.s
- 4300° Kelvin Color Temperature
- 65,000 Lux (6,000 foot candles)
- Available in a Variety of Mounting Options
- Patented 5-Stage Dimming and On/Off Switch, Controlled by the Removable Sterilizable Handle
- Additional Positioning Handle on Perimeter of Reflector
- EXCELED Minor Surgery Lights have auto seeking power supply (100V-240V)



Portable

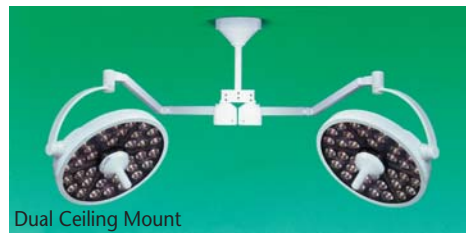


Dual Ceiling Mount

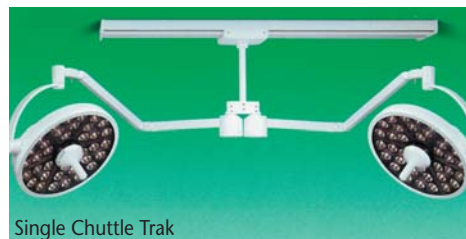

Single Chuttle Trak



Ceiling Mounted

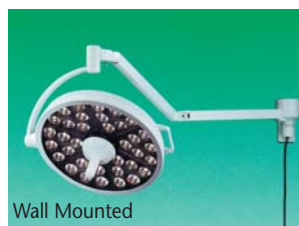

Wall Mounted



Dual Chuttle Trak

P/N	DESCRIPTION
MI000000-061512	EXCELED Wall Mounted (100V-240V)
MI000000-061513	EXCELED Portable Floor with Casters (100V-240V)
MI000000-061514	EXCELED Ceiling Mounted One (100V-240V)
MI000000-061515	EXCELED Ceiling Mounted Two (100V-240V)
<i>Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering.</i>	
CEILING ROD KITS	
Standard Ceiling Rod kit includes: Casted Mounting Ring, Rod and Wiring Harness. Mounted Ring must be attached to structural ceiling mount.	
MI00000-1001186	16" Ceiling Rod Kit for structural ceiling mount 7'6" - 8'6" above the floor.
MI00000-1001187	25" Ceiling Rod Kit for structural ceiling mount 8'7" - 9'6" above the floor.
MI00000-1001188	35" Ceiling Rod Kit for structural ceiling mount 9'7" - 10'6" above the floor.
Extended Ceiling Rod kit includes: H-D Casted Mounting Ring, Rod and Wiring Harness. Mounted Ring must be attached to structural ceiling mount.	
MI00-1000976-44	Extended Ceiling Rod Kit for structural ceiling mount 10'2" - 10'10" above the floor.
MI00-1000976-53	Extended Ceiling Rod Kit for structural ceiling mount 10'11" - 11'8" above the floor.
MI00-1000976-63	Extended Ceiling Rod Kit for structural ceiling mount 11'9" - 12'0" above the floor.
MI00-1000976-67	Extended Ceiling Rod Kit for structural ceiling mount 12'1" - 12'10" above the floor.
MI00-1000976-77	Extended Ceiling Rod Kit for structural ceiling mount 12'11" - 14'0" above the floor.
CEILING TRACKS	
ChuttleTrak™ offers exceptional versatility since it allows the horizontal movement of the lights, in addition to the movement that the individual lights offer.	
MI000000-044011	CHUTTLE TRAK Single Trolley 70-in. (178 cm) holds one Ceiling Mount or one Dual Ceiling Mounts *
MI000000-044012	CHUTTLE TRAK Dual Trolley 70-in. (178 cm) holds two Single Ceiling Mounts or a maximum of two Dual Ceiling Mounts *
* Ceiling Rod kit is required -- only 16" or 25" Rods can be used.	
ACCESSORIES	
MI00000-1000616	Extra Sterilizable Handle

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com

SYSTEM TWO Halogen

SYSTEM TWO HALOGEN Series

SYSTEM TWO offers great performance and high value, with no true competitor in the market. SYSTEM TWO lights are available in 16" or 20" diameter multifaceted reflectors. The SYSTEM TWO is available in a SOLO, DUO or TRIO Articulating Arm Assembly, which features continuous rotation of all key pivot points.

- Halogen Light Source
- 4000° Kelvin Color Temperature
- S1 Satellite Model: 16" diameter Multifaceted Reflector
72,000 Lux (6,700 foot candle output at 1 meter)
- S2 Surgery Light Model: 20" diameter Multifaceted Reflector
102,000 Lux (9,500 foot candle output at 1 meter)
- Available in a Variety of Mounting Options
- 4-Stage Dimming and On/Off Switch, Controlled by the Removable Sterilizable Handle
- Additional Positioning Handle on Perimeter of Reflector
- Video Integration Capabilities



Video Camera



Arm for Monitor
(monitor not included)



S1 Satellite

S2 Surgery Light



SOLO Ceiling Mounted

SOLO Ceiling Mounted is the Ideal Surgery Light.

P/N		DESCRIPTION
115V	220V	
MI00000-SS1120	MI00000-SS1230	One S1 Satellite
MI00000-SS2120	MI00000-SS2230	One S2 Surgery Light

Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 7 for Ceiling Rod Kits.

DUO Ceiling Mounted Two articulating arms rotates fully independent of the other, allowing limitless positioning of each arm assembly.

P/N		DESCRIPTION
115V	220V	
MI0000-SS1S1120	MI0000-SS1S1230	Two S1 Satellite
MI0000-SS2S2120	MI0000-SS2S2230	Two S2 Surgery Light
MI0000-SS2S1120	MI0000-SS2S1230	One S2 Surgery Light & One S1 Satellite
MI0000-SS2M1120	MI0000-SS2M1230	One S2 Surgery Light & One Monitor Arm
MI0000-SS1M1120	MI0000-SS1M1230	One S1 Satellite & One Monitor Arm
MI0000-SS2P1120	MI0000-SS2P1230	One S2 Surgery Light & One Video Camera
MI0000-SS1P1120	MI0000-SS1P1230	One S1 Satellite & One Video Camera

Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 7 for Ceiling Rod Kits.



DUO Ceiling Mounted



TRIO Ceiling Mounted Three articulating arms provide surgery lighting with full video integration for the state-of-the-art surgery room.



TRIO Ceiling Mounted

P/N		DESCRIPTION
115V	220V	
MI00-SS2S2S2120	MI00-SS2S2S2230	Three S2 Surgery Lights
MI00-SS2S2S1120	MI00-SS2S2S1230	Two S2 Surgery Lights & One S1 Satellite
MI00-SS2S1S1120	MI00-SS2S1S1230	One S2 Surgery Light & Two S1 Satellites
MI00-SS1S1S1120	MI00-SS1S1S1230	Three S1 Satellites
MI00-SS2S2M1120	MI00-SS2S2M1230	Two S2 Surgery Light & One Monitor Arm
MI00-SS2S1M1120	MI00-SS2S1M1230	One S2 Surgery Light, One S1 Satellite & One Monitor Arm
MI00-SS1S1M1120	MI00-SS1S1M1230	Two S1 Satellites & One Monitor Arm
MI00-SS1S1P1120	MI00-SS1S1P1230	Two S1 Satellites & One Video Camera
MI00-SS2S2P1120	MI00-SS2S2P1230	Two S2 Surgery Lights & One Video Camera
MI00-SS2S1P1120	MI00-SS2S1P1230	One S2 Surgery Light, One S1 Satellite & One Video Camera


Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 7 for Ceiling Rod Kits.



Portable



Wall Mounted

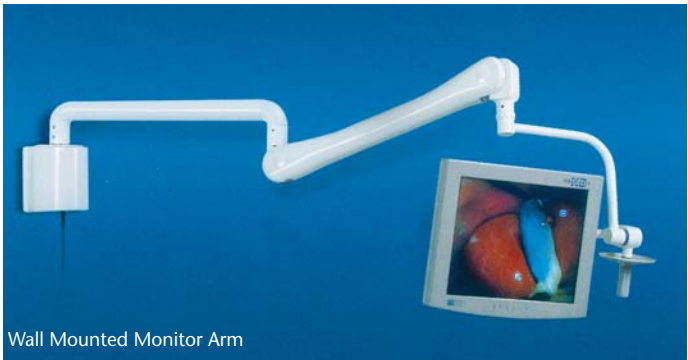

Wall Mounted Monitor Arm

P/N		DESCRIPTION
PORTABLE		
115V	220V	
MI000000-FS2120	MI000000-FS2230	Portable Floor with Casters One S2 Surgery Light
MI000000-FS1120	MI000000-FS1230	Portable Floor with Casters One S1 Satellite
WALL MOUNTED		
115V	220V	
MI000000-WS2120	MI000000-WS2230	One S2 Surgery Light
MI000000-WS1120	MI000000-WS1230	One S1 Satellite
MI000000-WM1120	MI000000-WM1230	One Monitor Arm, with S Video Cable and a 120V Cable (light or monitor are not included)

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com

SYSTEM TWO Metal Halide

SYSTEM TWO MH Series

SYSTEM TWO METAL HALIDE Series is the next level in the evolution of surgery lighting. This all new design heralds a notable increase in efficiency over traditional halogen light sources, producing over 100,000 lux output with only a 39 watt H.I.D. bulb, versus a 150 watt halogen bulb. Heat has been virtually eliminated, white color temperature has been significantly increased to over 4,300 Kelvin. Additionally, the Metal Halide light boasts a rated bulb life of over 10,000 hours and includes complete mechanical dimming, the most accurate method.

- Most efficient surgery lighting source on the market
- Extremely cool output -- lowest possible heat
- 4300° Kelvin Color Temperature
- 6.5" - 9" adjustable pattern
100,000 Lux (output at 1 meter)
- Over 10,000 hour rated bulb life for primary bulbs
- Available in a Variety of Mounting Options
- Mechanical On/Off Switch, wide range (50-100%) dimming
- Additional Positioning Handle on Perimeter of Reflector
- Video Integration Capabilities



Video Camera



Arm for Monitor
(monitor not included)



S1 Satellite

S2 Surgery Light



SOLO Ceiling Mounted

SOLO Ceiling Mounted is the Ideal Surgery Light.

P/N	DESCRIPTION
115V MI0000-SS2120MH	220V MI0000-SS2230MH
One S2 Surgery Light	
<i>Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 7 for Ceiling Rod Kits.</i>	

DUO Ceiling Mounted Two articulating arms rotates fully independent of the other, allowing limitless positioning of each arm assembly.

P/N	DESCRIPTION
115V MI00-SS2S2120MH	220V MI00-SS2S2230MH
Two S2 Surgery Lights	
MI00-SS2M1120MH	MI00-SS2M1230MH
One S2 Surgery Light & One Monitor Arm	
MI00-SS2P1120MH	MI00-SS2P1230MH
One S2 Surgery Light & One Video Camera	
<i>Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 7 for Ceiling Rod Kits.</i>	



DUO Ceiling Mounted



SYSTEM TWO Metal Halide

TRIO Ceiling Mounted Three articulating arms provide surgery lighting with full video integration for the state-of-the-art surgery room.



TRIO Ceiling Mounted

P/N		DESCRIPTION
115V	220V	
MI-SS2S2S2120MH	MI-SS2S2S2230MH	Three S2 Surgery Lights
MI-SS2S2M1120MH	MI-SS2S2M1230MH	Two S2 Surgery Lights & One Monitor Arm
MI-SS2S2P1120MH	MI-SS2S2P1230MH	Two S2 Surgery Lights & One Video Camera


Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 7 for Ceiling Rod Kits.



Portable



Wall Mounted

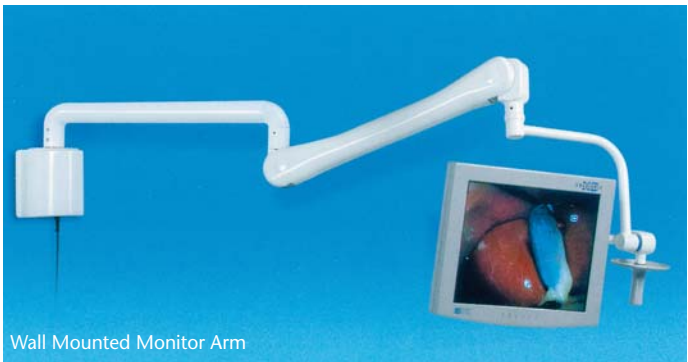

Wall Mounted Monitor Arm

P/N	DESCRIPTION	
PORTABLE		
115V	220V	
MI0000-FS2120MH	MI0000-FS2230MH	Portable Floor with Casters One S2 Surgery Light
WALL MOUNTED		
115V	220V	
MI0000-WS2120MH	MI0000-WS2230MH	One S2 Surgery Light
MI000000-WM1120	MI000000-WM1230	One Monitor Arm, with 5 Video Cable and a 120V Cable (light or monitor are not included)

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com

Accessories for SYSTEM TWO

P/N	DESCRIPTION
CEILING ROD KITS	
SOLO	
NA	Ceiling Rod Kit for structural ceiling mount 96" above the floor. (Headroom to the bottom of lights extension arm is approx 77.75".)
MI000-1000360-1	9.4" Ceiling Rod Kit for structural ceiling mount 108" above the floor. (Headroom to the bottom of lights extension arm is approx 80.35".)
MI000-1000369-1	13" Ceiling Rod Kit for structural ceiling mount 108" above the floor. (Headroom to the bottom of lights extension arm is approx 76.75".)
MI000-1000361-1	19" Ceiling Rod Kit for structural ceiling mount 120" above the floor. (Headroom to the bottom of lights extension arm is approx 82.75".)
MI000-1000388-1	25" Ceiling Rod Kit for structural ceiling mount 120" above the floor. (Headroom to the bottom of lights extension arm is approx 76.75".)
MI000-1000362-1	31" Ceiling Rod Kit for structural ceiling mount 132" above the floor. (Headroom to the bottom of lights extension arm is approx 82.75".)
MI000-1000363-1	43" Ceiling Rod Kit for structural ceiling mount 144" above the floor. (Headroom to the bottom of lights extension arm is approx 82.75".)
DUO	
NA	Ceiling Rod Kit for structural ceiling mount 96" above the floor. (Headroom to the bottom of lights extension arm is approx 72.75".)
MI000-1000360-1	9.4" Ceiling Rod Kit for structural ceiling mount 108" above the floor. (Headroom to the bottom of lights extension arm is approx 75.35".)
MI000-1000369-1	13" Ceiling Rod Kit for structural ceiling mount 108" above the floor. (Headroom to the bottom of lights extension arm is approx 71.75".)
MI000-1000361-1	19" Ceiling Rod Kit for structural ceiling mount 120" above the floor. (Headroom to the bottom of lights extension arm is approx 77.75".)
MI000-1000388-1	25" Ceiling Rod Kit for structural ceiling mount 120" above the floor. (Headroom to the bottom of lights extension arm is approx 71.75".)
MI000-1000362-1	31" Ceiling Rod Kit for structural ceiling mount 132" above the floor. (Headroom to the bottom of lights extension arm is approx 77.75".)
MI000-1000363-1	43" Ceiling Rod Kit for structural ceiling mount 144" above the floor. (Headroom to the bottom of lights extension arm is approx 77.75".)
TRIO	
NA	Ceiling Rod Kit for structural ceiling mount 96" above the floor. (Headroom to the bottom of lights extension arm is approx 67.75".)
MI000-1000360-1	9.4" Ceiling Rod Kit for structural ceiling mount 108" above the floor. (Headroom to the bottom of lights extension arm is approx 70.35".)
MI000-1000369-1	13" Ceiling Rod Kit for structural ceiling mount 108" above the floor. (Headroom to the bottom of lights extension arm is approx 66.75".)
MI000-1000361-1	19" Ceiling Rod Kit for structural ceiling mount 120" above the floor. (Headroom to the bottom of lights extension arm is approx 72.75".)
MI000-1000388-1	25" Ceiling Rod Kit for structural ceiling mount 120" above the floor. (Headroom to the bottom of lights extension arm is approx 66.75".)
MI000-1000362-1	31" Ceiling Rod Kit for structural ceiling mount 132" above the floor. (Headroom to the bottom of lights extension arm is approx 72.75".)
MI000-1000363-1	43" Ceiling Rod Kit for structural ceiling mount 144" above the floor. (Headroom to the bottom of lights extension arm is approx 72.75".)



Disposable Handle Cover

P/N	DESCRIPTION
ACCESSORIES	
MI00000-0001708	Disposable Handle Covers Carton of 40
MI00000-0001709	Disposable Handle Covers Carton of 160

P/N	DESCRIPTION
ACCESSORIES	
MI00000-0001282	Extra Sterilizable Handle for Halogen Lights
MI00000-1000563	Extra Sterilizable Handle for Metal Halide Lights





Sterilizable Handle

P/N	DESCRIPTION
BULBS	
MI00000-0003042	S2 Halogen Surgery Light Bulb (24V/150 watts)
MI00000-0003224	S1 Halogen Surgery Satellite Bulb (24V/150 watts)
MI00000-0003473	Metal Halide Surgery Light Backup Bulb (Halogen 110V/39 watts)
MI00000-0003479	Metal Halide Surgery Light Primary Bulb (Metal Halide 110V/39 watts)



Centurion Standard is a high-quality Minor Surgery light offers features and performance of much more expensive lights. Take the highly efficient single reflector design that assures an optimal, shadow-free surgical environment. All have fully adjustable articulated extension arms allow you to position the light head exactly where you need it. The pre-focused lens, along with an elliptical reflector and an infrared absorbing cylinder, creates a bright, 6-in. (150mm) diameter field of illumination that's cool, color-corrected and shadow-free.

- Halogen Light Source
- Illumination that's cool
- 4000° Kelvin Color Temperature
- 15" Diameter Single Reflector
43,000 Lux (4,000 foot candle output at 36-in/914 mm)
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- Lamphead and arm both rotate 540°
- Removable sterilizable positioning handle


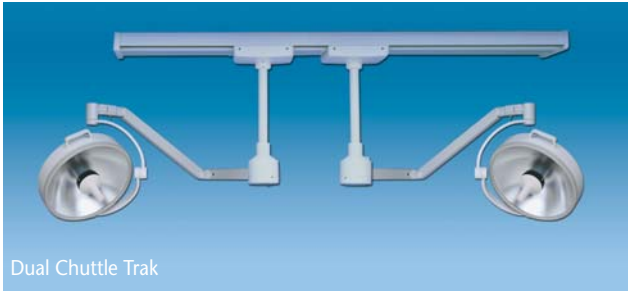
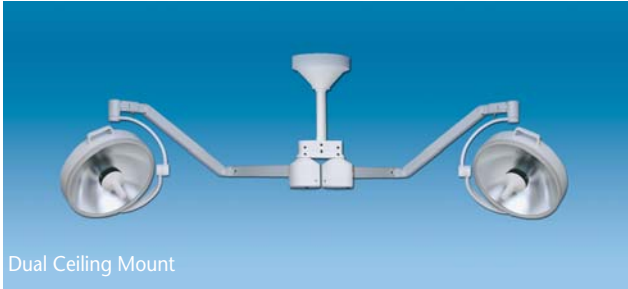


P/N		DESCRIPTION
115V	220V	
MI000000-021512	MI000000-021542	Wall Mounted
MI000000-021513	MI000000-021543	Portable Floor with Casters
MI000000-021514	MI000000-021544	Ceiling Mounted One
MI000000-021515	MI000000-021545	Ceiling Mounted Two
<i>Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 13 for Ceiling Rod Kits and Chuttle Trak.</i>		
BULBS		
MI000000-0001256	Tungsten Halogen Bulb (12V/90 watts)+B152	

21ST CENTURY

Centurion Excel has unparalleled positioning flexibility combined with high intensity illumination and quality manufacturing has made this surgery light a legend in the medical field. Pre-focused lens, plus an multi-faceted reflector and an IR absorbing cylinder, creates a bright, 6-in. (150mm) diameter field of illumination that's cool, color-corrected and shadow-free. The illumination is equal to 6,000 ft. candle intensity (64,500 lux) at a distance of 39-in. (1 meter) with a color temperature of 4,000° Kelvin to give you an accurate view of tissues in the operating field.

- Halogen Light Source
- Illumination that's cool
- 4000° Kelvin Color Temperature
- 15" Diameter Multi-Faceted Reflector
64,500 Lux (6,000 foot candle output at 39-in/1-m)
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- Lamphead and arm both rotate 540°
- Removable sterilizable positioning handle

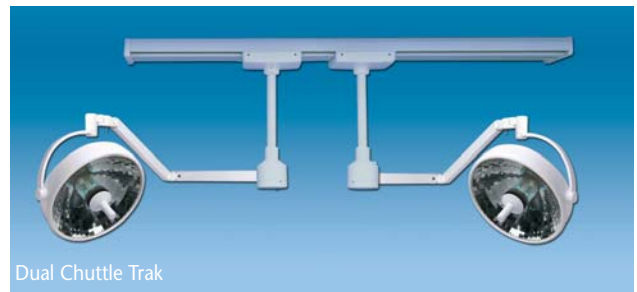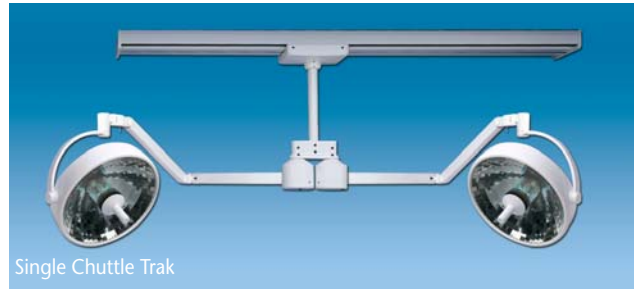
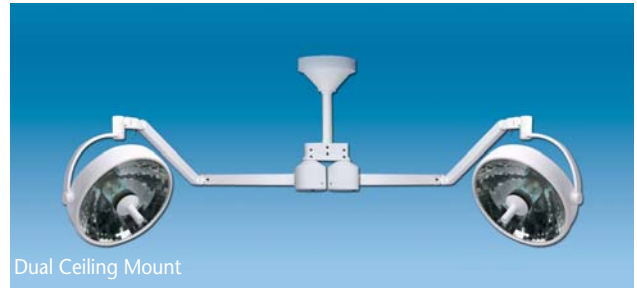


P/N		DESCRIPTION
115V	220V	
MI000000-051512	MI000000-051542	Wall Mounted One 15" Multifaceted Reflector Halogen Minor Surgery Light
MI000000-051512B	MI000000-051542B	Wall Mounted One 15" Multifaceted Reflector Halogen Minor Surgery Light w/Battery Backup
MI000000-051513	MI000000-051543	Portable Floor with Casters One 15" Multifaceted Reflector Halogen Minor Surgery Light
MI000000-051513B	MI000000-051543B	Portable Floor with Casters One 15" Multifaceted Reflector Halogen Minor Surgery Light w/Battery Backup
MI000000-051514	MI000000-051544	Ceiling Mounted One 15" Multifaceted Reflector Halogen Minor Surgery Light
MI000000-051515	MI000000-051545	Ceiling Mounted Two 15" Multifaceted Reflector Halogen Minor Surgery Light
MI000000-051514B	MI000000-051544B	Ceiling Mounted One 15" Multifaceted Reflector Halogen Minor Surgery Light w/Battery Backup
<i>Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering.</i>		
<i>See page 13 for Ceiling Rod Kits and Chuttle Trak.</i>		
BULBS		
MI000000-0001256	Tungsten Halogen Bulb (12V/90 watts)+B152	



Centura is a highly versatile full-featured Spot Light features a long-life halogen bulb and advanced optics to provide a tightly focused, high-intensity spot of light to meet a variety of medical needs. Add the optional built-in fiberoptic light generator with a 4-position turret to power fiberoptic headlamps or instruments for even greater versatility. Spot adjusts anywhere from 3-in. (76mm) to 16-in. (406mm) in diameter providing 2,000 ft. candles (21,500 lux) of illumination at a distance of 36-in. (914mm).

- Halogen Light Source
- Optional Built-in Fiber Optics
- 3500° Kelvin Color Temperature
- 2" Diameter Dichronic Glass Reflector
21,500 Lux (2000 foot candle output at 36-in/914-mm)
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- Lamphead and arm both rotate 540°
- Removable sterilizable positioning handle

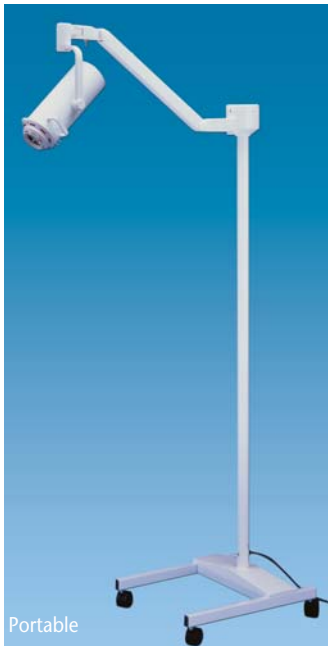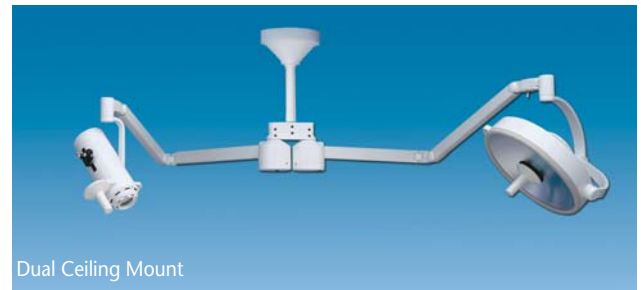

P/N	DESCRIPTION
115V	220V
MI000000-022812	MI000000-022842 Wall Mounted One Halogen Spot Light
MI000000-022813	MI000000-022843 Portable Floor with Casters One Halogen Spot Light
MI000000-022814	MI000000-022844 Ceiling Mounted One Halogen Spot Light *
MI000000-022815	MI000000-022845 Ceiling Mounted Two Halogen Spot Lights *
MI000000-022912	MI000000-022942 Wall Mounted One Halogen Spot Light w/Fiber Optic Generator
MI000000-022913	MI000000-022943 Portable Floor with Casters One Halogen Spot Light w/Fiber Optic Generator
MI000000-022914	MI000000-022944 Ceiling Mounted One Halogen Spot Light w/Fiber Optic Generator *
MI000000-022915	MI000000-022945 Ceiling Mounted Two Halogen Spot Lights w/Fiber Optic Generator *

*Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering.
See page 13 for Ceiling Rod Kits and Chuttle Trak.*



BULBS

MI000000-0001254	Fiber Optic Rear Bulb (24V/150 watts)
MI000000-0001255	Front Tungsten Halogen Bulb (21V/150 watts)

21ST CENTURY

Centry is an exceptionally compact light has a long-life halogen bulb to provide a large flood pattern of light with an intensity of 4,000 ft. candles (43,000 lux) and a color temperature of 4,300° Kelvin. The 5-in. (127mm) glass dichroic reflector creates a brilliant, prefocused 6-in. (150mm) spot at a distance of 24-in. (610mm) from lamphead. The **Centry** Light has two handle styles, fixed or removable.

- Halogen Light Source
- 4300° Kelvin Color Temperature
- Lens Filter System
- 5" Diameter Glass Dichroic Reflector
43,000 Lux (4000 foot candle output at 24-in/610-mm)
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- Lamphead and arm both rotate 540°
- Handle mounted to the perimeter of light head

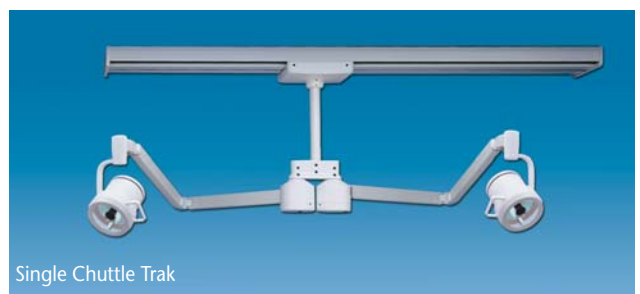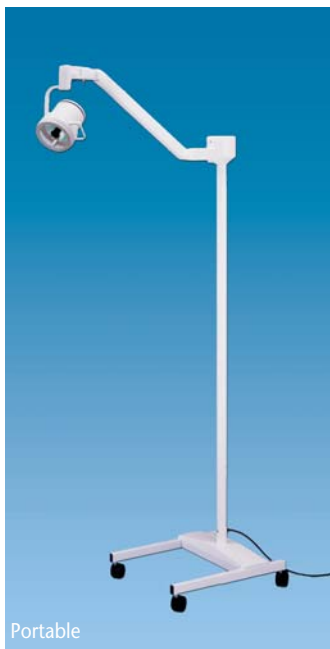

P/N		DESCRIPTION
115V	220V	
MI000000-021312	MI000000-021342	Wall Mounted with Fixed Handle
MI00-021312-SHO	MI00-021342-SHO	Wall Mounted with Sterilizable Handle
MI000000-021313	MI000000-021343	Floor Standing with Fixed Handle
MI00-021313-SHO	MI00-021343-SHO	Floor Standing with Sterilizable Handle
MI000000-021313D	MI000000-021343D	O.B. Portable Floor, Casters with Fixed Handle
MI0-021313D-SHO	MI0-021343D-SHO	O.B. Portable Floor, Casters with Sterilizable Handle
MI000000-021314	MI000000-021344	Ceiling Mounted with Fixed Handle
MI00-021314-SHO	MI00-021344-SHO	Ceiling Mounted with Sterilizable Handle
MI000000-021315	MI000000-021345	Ceiling Mounted Two with Fixed Handle
MI00-021315-SHO	MI00-021345-SHO	Ceiling Mounted Two with Sterilizable Handle
<i>Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 13 for Ceiling Rod Kits and Chuttle Trak.</i>		
BULBS		
MI00000-0001257	Tungsten Halogen Bulb (12V/75 watts)	
MI00000-0001717	O.B. Tungsten Halogen Bulb (12V/75 watts)	



21ST CENTURY COMBINATIONS


Combination Lights have been designed to provide you ease of ordering for our most popular combinations. All of these combinations can be used on the Single or Dual ChuttleTrak. See page 13 for more information on Ceiling Rod Kits and ChuttleTrak Ceiling Track.



Centurion Excel and Centura



Centurion Standard and Centura

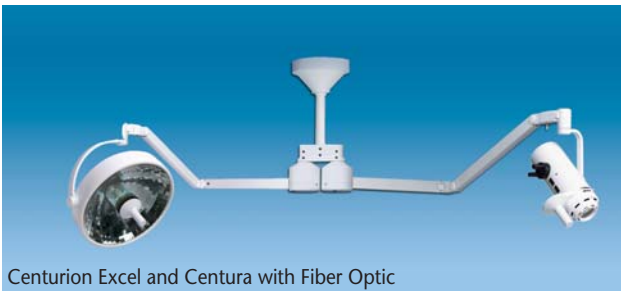


Centurion Excel and Centura with Fiber Optic



Centurion Standard and Centura with Fiber Optic



Centura and Centura with Fiber Optic

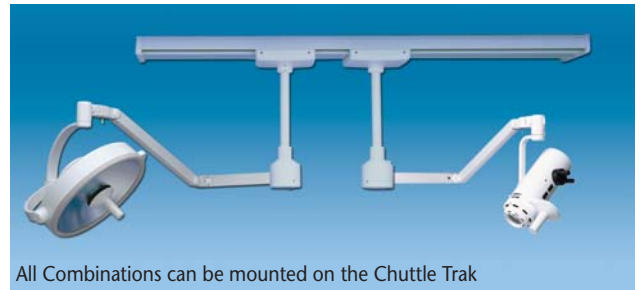

All Combinations can be mounted on the Chuttle Trak

P/N		DESCRIPTION
115V	220V	
MI000000-058515	MI000000-058545	CENTURION EXCEL Ceiling Mounted One 15" Multifaceted Reflector Halogen Minor Surgery Light & CENTURA Ceiling Mounted One Halogen Spot Light
MI000000-059515	MI000000-059545	CENTURION EXCEL Ceiling Mounted One 15" Multifaceted Reflector Halogen Minor Surgery Light & CENTURA Ceiling Mounted One Halogen Spot Light w/Fiber Optic Generator
MI000000-023015	MI000000-023045	CENTURA Ceiling Mounted One Halogen Spot Light & CENTURA One Halogen Spot Light w/Fiber Optic Generator
MI000000-028515	MI000000-028545	CENTURION STANDARD Ceiling Mounted One 15" Reflector Halogen Minor Surgery Light & CENTURA One Halogen Spot Light
MI000000-029515	MI000000-029545	CENTURION STANDARD Ceiling Mounted One 15" Reflector Halogen Minor Surgery Light & CENTURA One Halogen Spot Light w/Single Fiber Optic Light Generator

Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering. See page 13 for Ceiling Rod Kits and Chuttle Trak.

ACCESSORIES

MI000000-0001282	Extra Sterilizable Handle
MI000000-0001708	Disposable Sterile Handle Covers - Carton of 40
MI000000-0001709	Disposable Sterile Handle Covers - Carton of 160

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com



21ST CENTURY CHUTTLE TRAK & CEILING ROD KITS

Ceiling Rod Kits can be used on all models of the 21ST CENTURY Medical Lights. Correct height from the floor to the structural ceiling mount must be known at time of ordering.

P/N	DESCRIPTION
CEILING ROD KITS	
Standard Ceiling Rod kit includes: Casted Mounting Ring, Rod and Wiring Harness. Mounted Ring must be attached to structural ceiling mount.	
MI00000-1001186	16" Ceiling Rod Kit for structural ceiling mount 7'6" - 8'6" above the floor.
MI00000-1001187	25" Ceiling Rod Kit for structural ceiling mount 8'7" - 9'6" above the floor.
MI00000-1001188	35" Ceiling Rod Kit for structural ceiling mount 9'7" - 10'6" above the floor.
Extended Ceiling Rod kit includes: H-D Casted Mounting Ring, Rod and Wiring Harness. Mounted Ring must be attached to structural ceiling mount.	
MI00-1000976-49	Extended Ceiling Rod Kit for structural ceiling mount 10'2" - 10'10" above the floor.
MI00-1000976-53	Extended Ceiling Rod Kit for structural ceiling mount 10'11" - 11'8" above the floor.
MI00-1000976-63	Extended Ceiling Rod Kit for structural ceiling mount 11'9" - 12'0" above the floor.
MI00-1000976-67	Extended Ceiling Rod Kit for structural ceiling mount 12'1" - 12'10" above the floor.
MI00-1000976-77	Extended Ceiling Rod Kit for structural ceiling mount 12'11" - 14'0" above the floor.

ChuttleTrak™ Ceiling Track can be used on all models of the 21ST CENTURY Medical Lights. Ceiling Tracks offer exceptional versatility since they allow the horizontal movement of the lights, in addition to the movement that the individual lights offer. Tracks are available with a single trolley which can hold one Ceiling Mount or one Dual Ceiling Mount. The dual trolleys can accommodate a minimum of two Single Ceiling Mounts or a maximum of two Dual Ceiling Mounts. Standard track length is 70-in. (178cm). Other lengths are available, please call for details.

When ordering a ShuttleTrak™, correct height from the floor to the structural ceiling mount must be known at time of ordering. ShuttleTrak™ can only utilize the Rod Kits listed above and lengths up to a 25-in.

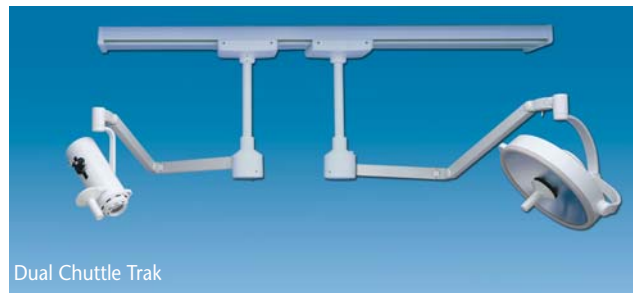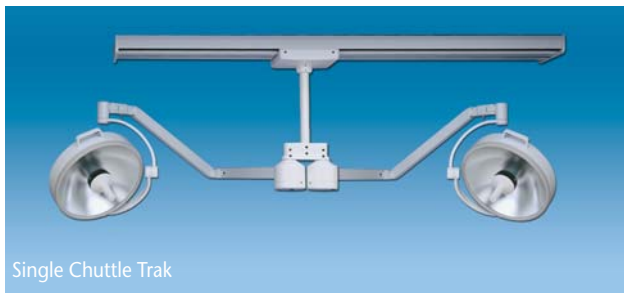


P/N	DESCRIPTION
CEILING TRACKS	
MI000000-044011	CHUTTLE TRAK Single Trolley Model *
MI000000-044012	CHUTTLE TRAK Dual Trolley Model *
* Ceiling Rod kit is required -- only 16" or 25" Rods can be used.	



Century "Woods" Diagnostic Light offers unparalleled versatility for diagnostic and general treatment procedures using two types of bulbs along with magnification. Hands free design in your choice of mounts allows you to position the light where you need it and rest assured it will stay in place. Illumination is provided by two ultraviolet bulbs for diagnostic applications and two daylight fluorescent bulbs with a correct color temperature of 6300° Kelvin for general examinations. Acrylic magnifying lens offers 2x power magnification and 6x power bifocal magnification. Features 24-in. flexible arm and 6-foot power cord.

- Ultraviolet & Fluorescent Light Source
- 6300° Kelvin Color Temperature
- 2x Standard Acrylic Magnifying Lens
- 6x Bifocal Lens
- Available in a Variety of Mounting Options
- Easy to Reach On/Off Switch
- 24" Flexible Arm

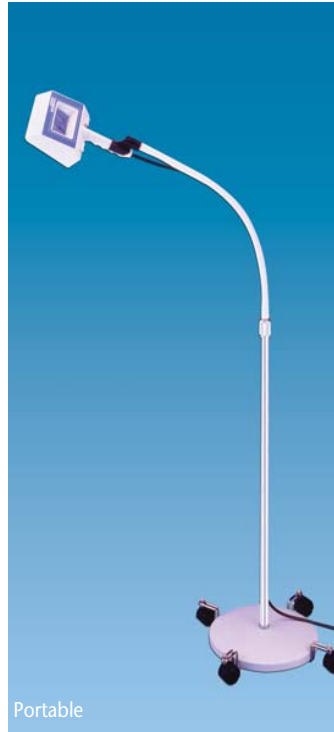

Portable



Hand Held



Table/Shelf Mounted



Wall Mounted

P/N		DESCRIPTION
115V	220V	
MI000000-015010	MI000000-015040	Hand-Held
MI000000-015011	MI000000-015041	Table/Shelf Mount
MI000000-015012	MI000000-015042	Wall Mounted
MI000000-015016	MI000000-015046	Floor Standing
MI000000-015017	MI000000-015047	Portable Floor with Casters
BULBS		
MI00000-0001481	U.V. Daylight Fluorescent Light (White Bulb) (24V/7 watts)	
MI00000-0001482	U.V. Black Light Blue Fluorescent Light (Black Bulb) (24V/7 watts)	

Versalite represents the best value today in general purpose exam lights. It features an amazing 10,000-hour rated bulb life, a break-resistant nylon light head housing, an adjustable upright rod and a cast iron base. It is energy-efficient, produces light output equal to a 130 watt bulb, while using only 20% of the same amount of energy. Color Temperature of over 6300° Kelvin.

- Compact Fluorescent Light Source
- 6300° Kelvin Color Temperature
- Specular Aluminum Reflector
- Lamphead rotates 340°




P/N		DESCRIPTION
115V	220V	
MI000000-021616	MI000000-021646	Floor Model
BULBS		
MI00000-0001579	VERSALITE Compact Fluorescent E26 Bulb (115V/26 watts)	
MI00000-0001587	VERSALITE Compact Fluorescent E27 Bulb (220V/26 watts)	

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061
Hours: 8:30am - 5:00pm CST Monday - Friday
 Visit our website at: www.suburbansurgical.com

MED-EXAM, MEDI-VUE

Med-Exam Light is perfect for examinations, minor procedures and diagnostics.

- Halogen Light Source
- 3700° Kelvin Color Temperature
- 11" Diameter Reflector
21,500 Lux (2000 foot candle output at 36-in/914-mm)
- Heat-Absorbing, Color-Correcting Glass Cylinder
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- 12" Flexible Neck

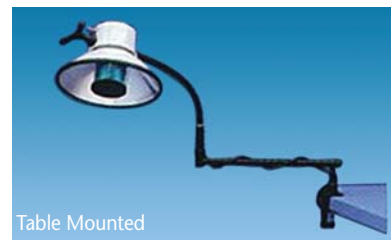

P/N	DESCRIPTION
115V	220V
MI000000-011111	MI000000-011141 Table/Shelf Clamp Mounted
MI000000-011112	MI000000-011142 Wall Mounted
MI000000-011116	MI000000-011146 Floor Standing
MI000000-011117	MI000000-011147 Portable Floor with Casters
ACCESSORIES/BULBS	
MI00000-1001201	MED-EXAM Tee Handle, NOT REMOVABLE, NOT AUTOCLAVABLE
MI00000-0001727	Medium Base Quartz/Halogen Bulb (115V/150 watt)
MI00000-0001728	Medium Base Quartz/Halogen Bulb (220V/150 watt)

Medi-Vue Magnifying Light Imagine being able to magnify and illuminate your treatment field while keeping your hands free to perform the procedure.

- 9-Watt Compact Fluorescent Light Source
- Three diopter lens is standard.
- Optional add-on 4D and 4D+12D (bifocal) lenses
- Lens pivots to give the best view of the treatment field
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- Toggle switch lets you choose one of two light intensities.
- Strong, flexible positioning arm allows you to place magnifying light where you need it
- Handle mounted to the perimeter of light head




P/N	DESCRIPTION
MI000000-045011	MEDI-VUE Table/Shelf Clamp Mounted One Fluorescent Magnifying Light
MI000000-045012	MEDI-VUE Wall Mounted One Fluorescent Magnifying Light
MI000000-045017	MEDI-VUE Portable Floor with Casters One Fluorescent Magnifying Light
ACCESSORIES/BULBS	
MI00000-0001702	MEDI-VUE 4d "Add-on" Lens
MI00000-0001703	MEDI-VUE 4d + 12d Bifocal "Add-on" Lens
MI00000-0001704	Compact Fluorescent Bulb (120V/9 watts)



Medi-Lite Light represents all the advantages of the new Medi-Light series. Considered one of the best values in the industry.

- Halogen Light Source
- 3700° Kelvin Color Temperature
- 13" Diameter Reflector
- 37,700 Lux (3500 foot candle output at 36-in/914-mm)
- Heat-Absorbing, Color-Correcting Glass Cylinder
- Available in a Variety of Mounting Options
- Easy to reach On/Off Switch
- 42" Arm Extension

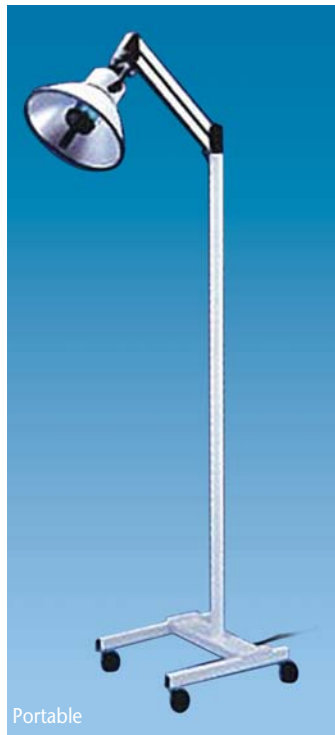

Portable



Dual Ceiling Mount



Ceiling Mounted



Wall Mounted

P/N	DESCRIPTION
115V	220V
MI000000-011312	MI000000-011342 Wall Mounted
MI000000-011313	MI000000-011343 Portable Floor with Casters
MI000000-011314	MI000000-011344 Ceiling Mounted One
MI000000-011315	MI000000-011345 Ceiling Mounted Two

Ceiling Rod kit is required. Correct height from the floor to the structural ceiling mount must be known at time of ordering.

CEILING ROD KITS

Basic Kit Includes: Ceiling Casting, Cover, Ceiling Rod & Wiring Enclosure. Mounted Ring must be attached to structural ceiling mount.

MI00-1000908-10	10" Ceiling Rod Kit for structural ceiling mount 7'6" - 7'11" above the floor.
MI00-1000908-15	15" Ceiling Rod Kit for structural ceiling mount 8'0" - 8'6" above the floor.
MI00-1000908-25	25" Ceiling Rod Kit for structural ceiling mount 8'7" - 9'6" above the floor.
MI00-1000908-35	35" Ceiling Rod Kit for structural ceiling mount 9'7" - 10'6" above the floor.
MI00-1000908-45	45" Ceiling Rod Kit for structural ceiling mount 10'7" - 11'6" above the floor.
MI00-1000908-55	55" Ceiling Rod Kit for structural ceiling mount 11'7" - 12'6" above the floor.

Single Ceiling Kit Includes: Ceiling Casting, Cover, Ceiling Rod, Wiring Enclosure, Extension Arm, Hinge Pin, Caps, Washers, Ceiling Rod Collar. Mounted Ring must be attached to structural ceiling mount.

MI00-1000019-10	10" Ceiling Rod Kit for structural ceiling mount 7'6" - 7'11" above the floor.
MI00-1000019-15	15" Ceiling Rod Kit for structural ceiling mount 8'0" - 8'6" above the floor.
MI00-1000019-25	25" Ceiling Rod Kit for structural ceiling mount 8'7" - 9'6" above the floor.
MI00-1000019-35	35" Ceiling Rod Kit for structural ceiling mount 9'7" - 10'6" above the floor.
MI00-1000019-45	45" Ceiling Rod Kit for structural ceiling mount 10'7" - 11'6" above the floor.
MI00-1000019-55	55" Ceiling Rod Kit for structural ceiling mount 11'7" - 12'6" above the floor.

ACCESSORIES/BULBS

MI000-1001006-2	MEDI-LITE SERIES Sterilizable Handle
MI000000-0001727	Medium Base Quartz/Halogen Bulb (115V/150 watt)
MI000000-0001728	Medium Base Quartz/Halogen Bulb (220V/150 watt)