

"LEADER IN ANIMAL CARE PRODUCTS FOR OVER 60 YEARS"



SSCI


SUBURBAN SURGICAL Co., Inc.

275 Twelfth Street, Wheeling, IL 60090

All intellectual material, including the text, graphics and other components of this catalog are Copyrighted by SSCI.
Any use or duplication of these intellectual materials, without written consent by SSCI, is strictly forbidden.

LAMINATED CAGES


ARISTOCRAT: DEVELOPED AND AVAILABLE FIRST AT SSCI — INDIVIDUAL MODULE DESIGN OFFERS ARRANGEMENT FLEXIBILITY

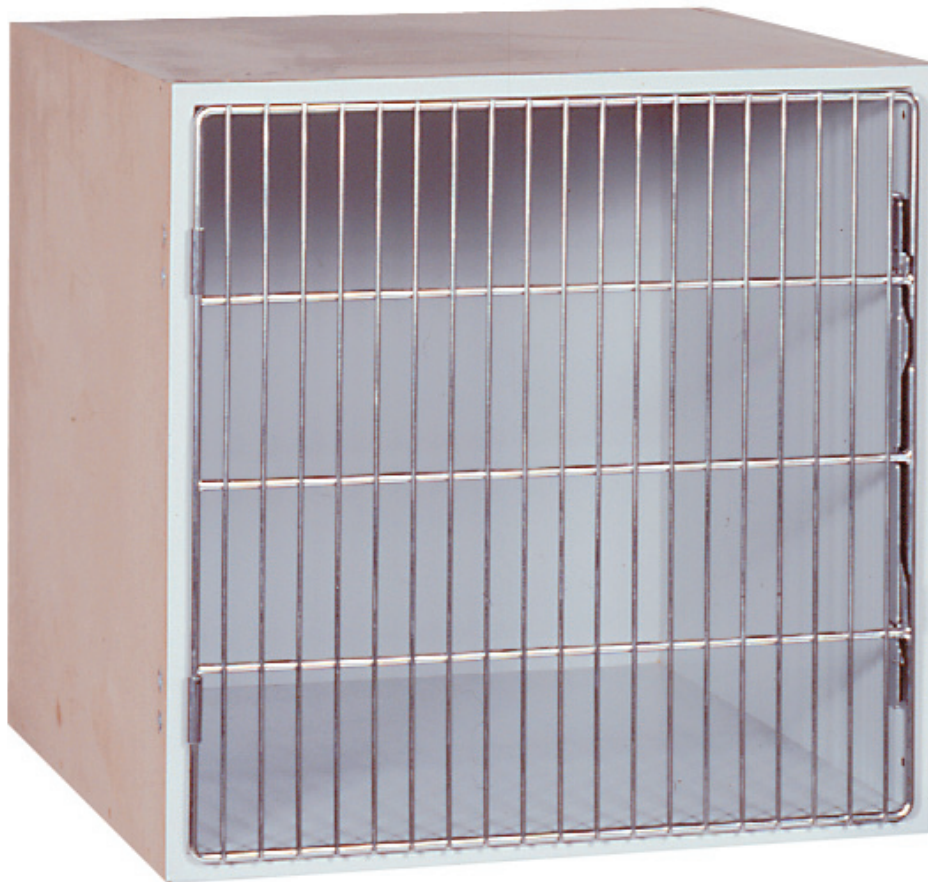

ARISTOCRAT CAGES ARE AVAILABLE IN A VARIETY OF STYLES:



Double Door Aristocrat Cage

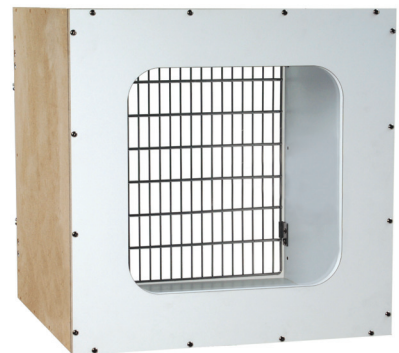

Full-View Back Aristocrat Cage

AN ATTRACTIVE ALTERNATIVE TO STAINLESS STEEL

If your facility needs do not require the ultimate durability of stainless steel, or you prefer to add the choice of color, our Aristocrat Cages are the obvious choice. Aristocrat Cages beat the competition by offering many of the unique, high-quality design features of our Stainless Steel Regal Cages in virtually any color you desire.

Each Aristocrat Laminated Cage begins with specially designed water-resistant materials. All interior surfaces and the front edge of each cage are covered with high-pressure laminate before being assembled. Once each panel is fastened together to form the final cage, each joint is epoxy sealed to create a long-lasting fluid barrier. This unique construction design makes our cages extra resistant to fluid damage.

The doors of our Aristocrat Cages are made entirely of stainless steel, electro-polished, the frame, crossmembers and vertical spokes are welded together at every intersection for exceptional structural integrity and strength. Our cage door mounts at the front of the cage, fitting close along the top and sides and overlapping the bottom of the cage opening to keep the animal's paws from becoming trapped between the door and the cage. Our distinctive door design also prevents paper liners from slipping out under the door.

Our latches and hinges are made from 14-gauge stainless steel. The adjustable, self-latching door latch positively secures the door at two points without creating crevices for debris to hide. Latch disassembles easily for adjustment and cleaning. Hinges and latch brackets mount inside the cage opening and are bolted through the side of the cage into a concealed threaded connector for greater structural integrity. Our cage doors are ambidextrous and are shipped with doors hinged left but may be easily changed to be hinged right if desired without additional parts.

All Aristocrat Cages can be used alone or readily assembled together to form the right cage arrangement for your facility. We ship each cage individually to keep shipping costs low. All installation hardware is included with your cage purchase. Better yet, Aristocrat Cage Arrangements can be easily added to or rearranged at any time.

Aristocrat Cage Arrangements may be made with any style Aristocrat Cage.



Aristocrat Cage Shelf Options

Make your cages more versatile with our laminated Shelves. Shelves are available for cages up to 36 in. wide. Shelves can be used to make room for multiple animals in one cage, or simply make cats more comfortable by allowing them to climb and rest higher than the floor. Cages can be purchased with factory-installed shelves or shelves can be purchased separately for installation by customer on new or existing cages. Shelves are mechanically-fastened thru the sides of the cage for stability. This prevents the shelf from ever tipping out of position. See Aristocrat Cage Finishing Accessories.

Aristocrat Cages with factory-installed shelves are positioned tight to the back spanning the width of the cage and the shelf is midway between the floor and the ceiling of the cage. *Select which style Aristocrat Cage you want from the following pages, then see the footnote below each chart to specify cages with shelves.*



Corner Aristocrat Cage



Full-View Back Corner Aristocrat Cage

Aristocrat Cages

ALL SSCI LAMINATED ARISTOCRAT CAGES ARE DESIGNED WITH THESE QUALITY FEATURES:

Door

Stainless steel door spokes are deburred smooth to eliminate sharp ends.

Stainless steel, electropolished removable ambidextrous door may be remounted with hinges on the right side with no additional parts.

Door extends below cage floor to help keep paper liners in place.

Door fits snug to the cage body to eliminate areas where an animal's paws could be pinched.

.1875-in. (.48cm) diameter stainless steel vertical spokes are spaced 1-in. apart and welded at every intersection to crossmembers and frame.

.3125-in. (.79cm) diameter stainless steel door frame and crossmembers are welded together for strength.

Latch

Stainless Steel Latch easily disassembles for sanitizing and adjusting.

.3125-in. (.79cm) diameter stainless steel latch keeper latches door securely.

Heavy 12-gauge stainless steel latch brackets.

Latch connector bar is notched for easy opening.

Heavy 12-gauge stainless steel self-latching door latch connector is designed to accommodate a padlock for greater security.

Hinges

Heavy 12-gauge stainless steel hinges allows the door to operate smoothly.

Hinges are thru-bolted with a threaded connector to sides of cage using a 'T'-nut

All holes are precision-machined with a computer-controlled boring machine. Components are precision-milled in a computer-controlled milling machine.


Floor seams are sealed with an epoxy sealant. The sealant extends up the side joints for a few inches.

Distinctive floor design with fluid-retaining lip is standard on all Aristocrat Cages.

Plastic laminate is bonded to interior and front surfaces before assembly with a high-pressure press to create a superior moisture barrier.

Body

Constructed of specially designed water-resistant .625 in. thick substrate. Glued and mechanically-fastened at all joints for extra strength. All joints are bonded with a glue that is impervious to liquid.

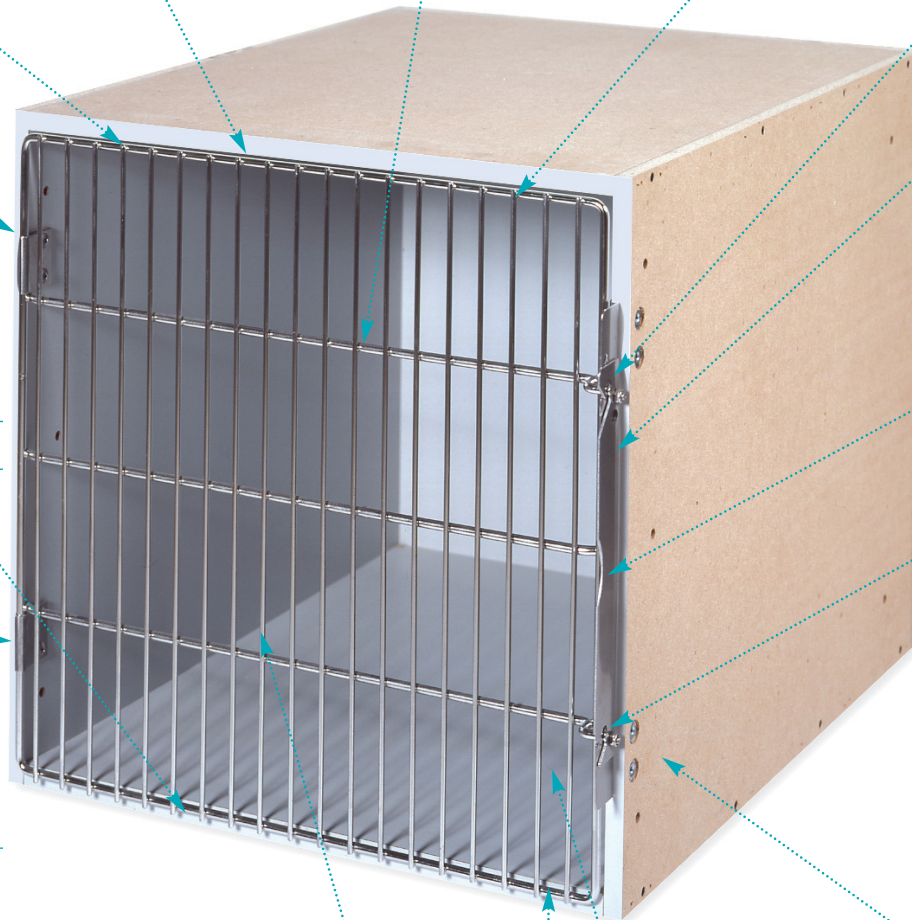


Aristocrat Cages

Two smaller doors allow you to have a large cage in a smaller space — less clearance space is needed to open doors all the way.

14-gauge stainless steel latch: top and bottom plate latches door in two places.

Stainless steel channel holds Divider Panel in place.

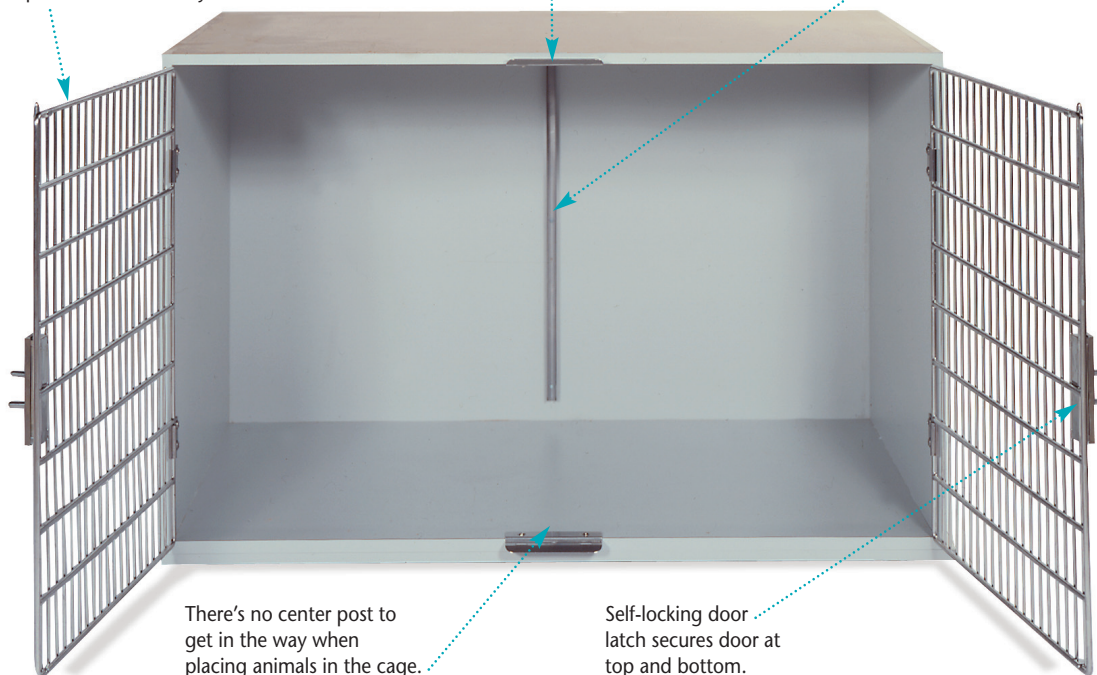

Aristocrat Cages

P/N	CAGE SIZE	
	Inches	Centimeters
14000-00-CTCTEI-__-†	18W x 18H	45.72W x 45.72H
14000-00-CTDREI-__-†	24W x 18H	60.96W x 45.72H
14000-00-DRCTEI-__-†	18W x 24H	45.72W x 60.96H
14000-00-DRDREI-__-†	24W x 24H	60.96W x 60.96H
14000-00-DREPEI-__-†	30W x 24H	76.20W x 60.96H
14000-00-DRFNEI-__-†	36W x 24H	91.44W x 60.96H
14001-00-DRHJEI-__-†	48W x 24H DD*	121.92W x 60.96H DD*
14000-00-EPDREI-__-†	24W x 30H	60.96W x 76.20H
14000-00-EPEPEI-__-†	30W x 30H	76.20W x 76.20H
14000-00-EPFNEI-__-†	36W x 30H	91.44W x 76.20H
14000-00-EPGLEI-__-†	42W x 30H	106.68W x 76.20H
14000-00-EPHJEI-__-†	48W x 30H	121.92W x 76.20H
14001-00-EPHJEI-__-†	48W x 30H DD*	121.92W x 76.20H DD*
14000-00-FNDREI-__-†	24W x 36H	60.96W x 91.44H
14000-00-FNEPEI-__-†	30W x 36H	76.20W x 91.44H
14000-00-FNFNEI-__-†	36W x 36H	91.44W x 91.44H
14000-00-FNGLEI-__-†	42W x 36H	106.68W x 91.44H
14000-00-FNHJEI-__-†	48W x 36H	121.92W x 91.44H
14001-00-FNHJEI-__-†	48W x 36H DD*	121.92W x 91.44H DD*

*DD after the height indicates a double door cage.

All cage bodies measure 28.25-in. (71.75cm) front to back.

Cages measure 30.188-in. (76.68cm) deep overall, including latch.

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.

Double Door (DD) Aristocrat Cages

Our double door cages are designed for easy use and animal comfort with no center door-latching post. Large animals can be easily placed inside without having to maneuver around a center post. Doors have self-latching pins that slide into top and bottom latches. An optional stainless steel divider panel can be used when your needs require two smaller cage compartments. Channel included with each Double Door Cage attaches to back of cage to stabilize Divider Panel. See Aristocrat Cage Finishing Accessories.

**AVAILABLE IN 19 DIFFERENT
SIZES—MORE THAN
ANYWHERE ELSE.**

To order cages with factory-installed Resting Shelf, change the cage P/N 7th position 0, to: A Straight resting shelf

B Corner resting shelf in left back corner

C Corner resting shelf in right back corner

D Corner resting shelves in left and right back corners

Specify the height of the resting shelves above the cage floor, change the cage

P/N, 8th and 9th position 0- to: BF approximately midway in 18-in high cages

BR approximately midway in 24-in high cages

CD approximately midway in 30-in high cages

CP approximately midway in 36-in high cages

See dimension chart for other Resting Shelf heights



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.


For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061

Hours: 8:30am - 5:00pm CST Monday - Friday

Visit our website at: www.suburbansurgical.com

Corner Aristocrat Cages

Designed to fit snugly into corner to give your Arrangement a built-in appearance.

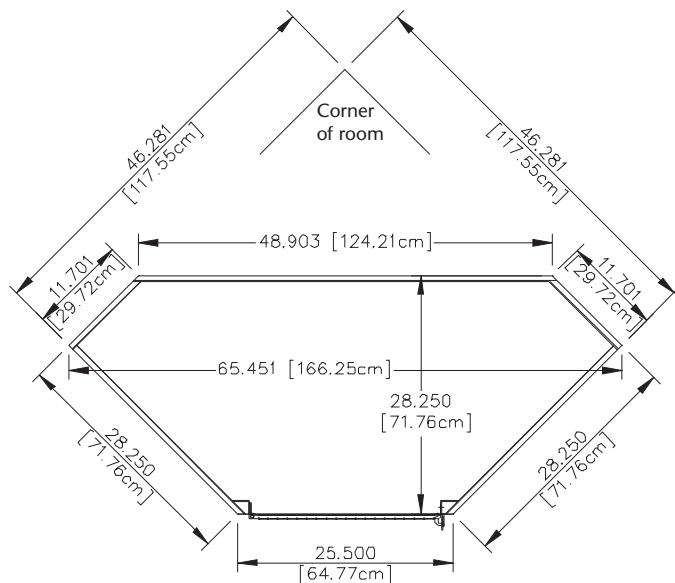

Corner Aristocrat Cages

These versatile, attractive cages are designed to fit perfectly in a corner. Allows you to create an L-shaped arrangement of cages when you place two straight banks against both sides of the Corner Cage.

P/N	DOOR
15000- 00-CTDREI-___†	24W x 18H
15000- 00-DRDREI-___†	24W x 24H
15000- 00-EPDREI-___†	24W x 30H
15000- 00-FNDREI-___†	24W x 36H

All cage bodies measure 28.250-in. (71.76cm) front to back.
Modules measure 30.182-in. (76.66cm) deep overall, including latch.

To order cages with factory-installed Resting Shelf, change the cage P/N 7th position 0, to: A Straight resting shelf
B Corner resting shelf in left back corner
C Corner resting shelf in right back corner
D Corner resting shelves in left and right back corners
Specify the height of the resting shelves above the cage floor, change the cage P/N, 8th and 9th position 0- to: BF approximately midway in 18-in high cages
BR approximately midway in 24-in high cages
CD approximately midway in 30-in high cages
CP approximately midway in 36-in high cages

See dimension chart for other Resting Shelf heights




SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:
†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.

Full-View Back Aristocrat Cages

Clear polycarbonate panel.

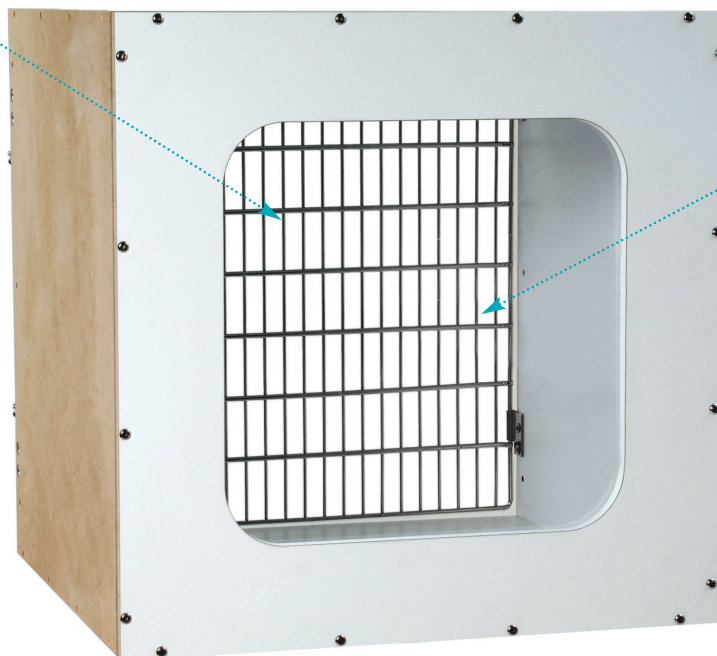

Full-width and height clear back panel offers an unobstructed view to allow observation without direct interaction in post-op, boarding and adoption areas.

P/N	CAGE SIZE	
	Inches	Centimeters
14002-00-CTCTEI-___†	18W x 18H	45.72W x 45.72H
14002-00-CTDREI-___†	24W x 18H	60.96W x 45.72H
14002-00-DRCTEI-___†	18W x 24H	45.72W x 60.96H
14002-00-DRDREI-___†	24W x 24H	60.96W x 60.96H
14002-00-DREPEI-___†	30W x 24H	76.20W x 60.96H
14002-00-DRFNEI-___†	36W x 24H	91.44W x 60.96H
14003-00-DRHJEI-___†	48W x 24H DD*	121.92W x 60.96H DD*
14002-00-EPDREI-___†	24W x 30H	60.96W x 76.20H
14002-00-EPEPEI-___†	30W x 30H	76.20W x 76.20H
14002-00-EPFNEI-___†	36W x 30H	91.44W x 76.20H
14002-00-EPGLEI-___†	42W x 30H	106.68W x 76.20H
14002-00-EPHJEI-___†	48W x 30H	121.92W x 76.20H
14003-00-EPHJEI-___†	48W x 30H DD*	121.92W x 76.20H DD*
14002-00-FNDREI-___†	24W x 36H	60.96W x 91.44H
14002-00-FNEPEI-___†	30W x 36H	76.20W x 91.44H
14002-00-FNFNEI-___†	36W x 36H	91.44W x 91.44H
14002-00-FNGLEI-___†	42W x 36H	106.68W x 91.44H
14002-00-FNHJEI-___†	48W x 36H	121.92W x 91.44H
14003-00-FNHJEI-___†	48W x 36H DD*	121.92W x 91.44H DD*

*DD after the height indicates a double door cage.

All cage bodies measure 28.25-in. (71.75cm) front to back.

Cages measure 30.188-in. (76.68cm) deep overall, including latch.

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.

Full-View Back Aristocrat Cages

These versatile cages are ideal for medical observation or for facilities where animals are on display to the public. Clear polycarbonate back panel allows animals to be seen from either side of the cage. These cages have all the fine features of our Aristocrat Cages with the addition of our uniquely designed resilient polycarbonate back panel.

To order cages with factory-installed Resting Shelf, change the cage P/N 7th position 0, to:

A Straight resting shelf

B Corner resting shelf in left back corner

C Corner resting shelf in right back corner

D Corner resting shelves in left and right back corners

Specify the height of the resting shelves above the cage floor, change the cage

P/N, 8th and 9th position 0- to: BF approximately midway in 18-in high cages

BR approximately midway in 24-in high cages

CD approximately midway in 30-in high cages

CP approximately midway in 36-in high cages

See dimension chart for other Resting Shelf heights



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061


Hours: 8:30am - 5:00pm CST Monday - Friday

Visit our website at: www.suburbansurgical.com

Full-View Back Corner Aristocrat Cage

Clear polycarbonate panel offers an unobstructed view of the animal inside.

Designed to fit snugly into corner to give your Arrangement a built-in appearance.

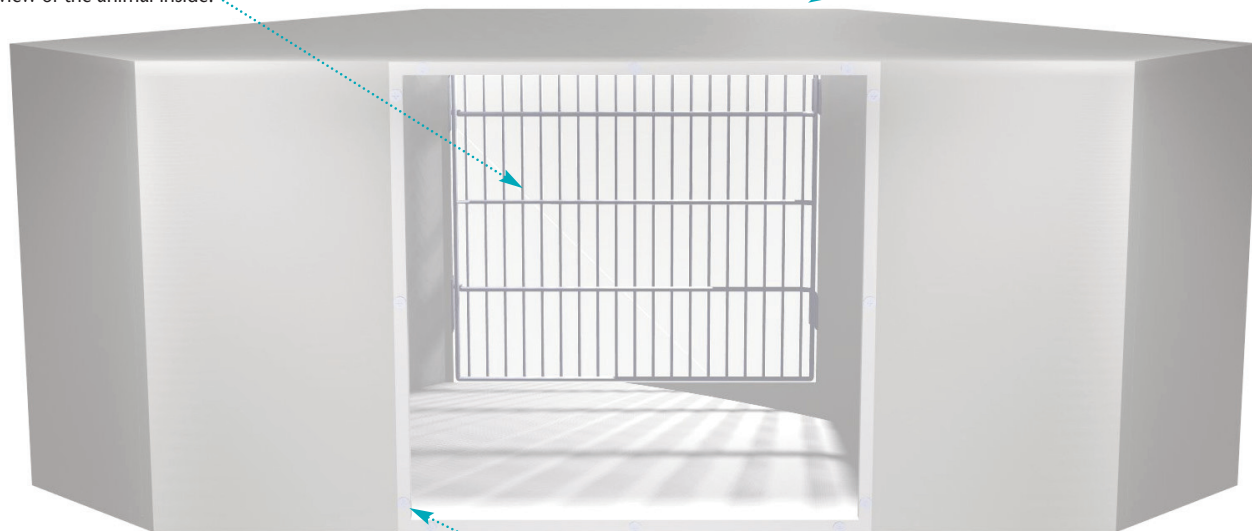

Sealer cushions panel, providing a tight seal.

Full-View Back Corner Aristocrat Cages

These versatile, attractive cages are designed to fit perfectly in a corner. Cage features a clear back panel to allow animals to be seen from either side of the cage. Allows you to create an L-shaped arrangement of cages when you place two straight banks against both sides of the Corner Cage. These cages have all the fine features of our Aristocrat Corner Cages with the addition of our uniquely designed polycarbonate back panel.

P/N	DOOR
15002- 00-CTDREI-___†	24W x 18H
15002- 00-DRDREI-___†	24W x 24H
15002- 00-EPDREI-___†	24W x 30H
15002- 00-FNDREI-___†	24W x 36H

To order cages with factory-installed Resting Shelf, change the cage P/N 7th position 0, to: A Straight resting shelf

B Corner resting shelf in left back corner

C Corner resting shelf in right back corner

D Corner resting shelves in left and right back corners

Specify the height of the resting shelves above the cage floor, change the cage


P/N, 8th and 9th position 0- to: BF approximately midway in 18-in high cages

BR approximately midway in 24-in high cages

CD approximately midway in 30-in high cages

CP approximately midway in 36-in high cages

See dimension chart for other Resting Shelf heights

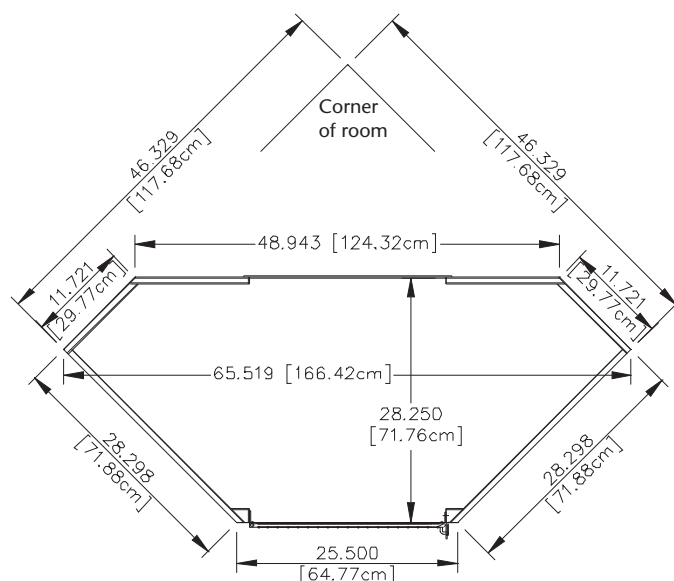

All cage bodies measure 28.25-in. (71.75cm) front to back.

Cages measure 30.188-in. (76.68cm) deep overall, including latch.

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Aristocrat Cage Arrangements

Shown with optional End and Top Panels.

Individual cages fit tightly together. Cages are shipped separately to reduce your shipping costs.

All hardware is included for quick and easy installation.



*For
Bases, Top and
End Panels, see
Aristocrat Cage
Finishing
Accessories.*

**Aristocrat Cage Arrangement
shown with Top and End Panels**

Each and every Aristocrat Cage can be combined with other Aristocrat Cages to create the ideal Arrangement for your facility. Better yet, if your needs or available space change, the Arrangement can be added to or rearranged to fit the new requirements.

You can reverse the ambidextrous doors to open the other way with no additional hardware if you find it would work better for your needs. And when you add some of our finishing accessories, you'll find that Aristocrat Cages are not only the best, but also the most versatile laminated cages you can buy.

Arrangements may be mounted on a SSCI Laminated Toe-Kick Base, stainless steel Base with or without casters, or on your wooden or concrete platform.

A variety of Cage Arrangements are shown on the following pages. One may be ideal for your needs. However if you don't see an Arrangement that fits your facility, we will work with you to create the perfect configuration of cage types and sizes for your needs.

Design Service For Aristocrat Cage Arrangement

THE PERFECT ARISTOCRAT CAGE ARRANGEMENT IS JUST A PHONE CALL AWAY

SSCI's Design Service is offered to customers after discussing and agreeing on a viable project with an SSCI Representative. After a methodical review of the customer's needs for their facility, the SSCI Representative will present the project to our staff of designers to create a set of approval drawings. The drawings ensure that SSCI will supply exactly what the customer has discussed with the

SSCI Representative. Once the customer has committed to purchase from SSCI and, if the drawings require revisions, the Final Approval Drawings will be created and presented to the customer for their approval. Other Customer Requirements may be required before manufacturing time can begin.

CREATING THE PERFECT ARISTOCRAT CAGE ARRANGEMENT TAKES JUST 5 EASY STEPS:

1 Measure your space.

Obtain overall width, height and depth of the available space in your facility where you plan to place your Aristocrat Cage Arrangement. To allow for the cage, the minimum depth of the area must be 29-in. (73.66cm) plus the width of the widest cage door in the arrangement. This will allow cage doors to open fully. Make sure you have sufficient space in the front of the cage arrangement for your staff to handle patients safely.

2 Choose the style of the Aristocrat Cages in your Arrangement.

Our Aristocrat Cage styles give you two options.

3 Identify the Aristocrat Cage sizes you require.

A variety of sizes can be used in a single Arrangement, as long as the overall width and height of each tier is the same. An idea of the size of animals you will house in your Cages will help you select the sizes you need.

Determine how your Cages will be arranged.

4

Some Aristocrat Cage Arrangements are shown on the following pages. If you don't see an Arrangement that will work for you, call us with your space measurements, and we will help you design a Aristocrat Cage Arrangement that will include the largest cages for the total number of animals you wish to accommodate.

Aristocrat Cage Arrangements should be set on a base.

5

Optional SSCI Stainless Steel Base:

(add casters to create a mobile base).

Customer-Supplied Bases:


Wood Platform Bases should be constructed level front to back and right to left of 2x4s measuring the width of your Arrangement and 29-in. (73.66cm) front to back.

Trowel-finished Concrete Pad should be poured level 3.5-in. (8.90cm) high and 29-in. (73.66cm) front to back.

If one of our illustrated Arrangements shown on the following pages won't fit your facility perfectly, let our designers help you create the ideal Arrangement for your needs. SSCI's design team will work with your specifications to create a Cage Arrangement suited to your unique needs, space requirements and budget. And if you already own some cages, even some of our competitors' cages, we can often incorporate those into your Cage Arrangement design.


It all starts with a phone call. Contact one of our repre-

sentatives to discuss your requirements. We'll create a CAD-generated drawing of your Arrangement complete with all hardware needed for installation.

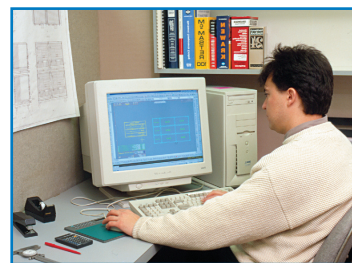

Our team of experienced designers will work with your measurements, cage style and size and budget to create the perfect Cage Arrangement for your requirements.


FOR EVERY ARISTOCRAT ARRANGEMENT PURCHASE SSCI CREATES A DRAWING TO ENSURE CUSTOMER'S SATISFACTION

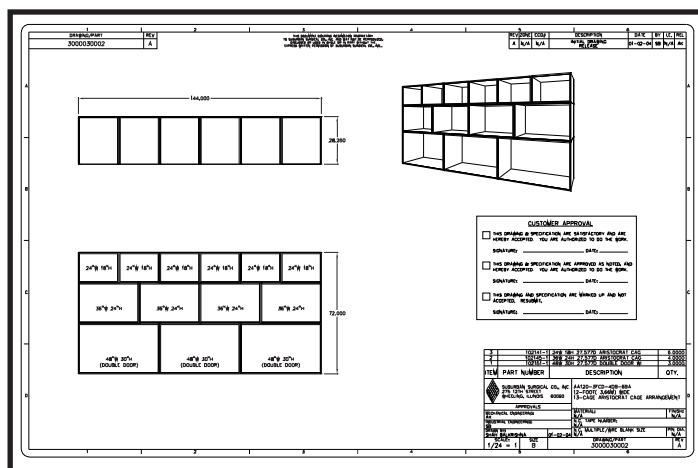

SSCI creates an Approval Drawing for each cage arrangement ordered. These detailed drawings of your cage arrangement list cage sizes, quantities, laminate choice and all appropriate hardware and illustrate how the arrangement is to be assembled. If you choose to alter the positions of the cages in the arrangement or your choice of laminate, do not sign the existing Approval Drawing. Contact your SSCI representative, so that a new Approval Drawing can be created.

All laminated products are non-stocked, made specifically per your specifications. Once SSCI has received the signed Approval Drawing the order is non-cancelable. In addition to other customer requirements, the Final Approval Drawings must be signed and received by SSCI prior to the beginning of manufacturing lead time.

Aristocrat Cage Arrangements

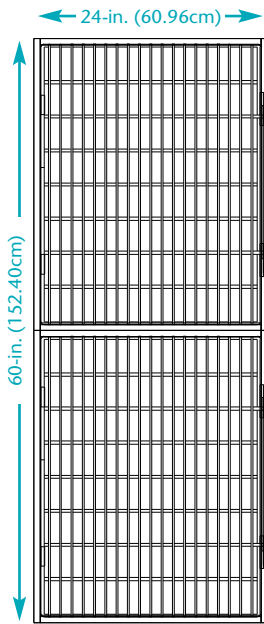

2-Foot (.61m) Wide, 2-Cage Aristocrat Cage Arrangement


Includes:

Top:

One 24-in.W x 30-in.H (60.96cmW x 76.20cmH) Cage

Bottom:

One 24-in.W x 30-in.H (60.96cmW x 76.20cmH) Cage

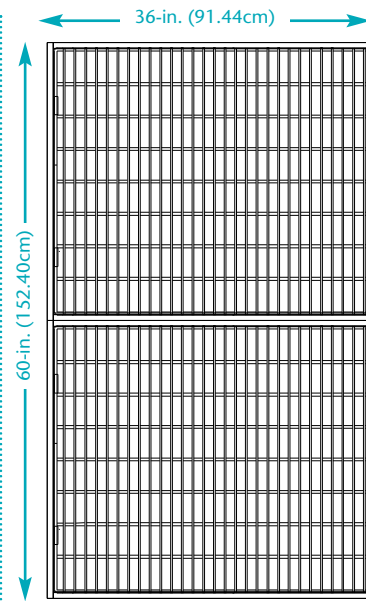

3-Foot (.91m) Wide, 2-Cage Aristocrat Cage Arrangement

Includes:

Top:


One 36-in.W x 30-in.H (91.44cmW x 76.20cmH) Cage

Bottom:

One 36-in.W x 30-in.H (91.44cmW x 76.20cmH) Cage

P/N	Description
3000030020-___†	Aristocrat Cage Arrangement
3000033000-___†	Full-View Polycarbonate Back Aristocrat Cage Arrangement

P/N	Description
3000030021-___†	Aristocrat Cage Arrangement
3000033001-___†	Full-View Polycarbonate Back Aristocrat Cage Arrangement

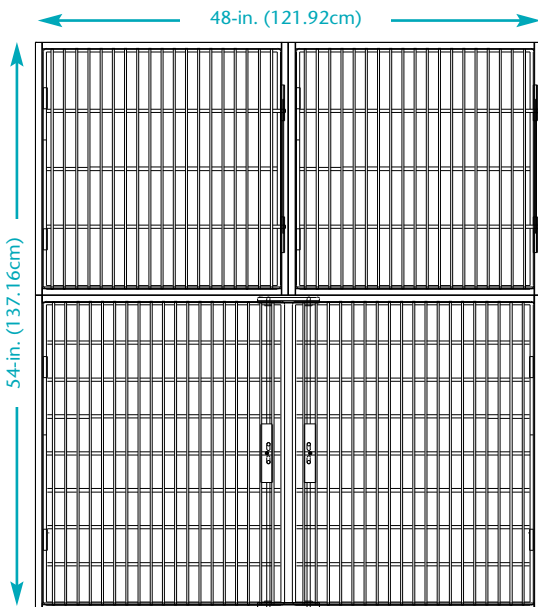

4-Foot (1.22m) Wide, 3-Cage Aristocrat Cage Arrangement

Includes:

Top: Two 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages

Bottom: One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage

P/N	Description
3000030014-___†	Aristocrat Cage Arrangement
3000033002-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.




SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061

Hours: 8:30am - 5:00pm CST Monday - Friday

Visit our website at: www.suburbansurgical.com

Aristocrat Cage Arrangements

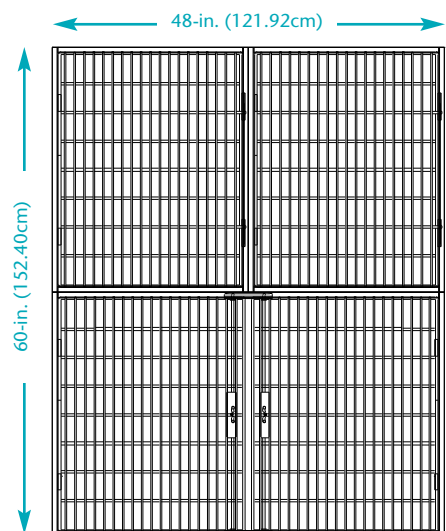


4-Foot (1.22m) Wide, 3-Cage Aristocrat Cage Arrangement

Includes:

Top: Two 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cages

Bottom: One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage

P/N	Description
3000030013-___†	Aristocrat Cage Arrangement
3000033003-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

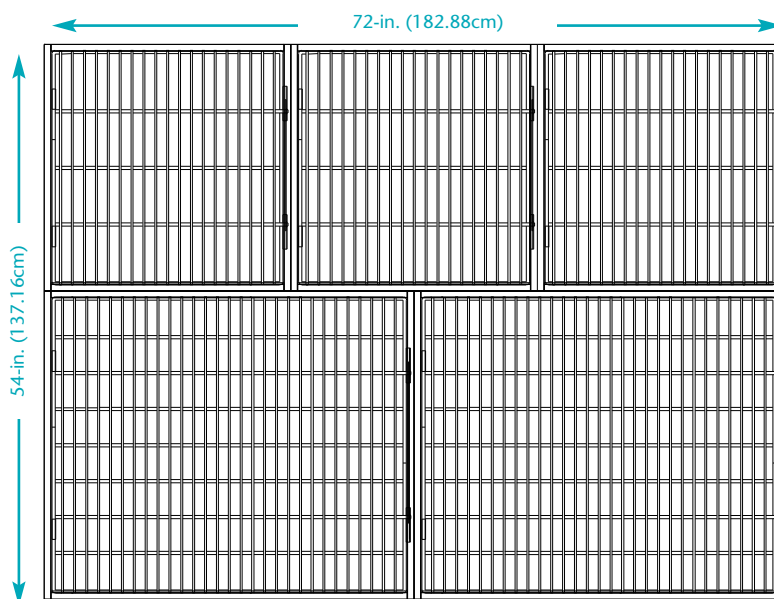

6-Foot (1.83m) Wide, 5-Cage Aristocrat Cage Arrangement

Includes:

Top: Three 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages

Bottom: Two 36-in.W x 30-in.H (91.44cmW x 76.20cmH) cages

P/N	Description
3000030000-___†	Aristocrat Cage Arrangement
3000033004-___†	Full-View Polycarbonate Aristocrat Cage Arrangement


Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Aristocrat Cage Arrangements

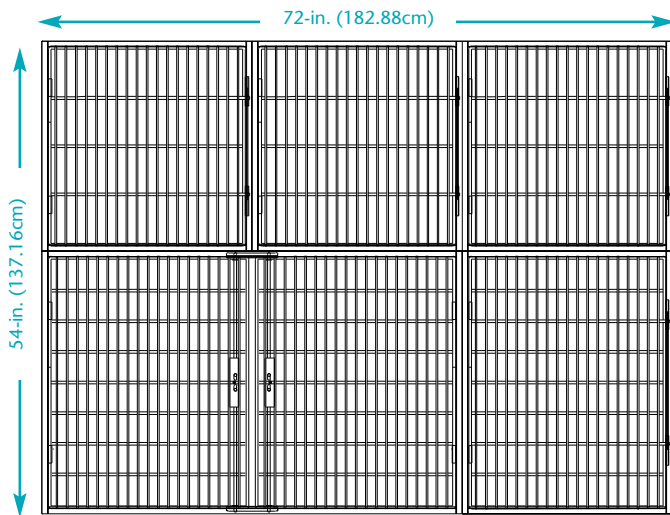

6-Foot (1.83m) Wide, 5-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Bottom:** One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage
One 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cage


P/N	Description
3000030022-___†	Aristocrat Cage Arrangement
3000033005-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

6-Foot (1.83m) Wide, 5-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Three 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cages
- Bottom:** Two 36-in.W x 30-in.H (91.44cmW x 76.20cmH) cages

P/N	Description
3000030018-___†	Aristocrat Cage Arrangement
3000033006-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

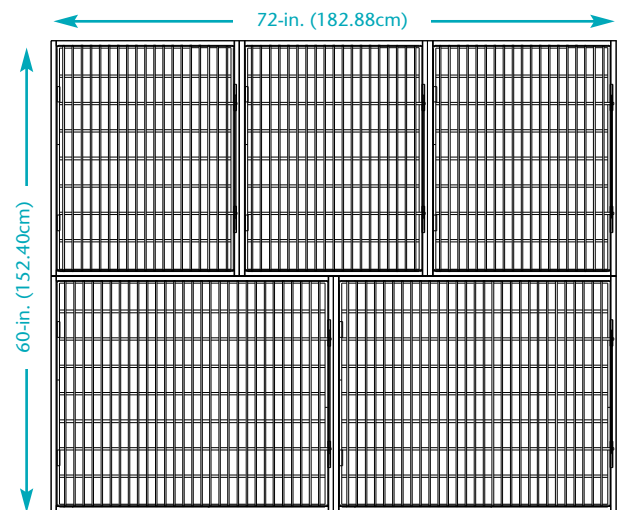


6-Foot (1.83m) Wide, 8-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cm) cages
- Middle:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cm) cages
- Bottom:** Two 36-in.W x 30-in.H (91.44cmW x 76.20cm) cages

P/N	Description
3000030023-___†	Aristocrat Cage Arrangement
3000033007-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

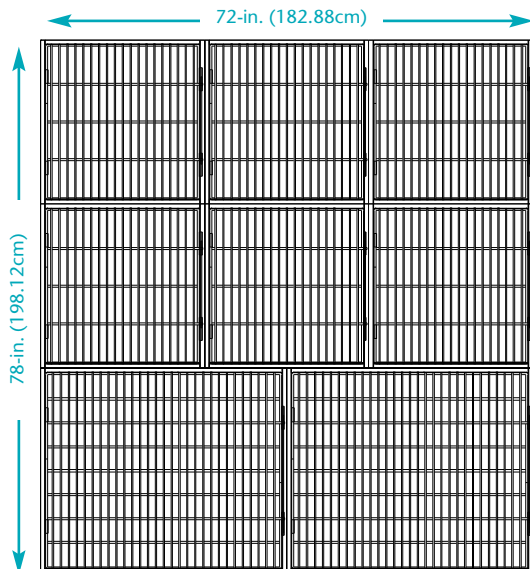


Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:
†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Aristocrat Cage Arrangements

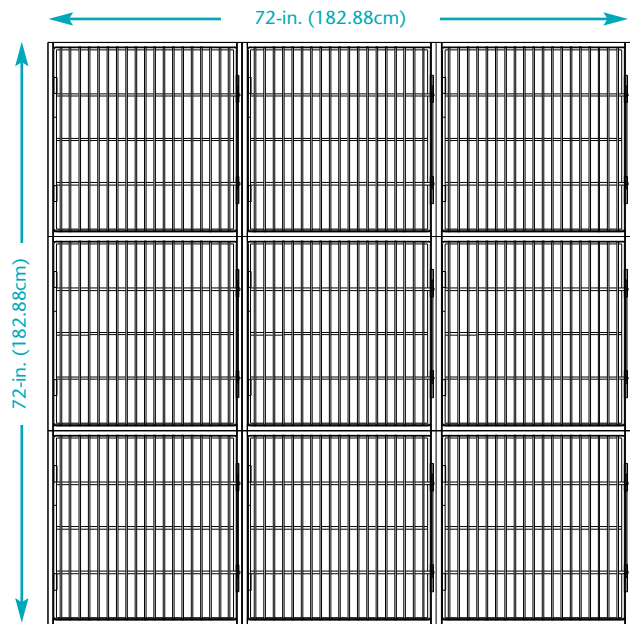

6-Foot (1.83m) Wide, 9-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Middle:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Bottom:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages


P/N	Description
3000030024-___†	Aristocrat Cage Arrangement
3000033008-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

8-Foot (2.44m) Wide, 9-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Four 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Middle:** Two 36-in.W x 30-in.H (91.44cmW x 76.20cmH) cages
One 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cage
- Bottom:** Two 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cages

P/N	Description
3000030025-___†	Aristocrat Cage Arrangement
3000033009-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

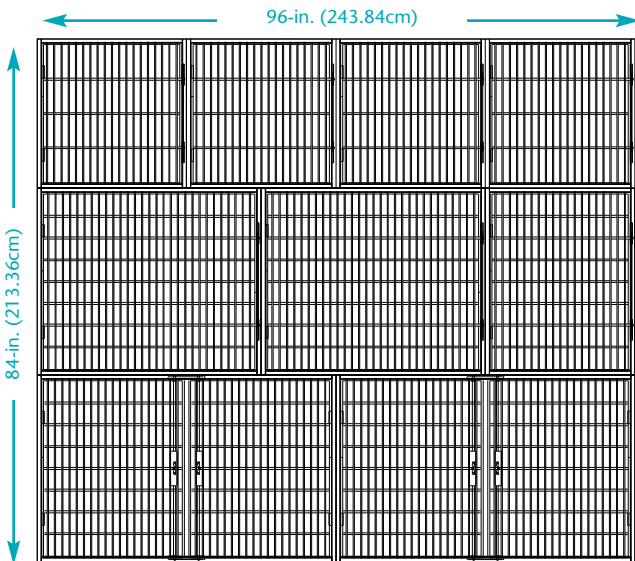

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:
†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.



Aristocrat Cage Arrangements

8-Foot (2.44m) Wide, 10-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Four 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
Middle: Four 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
Bottom: Two 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cages

P/N	Description
3000030026-___†	Aristocrat Cage Arrangement
3000033010-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

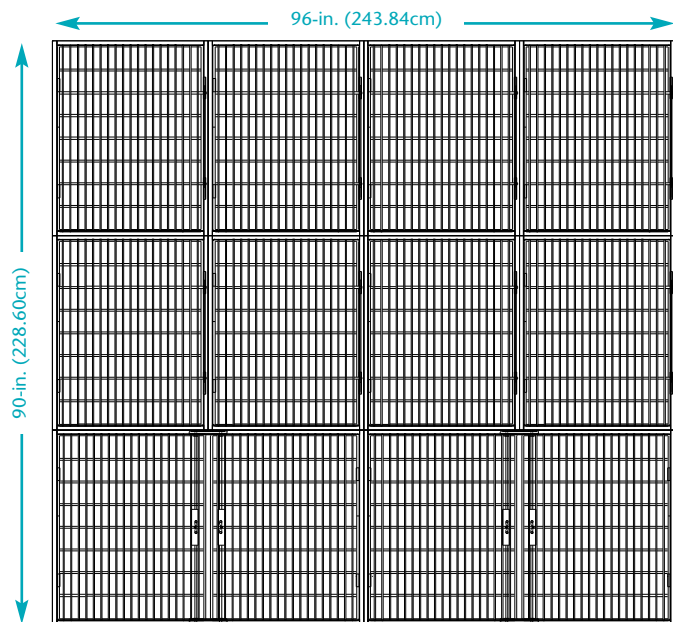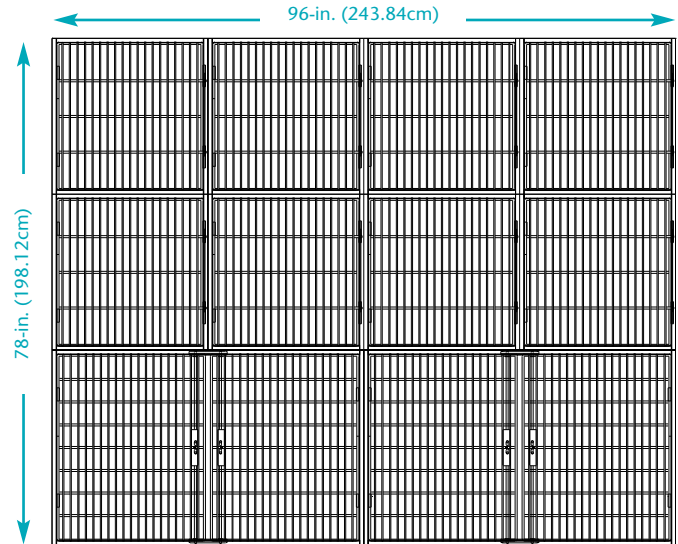

8-Foot (2.44m) Wide, 10-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Four 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cages
Middle: Four 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cages
Bottom: Two 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cages

P/N	Description
3000030027-___†	Aristocrat Cage Arrangement
3000033011-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061

Hours: 8:30am - 5:00pm CST Monday - Friday


Visit our website at: www.suburbansurgical.com

Aristocrat Cage Arrangements

10-Foot (3.05m) Wide, 12-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Five 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Middle:** Four 30-in.W x 30-in.H (76.20cmW x 60.96cmH) cages
- Bottom:** One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage
Two 36-in.W x 30-in.H (91.44cmW x 76.20cmH) cages

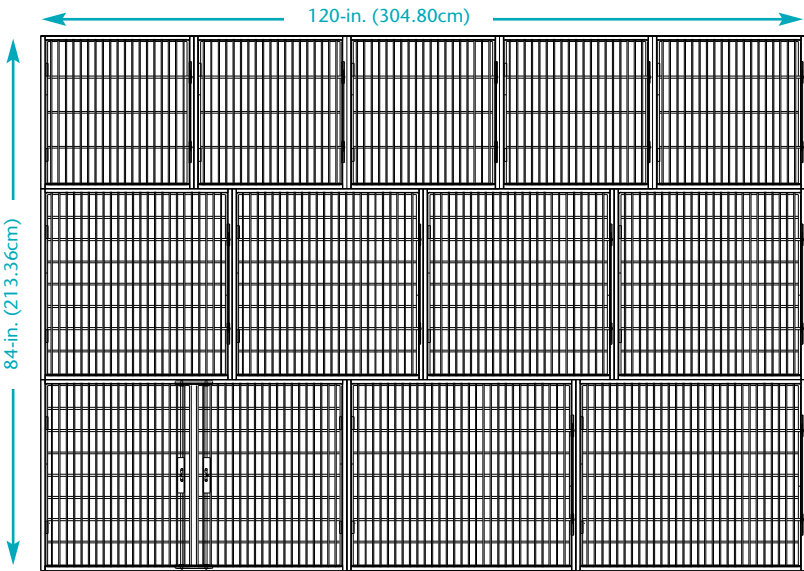


P/N	Description
3000030028-___†	Aristocrat Cage Arrangement
3000033012-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

10-Foot (3.05m) Wide, 12-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Five 24-in.W x 30-in.H (60.96cmW x 60.96cmH) cages
- Middle:** Four 30-in.W x 24-in.H (76.20cmW x 76.20cmH) cages
- Bottom:** One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage
Two 36-in.W x 30-in.H (91.44cmW x 76.20cmH) cages

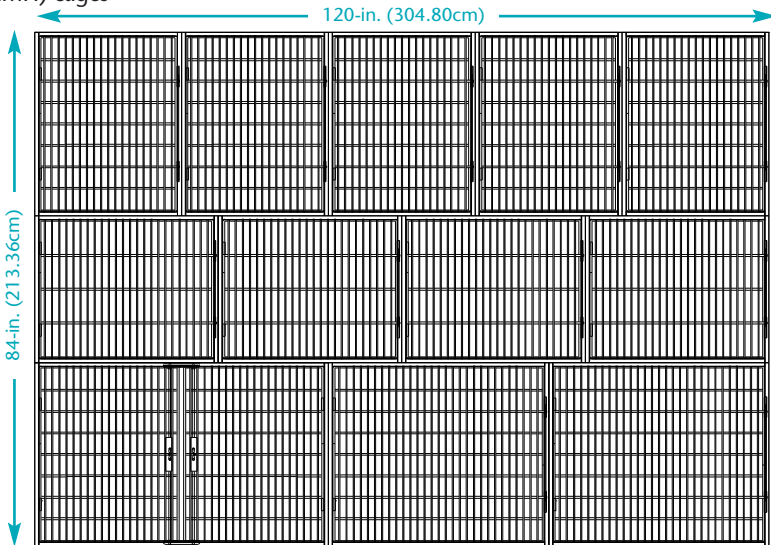

P/N	Description
3000030029-___†	Aristocrat Cage Arrangement
3000033013-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:
†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.


Aristocrat Cage Arrangements

12-Foot (3.66m) Wide, 14-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Six 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
Middle: Four 30-in.W x 30-in.H (76.20cmW x 76.20cmH) cages
 One 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cages
Bottom: Three 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage

P/N	Description
3000030030-___†	Aristocrat Cage Arrangement
3000033014-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

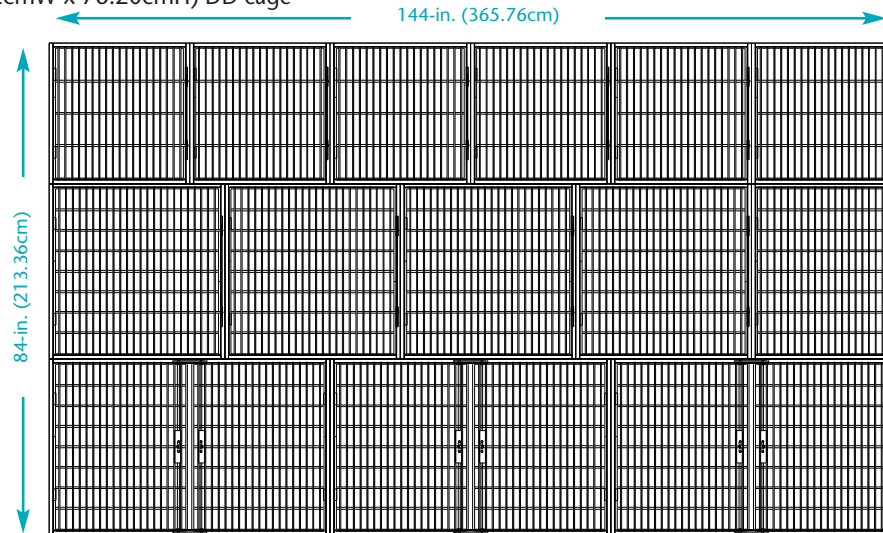


14-Foot (4.27m) Wide, 17-Cage Aristocrat Cage Arrangement

Includes:

- Top:** Seven 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
Middle: Four 30-in.W x 30-in.H (76.20cmW x 76.20cmH) cages
 Two 24-in.W x 30-in.H (60.96cmW x 76.20cmH) cages
Bottom: Two 36-in.W x 30-in.H (121.92cmW x 76.20cmH) cages
 Two 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cages

P/N	Description
3000030031-___†	Aristocrat Cage Arrangement
3000033015-___†	Full-View Polycarbonate Aristocrat Cage Arrangement

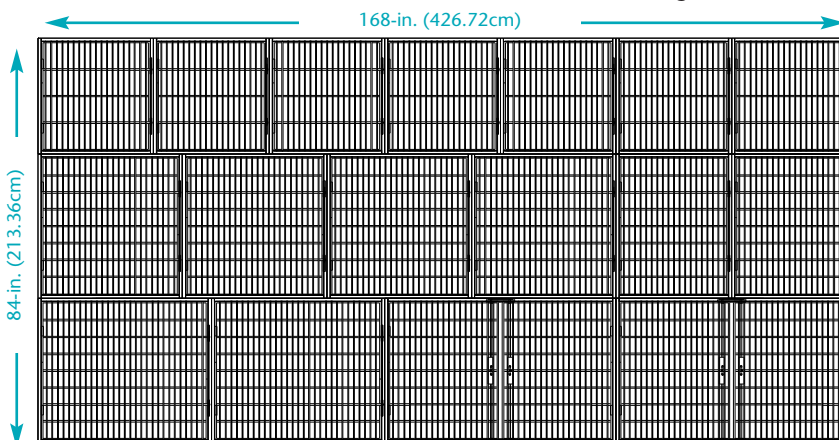

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.


Aristocrat Corner Cage Arrangements

4-Foot (2.44m), 2-Foot Corner (1.22m), 4-Foot (2.44m), 8-Cage Aristocrat Corner Cage Arrangement

Includes:

- Top:** Two 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
One 24-in.W x 24-in.H (60.96cmW x 60.96cmH) corner cage
Two 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Bottom:** One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage
One 24in.W x 30-in.H (60.96cmW x 76.20cmH) corner cage
One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage

P/N	Description
3000030011____†	Aristocrat Corner Cage Arrangement
3000033016____†	Full-View Polycarbonate Aristocrat Corner Cage Arrangement

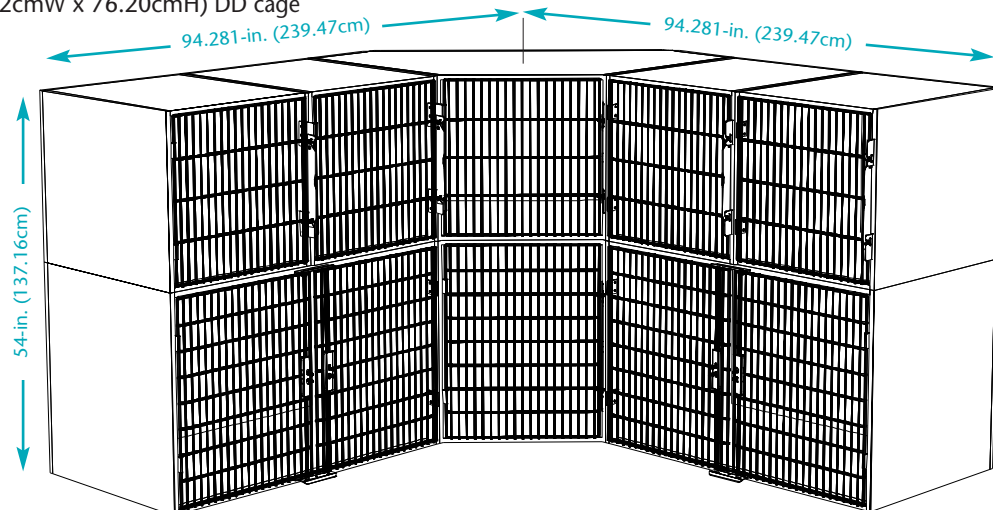


6-Foot (3.05m), 2-Foot (1.22m) Corner, 4-Foot (2.44m), 10-Cage Aristocrat Corner Cage Arrangement

Includes:

- Top:** Three 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
One 24-in.W x 24-in.H (60.96cmW x 60.96cmH) corner cage
Two 24-in.W x 24-in.H (60.96cmW x 60.96cmH) cages
- Bottom:** Two 36-in.W x 30-in.H (91.44cmW x 76.20cmH) cages
One 24in.W x 30-in.H (60.96cmW x 76.20cmH) corner cage
One 48-in.W x 30-in.H (121.92cmW x 76.20cmH) DD cage

P/N	Description
3000030032-____†	Aristocrat Cage Arrangement
3000033017-____†	Full-View Polycarbonate Aristocrat Cage Arrangement

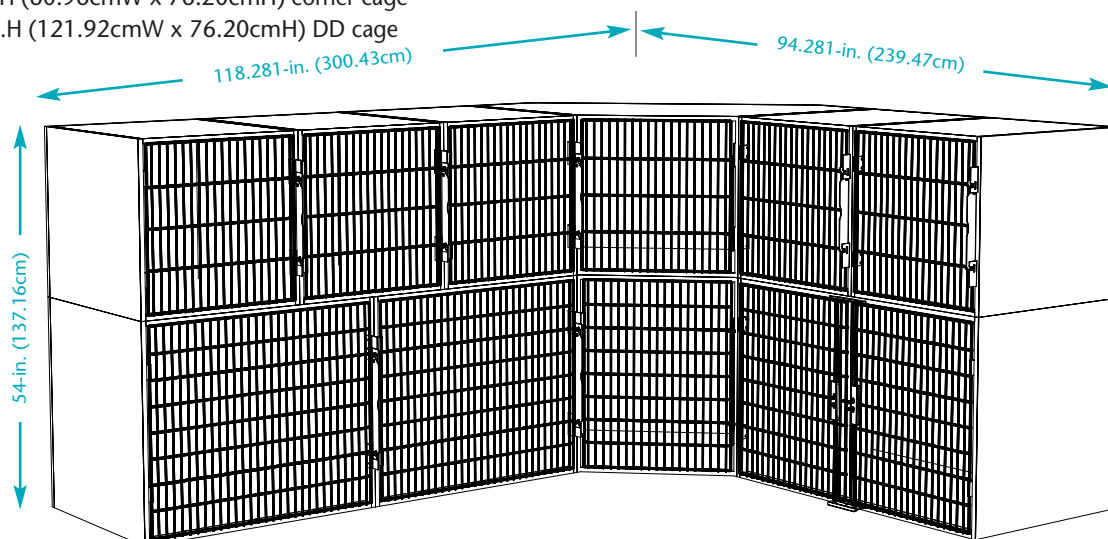

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Finishing Accessories For Aristocrat Cage

Finishing Accessories and Cage Accessories, give your Aristocrat Cage Arrangement a finished appearance and greater versatility. Aristocrat Cage Finishing Accessories will only fit SSCI Aristocrat Cages. SSCI Cage Accessories fit all types of SSCI Cages.



**Arrangement above shown
with Top and End Panels**




**Arrangement above shown
without Top and End Panels**

Finishing Accessories For Aristocrat Cage

When the Divider Panel is removed, there's no center post to get in the way when placing animals in the cage.

Stainless steel Divider Panel slides securely into rear channel.

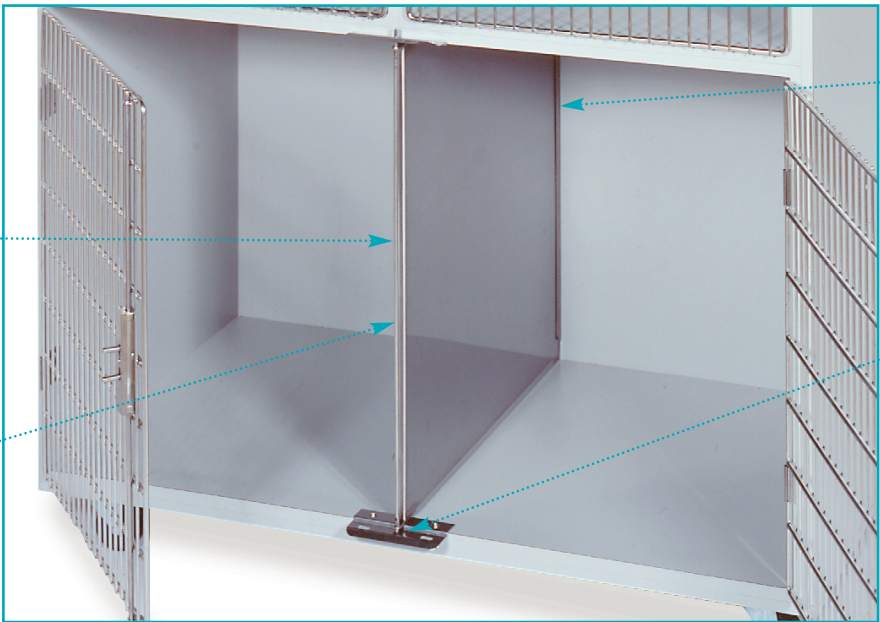

Stainless steel channel holds Divider Panel in place.

Front pin slips into specially designed upper and lower door latches.

Removable Divider Panel for Double Door Aristocrat Cages

Use this rugged stainless steel Divider to turn one large cage into two smaller ones. Panel slides into Channel included with all Double Door Aristocrat Cages, and the self-latching front pin fits into top and bottom cage-door latches to hold panel securely in place. Panel removes quickly and easily when not required.


P/N	FITS CAGE SIZE	
	Inches	Centimeters
12111-00-AADREI	24H	60.96H
12111-00-AAEPEI	30H	76.20H
12111-00-AAFNEI	36H	91.44H

Laminated Resting Shelves

Make your cages more versatile with our laminated Resting Shelves. Shelf can be used to make room for multiple animals in one cage, or simply to make cats more comfortable by allowing them to climb and rest higher than the floor. Shelves, factory-installed, are available when ordering new cages, or shelves may be ordered separately for customer installation in new or existing cages. Shelves are mechanically-fastened thru the sides of the cage for stability. This prevents the shelf from ever tipping out of position. Existing arrangements must be disassembled to install shelves. Three types of shelves are available.

- Across the back shelves, 6" deep, are positioned against the back spanning the width of the cage.
- Corner back shelves, 12" deep x 12" wide, are positioned in the back and side of the left or right corner of the cage.
- Side shelves, 6" wide x 27" deep, are positioned on the left or right side spanning the depth of the cage.

Typically, the shelves are installed approximately midway between the floor and the ceiling of the cage.

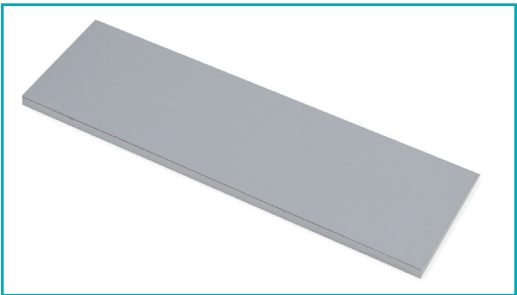

SIZE/P/N	DESCRIPTION	FITS CAGE	
		In.	Cm
14165-00-CTACAX-__-†	Back Resting Shelf	18W	45.72W
14165-00-DRACAX-__-†	Back Resting Shelf	24W	60.96W
14165-00-EPACAX-__-†	Back Resting Shelf	30W	76.20W
14165-00-FNACAX-__-†	Back Resting Shelf	36W	91.44W
14166-00-BVACBV-__-†	Corner shelf	Fits any cage	
14167-00-AXACED-__-†	Side Resting Shelf	Fits any cage	

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Finishing Accessories For Aristocrat Cage



Designed to attach easily to your Aristocrat Cage Arrangements, laminated Top and End Panels are laminated in your choice of colors to match or coordinate with your Arrangement. Without these panels the outside of the arrangements are unfinished. Top and End Panels are required to give cage arrangements a finished appearance.

A Top Panels

Order the combination of panels you need to span the entire top of your arrangement. For example: for a 48-in. (121.92cm) wide arrangement, you need to order two 24-in. (60.96cm) wide panels or one 48-in. (121.92cm) wide panel.

P/N	WIDTH	
	Inches	Centimeters
Top Panel for Aristocrat Cages		
14140-00-CTACEI-___†	18W	45.72W
14140-00-DRACEI-___†	24W	60.96W
14140-00-EPACEI-___†	30W	76.20W
14140-00-FNACEI-___†	36W	91.44W
14140-00-GLACEI-___†	42W	106.68W
14140-00-HJACEI-___†	48W	121.92W
14140-00-JFACEI-___†	60W	152.40W
14140-00-LBACEI-___†	72W	182.88W
Top Panel for Corner Aristocrat Cages		
14141-00-DRACEI-___†	Square Back Corner	
14142-00-DRACEI-___†	Mitered Back Corner*	

* Required for Kitty Suite Module Arrangements with exposed back.

B End Panels

Order the appropriate height End Panel for each exposed end of each tier of cages.

P/N	DIMENSIONS	
	Inches	Centimeters
End Panels for Arrangements Without Top Panel		
14120-00-ACCTEI-___†	18H	45.72H
14120-00-ACDREI-___†	24H	60.96H
14120-00-ACEPEI-___†	30H	76.20H
14120-00-ACFNEI-___†	36H	91.44H
Top Tier End Panels for Arrangements With Top Panel		
14121-00-ACCWEI-___†	18.625H	47.31H
14121-00-ACDUEI-___†	24.625H	62.55H
14121-00-ACESEI-___†	30.625H	77.79H
14121-00-ACFQEI-___†	36.625H	93.03H

Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

For Product Information Call TOLL-FREE 1-800-323-7366 • (847) 537-9320 • Fax (847) 537-9061

Hours: 8:30am - 5:00pm CST Monday - Friday

Visit our website at: www.suburbansurgical.com

Catalog No. 702488 12/20/07

Finishing Accessories For Aristocrat Cage

Toe-Kick Bases

Laminated Toe-Kick Bases give Aristocrat Cage Arrangement stability. Front Toe Kick Bases are used with arrangements with standard-style modules. Front & Back Toe Kick Bases are used with arrangements with Full-View Back Cages. Base's sides and toe-kicks are laminated.

Base adds just 4-in. (10.16cm) to Cage Arrangement's height .



Laminated Base

P/N Toe Kick Front	P/N Toe Kick Front & Back	Feet	Meters
14115- 00-CTAPEI- _ _ _ †	14117- 00-CTAPEI- _ _ _ †	1'6"W	.46W
14115- 00-DRAPEI- _ _ _ †	14117- 00-DRAPEI- _ _ _ †	2'W	.61W
14115- 00-EPAPEI- _ _ _ †	14117- 00-EPAPEI- _ _ _ †	2'6"W	.76W
14115- 00-FNAPEI- _ _ _ †	14117- 00-FNAPEI- _ _ _ †	3'W	.91W
14115- 00-GLAPEI- _ _ _ †	14117- 00-GLAPEI- _ _ _ †	3'6"W	1.07W
14115- 00-HJAPEI- _ _ _ †	14117- 00-HJAPEI- _ _ _ †	4'W	1.22W
14115- 00-IHAPEI- _ _ _ †	14117- 00-IHAPEI- _ _ _ †	4'6"W	1.37W
14115- 00-JFAPEI- _ _ _ †	14117- 00-JFAPEI- _ _ _ †	5'W	1.52W
14115- 00-KDAPEI- _ _ _ †	14117- 00-KDAPEI- _ _ _ †	5'6"W	1.68W
14115- 00-LBAPEI- _ _ _ †	14117- 00-LBAPEI- _ _ _ †	6'W	1.83W

P/N Toe Kick Front	P/N Toe Kick Front & Back	Feet	Meters
14115- 00-LZAPEI- _ _ _ †	14117- 00-LZAPEI- _ _ _ †	6'6"W	1.98W
14115- 00-MXAPEI- _ _ _ †	14117- 00-MXAPEI- _ _ _ †	7'W	2.13W
14115- 00-NVAPEI- _ _ _ †	14117- 00-NVAPEI- _ _ _ †	7'6"W	2.29W
14115- 00-OTAPEI- _ _ _ †	14117- 00-OTAPEI- _ _ _ †	8'W	2.44W
14115- 00-PRAPEI- _ _ _ †	14117- 00-PRAPEI- _ _ _ †	8'6"W	2.59W
14115- 00-QPAPEI- _ _ _ †	14117- 00-QPAPEI- _ _ _ †	9'W	2.74W
14115- 00-RNAPEI- _ _ _ †	14117- 00-RNAPEI- _ _ _ †	9'6"W	2.90W
14115- 00-SLAPEI- _ _ _ †	14117- 00-SLAPEI- _ _ _ †	10'W	3.05W
14116-00-DRAPEI- _ _ _ †	14118-00-DRAPEI- _ _ _ †	For Corner Cages	

Stainless Steel Stationary Bases

Stainless Steel Bases give your Cage Arrangement stability. Leg levelers ensure your Arrangement can be level even if your floor is slightly uneven.

Base adds just 1.5-in. (3.81 cm) to Cage Arrangement's height .



P/N STATIONARY BASE	FITS CAGE ARRANGEMENT Feet	Meters
12117-00-CTAHEI	1'6"W	.46W
12117-00-DRAHEI	2'W	.61W
12117-00-EPAHEI	2'6"W	.76W
12117-00-FNAHEI	3'W	.91W
12117-00-GLAHEI	3'6"W	1.07W
12117-00-HJAHEI	4'W	1.22W
12117-00-IHAHEI	4'6"W	1.37W
12117-00-JFAHEI	5'W	1.52W
12117-00-KDAHEI	5'6"W	1.68W
12117-00-LBAHEI	6'W	1.83W
12117-00-LZAHEI	6'6"W	1.98W
12117-00-MXAHEI	7'W	2.13W
12117-00-NVAHEI	7'6"W	2.29W
12117-00-OTAHEI	8'W	2.44W
12117-00-PRAHEI	8'6"W	2.59W
12117-00-QPAHEI	9'W	2.74W
12117-00-RNAHEI	9'6"W	2.90W
12117-00-SLAHEI	10'W	3.05W
12118-50-DRDREI for Corner Cage	2'1.500"W	.65W

Stainless Steel Mobile Bases

Stainless Steel Mobile Base give your Cage Arrangement mobility. Base adds 6.75-in. (17.15cm) to Cage Arrangement's height. Swivel casters with 5-in. (12.70cm) diameter wheels allow for easy movement throughout your facility.

Bases up to 6-ft. (1.83m) long have four casters, two with wheel brakes. Bases 6-ft. 6-in. (1.98m) long or longer have six casters, three with wheel brakes.



Please Note: Mobile bases are not recommended for L-shaped Aristocrat Cage Arrangements.

P/N MOBILE BASE	FITS CAGE ARRANGEMENT Feet	Meters
12117-50-CTAHEI	1'6"W	.46W
12117-50-DRAHEI	2'W	.61W
12117-50-EPAHEI	2'6"W	.76W
12117-50-FNAHEI	3'W	.91W
12117-50-GLAHEI	3'6"W	1.07W
12117-50-HJAHEI	4'W	1.22W
12117-50-IHAHEI	4'6"W	1.37W
12117-50-JFAHEI	5'W	1.52W
12117-50-KDAHEI	5'6"W	1.68W
12117-50-LBAHEI	6'W	1.83W
12117-50-LZAHEI	6'6"W	1.98W
12117-50-MXAHEI	7'W	2.13W
12117-50-NVAHEI	7'6"W	2.29W
12117-50-OTAHEI	8'W	2.44W
12117-50-PRAHEI	8'6"W	2.59W
12117-50-QPAHEI	9'W	2.74W
12117-50-RNAHEI	9'6"W	2.90W
12117-50-SLAHEI	10'W	3.05W
12118-50-DRDREI for Corner Cage	2'1.500"W	.65W


Important: These customer order requirements, in addition to other customer order requirements, must be satisfied prior to the beginning of manufacturing lead time:

†Specify 4-digit Color code. See Color Selections, Section A, pages 1 & 2 for Color Codes.



SSCI laminated products are available in SSCI standard laminates, SSCI optional laminates, or Customer selected laminates.

Cage Accessories will fit all SSCI Cages.

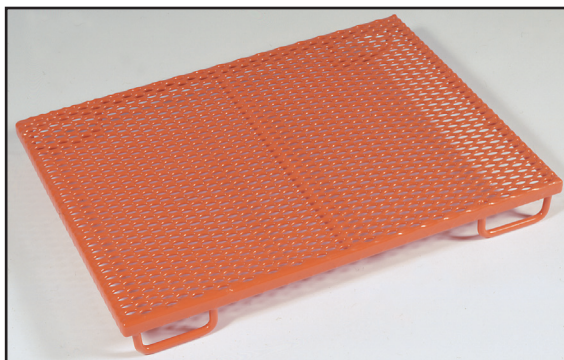


Plastisol-Coated Removable Cage Floors

Orange plastisol-coated expanded metal Cage Floor. Allows animals to rest 3 -in. (7.62cm) above the cage bottom to stay clean, dry and comfortable. For double-door cages, order 2 Cage Floors of the appropriate size.

P/N	FITS CAGES (W)		FLOOR SIZE (WxHxD)	
	In.	Cm	Inches	Centimeters
12385-00-CTALEI	18W	45.72W	15.438W x 3H x 28D	38.98W x 7.62H x 71.12D
12385-00-DRALEI	24W	60.96W	21.438W x 3H x 28D	54.45W x 7.62H x 71.12D
12385-00-EPALEI	30W	76.20W	27.438W x 3H x 28D	69.69W x 7.62H x 71.12D
12385-00-FNALEI	36W	91.44W	33.438W x 3H x 28D	84.93W x 7.62H x 71.12D

Stainless Steel Removable Cage Floors

Round rod stainless steel Cage Floor is welded at all intersections for superior strength and durability. Rods are deburred and spaced .750-in. (1.90cm) apart and have electropolished finish. Allows animals to rest 1.5 -in. (3.81cm) above the cage bottom to stay clean, dry and comfortable. For double-door cages, order 2 Cage Floors of the appropriate size.

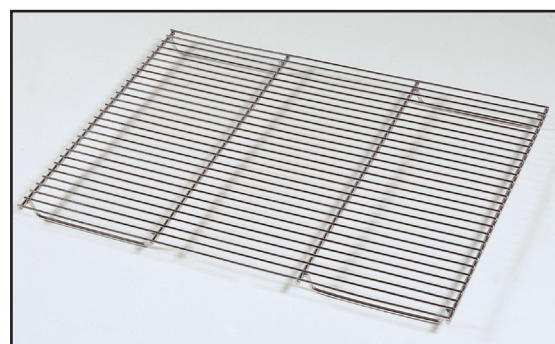

P/N	FITS CAGES (W)		FLOOR SIZE (WxHxD)	
	In.	Cm	Inches	Centimeters
12375-00-CTAFEI	18W	45.72W	15.188W x 1.412H x 27.75D	38.58W x 3.59H x 70.49D
12375-00-DRAFEI	24W	60.96W	21.188W x 1.412H x 27.75D	53.82W x 3.59H x 70.49D
12375-00-EPAFEI	30W	76.20W	27.188W x 1.412H x 27.75D	69.06W x 3.59H x 70.49D
12375-00-FNAFEI	36W	91.44W	33.188W x 1.412H x 27.75D	84.30W x 3.59H x 70.49D



Tapered Litter Pan

Our stainless steel Litter Pan is just the right size to fit neatly into your cage. Features tapered sides for nesting. Easily cleaned and disinfected.

P/N	DESCRIPTION	MEASUREMENTS (WxHxD)	
		Inches	Centimeters
12350-00-BZAPBJ	Tapered Litter Pan	13W x 4H x 9D	33.02W x 10.16H x 22.86D

Cage Accessories



Clamp-On Adjustable IV Hook

Stainless steel Clamp-on IV Hook attaches securely to vertical spokes in cage door. Height adjusts total amount shown below.

P/N	HEIGHT ADJUSTS	
	In.	Cm
12397-00-AACXAA	19H	48.26H
12397-00-AAFJAA	35H	88.90H



Food and Water Bowls

Stainless steel bowls are easy to clean and sanitize to ensure animal health.

P/N	CAPACITY	
	U.S. Liquid	Metric
12305-00-A2AAAA	1 Pint	0.47 Liters
12305-00-A4AAAA	1 Quart	0.95 Liters
12305-00-A6AAAA	2 Quarts	1.89 Liters



Food and Water Bowl Holders

Attach securely to cage door to hold stainless steel bowls sold above. Designed to keep bowl off the floor and eliminate spills. Made of stainless steel with electropolished finish.

P/N	FITS BOWL SIZE	
	U.S. Liquid	Metric
12335-00-A2AAAA	1 Pint	0.47 Liters
12335-00-A4AAAA	1 Quart	0.95 Liters
12335-00-A6AAAA	2 Quarts	1.89 Liters





Card Holders

Use to display animal's name, weight, medication needs, release date, or other pertinent information. Attaches securely to the cage door to hold identification cards on three sides. Made of durable, easy to sanitize stainless steel. Cards slide into holder from top.

P/N	DIMENSIONS	
	Inches	Centimeters
12300-00-ATALAA	5W x 3H	12.70W x 7.62H
12300-00-AXAPAA	6W x 4H	15.24W x 10.16H
12300-00-BFATAA	8W x 5H	20.32W x 12.70H

Clipboards

Holds information about the animal in your care. Yet lifts off quickly and easily to use as a writing surface to update the chart. Clip holds chart securely to keep all pertinent information at hand and easily accessible. Made of durable, easy to sanitize stainless steel.

P/N	DIMENSIONS	
	Inches	Centimeters
12302-00-BJAXAA	9W x 6H	22.86W x 15.24H
12302-00-BJBVAA	9W x 12H	22.86W x 30.48H



Clipboard with Bin

Bin/clipboard combination lifts on and off easily, lets you keep medicines, special foods or patient's personal belongings, leash, collar, etc. and dispensing schedule right on patient's cage. Helps to remind you when the next dose is due. Made of durable, easy to sanitize stainless steel.

P/N	BIN HEIGHT		OVERALL DIMENSIONS	
	In.	Cm	In.	Cm
12304-00-AXBIAL	4.188H	10.64H	6W x 9H x 3D	15.24W x 22.86H x 7.62D
12304-00-BJBVAL	5.188H	13.18H	9W x 12H x 3D	22.86W x 30.48H x 7.62D

